

# **Lab Manual for Physics 206/211**

Benjamin Crowell  
Fullerton College

[www.lightandmatter.com](http://www.lightandmatter.com)

Copyright (c) 1999-2011 by B. Crowell. This lab manual is subject to the Creative Commons CC-BY-SA license.

# Contents

1	Electricity . . . . .	6
2	Resistivity . . . . .	12
3	The Loop and Junction Rules . . . . .	14
4	Electric Fields . . . . .	16
5	Relativity . . . . .	20
6	Magnetism . . . . .	24
7	Electromagnetism . . . . .	28
8	The Charge to Mass Ratio of the Electron. . . . .	34
9	Radioactivity . . . . .	38
10	Two-Source Interference . . . . .	42
11	Refraction and Images . . . . .	44
12	Geometric Optics. . . . .	48
13	Wave Optics . . . . .	52
14	The Photoelectric Effect . . . . .	56
15	Electron Diffraction . . . . .	60
16a	Setup of the Spectrometer . . . . .	64
16b	The Mass of the Electron . . . . .	68
16c	The Nitrogen Molecule. . . . .	74
	Appendix 1: Format of Lab Writeups . . . . .	80
	Appendix 2: Basic Error Analysis. . . . .	82
	Appendix 3: Propagation of Errors . . . . .	88
	Appendix 4: Graphing . . . . .	90
	Appendix 5: Finding Power Laws from Data . . . . .	94
	Appendix 6: High Voltage Safety Checklist . . . . .	96
	Appendix 7: Laser Safety Checklist . . . . .	98





# 1 Electricity

## Apparatus

scotch tape  
rubber rod  
heat lamp  
fur  
bits of paper  
rods and strips of various materials  
30-50 cm rods, and angle brackets, for hanging charged rods  
power supply (Thornton), in lab benches ..1/group  
multimeter (PRO-100), in lab benches .... 1/group  
alligator clips  
flashlight bulbs  
spare fuses for multimeters — Let students replace fuses themselves.

## Goals

Determine the qualitative rules governing electrical charge and forces.

Light up a lightbulb, and measure the current through it and the voltage difference across it.

## Introduction

Newton's law of gravity gave a mathematical formula for the gravitational force, but his theory also made several important non-mathematical statements about gravity:

Every mass in the universe attracts every other mass in the universe.

Gravity works the same for earthly objects as for heavenly bodies.

The force acts at a distance, without any need for physical contact.

Mass is always positive, and gravity is always attractive, not repulsive.

The last statement is interesting, especially because it would be fun and useful to have access to some

negative mass, which would fall up instead of down (like the “upsydaisium” of Rocky and Bullwinkle fame).

Although it has never been found, there is no theoretical reason why a second, negative type of mass can't exist. Indeed, it is believed that the nuclear force, which holds quarks together to form protons and neutrons, involves three qualities analogous to mass. These are facetiously referred to as “red,” “green,” and “blue,” although they have nothing to do with the actual colors. The force between two of the same “colors” is repulsive: red repels red, green repels green, and blue repels blue. The force between two different “colors” is attractive: red and green attract each other, as do green and blue, and red and blue.

When your freshly laundered socks cling together, that is an example of an electrical force. If the gravitational force involves one type of mass, and the nuclear force involves three colors, how many types of electrical “stuff” are there? In the days of Benjamin Franklin, some scientists thought there were two types of electrical “charge” or “fluid,” while others thought there was only a single type. In the first part of this lab, you will try to find out experimentally how many types of electrical charge there are.

The unit of charge is the coulomb, C; one coulomb is defined as the amount of charge such that if two objects, each with a charge of one coulomb, are one meter apart, the magnitude of the electrical force between them is  $9 \times 10^9$  N. Practical applications of electricity usually involve an electric circuit, in which charge is sent around and around in a circle and recycled. Electric current,  $I$ , measures how many coulombs per second flow past a given point; a shorthand for units of C/s is the ampere, A. Voltage,  $V$ , measures the electrical potential energy per unit charge; its units of J/C can be abbreviated as volts, V. Making the analogy between electrical interactions and gravitational ones, voltage is like height. Just as water loses gravitational potential energy by going over a waterfall, electrically charged particles lose electrical potential energy as they flow through a circuit. The second part of this lab involves building an electric circuit to light up a lightbulb, and measuring both the current that flows through the bulb and the voltage difference across it.

## Observations

The following important rules serve to keep facts separate from opinions and reduce the chances of getting a garbled copy of the data:

(1) Take your raw data in pen, directly into your lab notebook. This is what real scientists do. The point is to make sure that what you're writing down is a first-hand record, without mistakes introduced by recopying it. (If you don't have your two lab notebooks yet, staple today's raw data into your notebook when you get it.)

(2) Everybody should record their own copy of the raw data. Do not depend on a "group secretary."

(3) If you do calculations during lab, keep them on a separate page or draw a line down the page and keep calculations on one side of the line and raw data on the other. This is to distinguish facts from inferences. (I will deduct 25% from your grade if you mix calculations and raw data.)

(4) Never write numbers without units. Without units, a number is meaningless. There is a big difference between "Johnny is six" and "Johnny is six feet." (I will deduct 25% from your grade if you write numbers without units.)

Because this is the first meeting of the lab class, there is no prelab writeup due at the beginning of the class. Instead, you will discuss your results with your instructor at various points.

### A Inferring the rules of electrical repulsion and attraction

Stick a piece of scotch tape on a table, and then lay another piece on top of it. Pull both pieces off the table, and then separate them. If you now bring them close together, you will observe them exerting a force on each other. Electrical effects can also be created by rubbing the fur against the rubber rod.

Your job in this lab is to use these techniques to test various hypotheses about electric charge. The most common difficulty students encounter is that the charge tends to leak off, especially if the weather is humid. If you have charged an object up, you should not wait any longer than necessary before making your measurements. It helps if you keep your hands dry.

To keep this lab from being too long, the class will pool its data for part A. Your instructor will organize the results on the whiteboard.

*i. Repulsion and/or attraction*

Test the following hypotheses. Note that they are mutually exclusive, i.e., only one of them can be true.

A) Electrical forces are always attractive.

R) Electrical forces are always repulsive.

AR) Electrical forces are sometimes attractive and sometimes repulsive.

Interpretation: Once the class has tested these hypotheses thoroughly, we will discuss what this implies about how many different types of charge there might be.

*ii. Are there forces on objects that have not been specially prepared?*

So far, special preparations have been necessary in order to get objects to exhibit electrical forces. These preparations involved either rubbing objects against each other (against resistance from friction) or pulling objects apart (e.g. overcoming the sticky force that holds the tape together). In everyday life, we do not seem to notice electrical forces in objects that have not been prepared this way.

Now try to test the following hypotheses. Bits of paper are a good thing to use as unprepared objects, since they are light and therefore would be easily moved by any force. *Do not* use tape as an uncharged object, since it can become charged a little bit just by pulling it off the roll.

U0) Objects that have not been specially prepared are immune to electrical forces.

UA) Unprepared objects can participate in electrical forces with prepared objects, and the forces involved are always attractive.

UR) Unprepared objects can participate in electrical forces with prepared objects, and the forces involved are always repulsive.

UAR) Unprepared objects can participate in electrical forces with prepared objects, and the forces involved can be either repulsive or attractive.

These four hypotheses are mutually exclusive.

Once the class has tested these hypotheses thoroughly, we will discuss what practical implications this has for planning the observations for part iii.

*iii. Rules of repulsion and/or attraction and the number of types of charge*

Test the following mutually exclusive hypotheses:

1A) There is only one type of electric charge, and the force is always attractive.

1R) There is only one type of electric charge, and the force is always repulsive.

2LR) There are two types of electric charge, call them X and Y. Like charges repel (X repels X and Y repels Y) and opposite charges attract (X and Y attract each other).

2LA) There are two types of electric charge. Like charges attract and opposite charges repel.

3LR) There are three types of electric charge, X, Y and Z. Like charges repel and unlike charges attract.

On the whiteboard, we will make a square table, in which the rows and columns correspond to the different objects you're testing against each other for attraction and repulsion. To test hypotheses 1A through 3LR, you'll need to see if you can successfully explain your whole table by labeling the objects with only one label, X, or whether you need two or three.

Some of the equipment may look identical, but not be identical. In particular, some of the clear rods have higher density than others, which may be because they're made of different types of plastic, or glass. This could affect your conclusions, so you may want to check, for example, whether two rods with the same diameter, that you think are made of the same material, actually weigh the same.


In general, you will find that some materials, and some combinations of materials, are more easily charged than others. For example, if you find that the mahogany rod rubbed with the weasel fur doesn't charge well, then don't keep using it! The white plastic strips tend to work well, so don't neglect them.

Once we have enough data in the table to reach a definite conclusion, we will summarize the results from part A and then discuss the following examples of incorrect reasoning about this lab.

- (1) "The first piece of tape exerted a force on the second, but the second didn't exert one on the first."
- (2) "The first piece of tape repelled the second, and the second attracted the first."
- (3) "We observed three types of charge: two that exert forces, and a third, neutral type."
- (4) "The piece of tape that came from the top was positive, and the bottom was negative."
- (5) "One piece of tape had electrons on it, and the other had protons on it."
- (6) "We know there were two types of charge, not three, because we observed two types of interactions, attraction and repulsion."

## B Measuring current and voltage

As shown in the figure, measuring current and voltage requires hooking the meter into the circuit in two completely different ways.

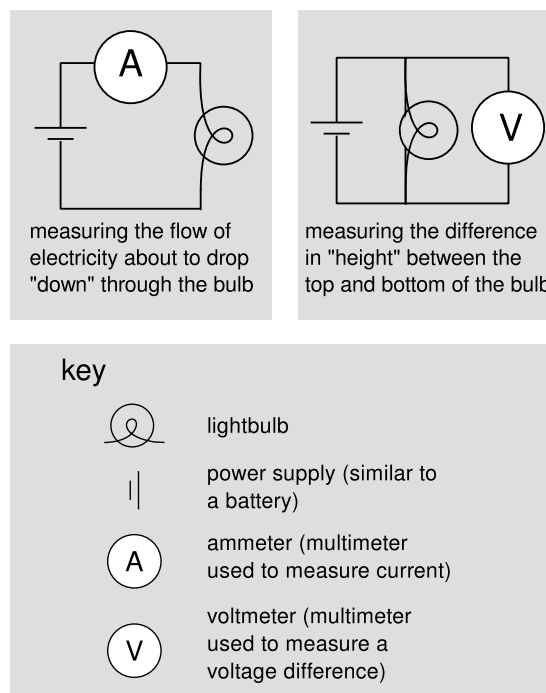

The arrangement for the ammeter is called a series circuit, because every charged particle that travels the circuit has to go through each component in a row, one after another. The series circuit is arranged like beads on a necklace.

The setup for the voltmeter is an example of a parallel circuit. A charged particle flowing, say, clockwise around the circuit passes through the power supply and then reaches a fork in the road, where it has a choice of which way to go. Some particles will pass through the bulb, others (not as many) through the meter; all of them are reunited when they reach the junction on the right.

Students tend to have a mental block against setting up the ammeter correctly in series, because it involves breaking the circuit apart in order to insert the meter. To drive home this point, we will act out the process using students to represent the circuit components. If you hook up the ammeter incorrectly, in parallel rather than in series, the meter provides an easy path for the flow of current, so a large amount of current will flow. To protect the meter from this surge, there is a fuse inside, which will blow, and the meter will stop working. This is not a huge tragedy; just ask your instructor for a replacement fuse and open up the meter to replace



it.

Unscrew your lightbulb from its holder and look closely at it. Note that it has two separate electrical contacts: one at its tip and one at the metal screw threads.

Turn the power supply's off-on switch to the off position, and turn its (uncalibrated) knob to zero. Set up the basic lightbulb circuit without any meter in it. There is a rack of cables in the back of the room with banana-plug connectors on the end, and most of your equipment accepts these plugs. To connect to the two brass screws on the lightbulb's base, you'll need to stick alligator clips on the banana plugs.

Check your basic circuit with your instructor, then turn on the power switch and *slowly* turn up the knob until the bulb lights. The knob is uncalibrated and highly nonlinear; as you turn it up, the voltage it produces goes zerozerozerozerozero*six*! To light the bulb without burning it out, you will need to find a position for the knob in the narrow range where it rapidly ramps up from 0 to 6 V.

Once you have your bulb lit, do not mess with the knob on the power supply anymore. You do not even need to switch the power supply off while rearranging the circuit for the two measurements with the meter; the voltage that lights the bulb is only about a volt or a volt and a half (similar to a battery), so it can't hurt you.

We have a single meter that plays both the role of the voltmeter and the role of the ammeter in this lab. Because it can do both these things, it is referred to as a multimeter. Multimeters are highly standardized, and the following instructions are generic ones that will work with whatever meters you happen to be using in this lab.

#### *Voltage difference*

Two wires connect the meter to the circuit. At the places where three wires come together at one point, you can plug a banana plug into the back of another banana plug. At the meter, make one connection at the "common" socket ("COM") and the other at the socket labeled "V" for volts. The common plug is called that because it is used for every measurement, not just for voltage.

Many multimeters have more than one scale for measuring a given thing. For instance, a meter may have a millivolt scale and a volt scale. One is used for measuring small voltage differences and the other for large ones. You may not be sure in advance what scale is appropriate, but that's not a big problem —

once everything is hooked up, you can try different scales and see what's appropriate. Use the switch or buttons on the front to select one of the voltage scales. By trial and error, find the most precise scale that doesn't cause the meter to display an error message about being overloaded.

Write down your measurement, with the units of volts, and stop for a moment to think about what it is that you've measured. Imagine holding your breath and trying to make your eyeballs pop out with the pressure. Intuitively, the voltage difference is like the pressure difference between the inside and outside of your body.

What do you think will happen if you unscrew the bulb, leaving an air gap, while the power supply and the voltmeter are still going? Try it. Interpret your observation in terms of the breath-holding metaphor.

#### *Current*

The procedure for measuring the current differs only because you have to hook the meter up in series and because you have to use the "A" (amps) plug on the meter and select a current scale.

In the breath-holding metaphor, the number you're measuring now is like the rate at which air flows through your lips as you let it hiss out. Based on this metaphor, what do you think will happen to the reading when you unscrew the bulb? Try it.

Discuss with your group and check with your instructor:

- (1) What *goes through* the wires? Current? Voltage? Both?
- (2) Using the breath-holding metaphor, explain why the voltmeter needs *two* connections to the circuit, not just one. What about the ammeter?

While waiting for your instructor to come around and discuss these questions with you, you can go on to the next part of the lab.

#### *Resistance*

The ratio of voltage difference to current is called the resistance of the bulb,  $R = \Delta V/I$ . Its units of volts per amp can be abbreviated as ohms,  $\Omega$  (capital Greek letter omega).

Calculate the resistance of your lightbulb. Resistance is the electrical equivalent of kinetic friction. Just as rubbing your hands together heats them up, objects that have electrical resistance produce heat when a current is passed through them. This is why the bulb's filament gets hot enough to heat up.

When you unscrew the bulb, leaving an air gap, what is the resistance of the air?

Ohm's law is a generalization about the electrical properties of a variety of materials. It states that the resistance is constant, i.e., that when you increase the voltage difference, the flow of current increases exactly in proportion. If you have time, test whether Ohm's law holds for your lightbulb, by cutting the voltage to half of what you had before and checking whether the current drops by the same factor. (In this condition, the bulb's filament doesn't get hot enough to create enough visible light for your eye to see, but it does emit infrared light.)

#### *List of materials for static electricity*

You don't have to know anything about what the various materials are in order to do this lab, but here is a list for use by instructors and the lab technician:

- scotch tape (used as two different objects, top and bottom)
- teflon fabric (brown, coarse)
- teflon rods (white, rigid, slippery, skinny)
- PVC pipe
- polyurethane rods (brown, flexible)
- nylon (?) fabric (blue)
- fur

## Notes For Next Week

(1) Next week, when you turn in your writeup for this lab, you also need to turn in a prelab writeup for the next lab. The prelab questions are listed at the end of the description of that lab in the lab manual. Never start a lab without understanding the answers to all the prelab questions; if you turn in partial answers or answers you're unsure of, discuss the questions with your instructor or with other students to make sure you understand what's going on.

(2) You should exchange phone numbers with your lab partners for general convenience throughout the semester. You can also get each other's e-mail addresses by logging in to Spotter and clicking on "e-mail."

## Rules and Organization

Collection of raw data is work you share with your lab partners. Once you're done collecting data, you need to do your own analysis. E.g., it is not okay for two people to turn in the same calculations, or on a lab requiring a graph for the whole group to make one graph and turn in copies.

You'll do some labs as formal writeups, others as informal "check-off" labs. As described in the syllabus, they're worth different numbers of points, and you have to do a certain number of each type by the end of the semester.

The format of formal lab writeups is given in appendix 1 on page 80. The raw data section must be contained in your bound lab notebook. Typically people word-process the abstract section, and any other sections that don't include much math, and stick the printout in the notebook to turn it in. The calculations and reasoning section will usually just consist of hand-written calculations you do in your lab notebook. You need two lab notebooks, because on days when you turn one in, you need your other one to take raw data in for the next lab. You may find it convenient to leave one or both of your notebooks in the cupboard at your lab bench whenever you don't need to have them at home to work on; this eliminates the problem of forgetting to bring your notebook to school.

For a check-off lab, the main thing I'll pay attention to is your abstract. The rest of your work for a check-off lab can be informal, and I may not ask to see it unless I think there's a problem after reading your abstract.



## 2 Resistivity

### Apparatus

power supply (Thornton), in lab benches ..1/group  
 multimeter (PRO-100), in lab benches ....1/group  
 digital multimeters (FlukeHP) .....1/group  
 resistors of various values  
 alligator clips  
 spare fuses for multimeters — Let students replace fuses themselves.  
 Play-Doh (Crayloa super soft dough, \$10/3 lb at school supply store on Commonwealth)  
 1 g aluminum slotted weights for use as electrical contacts

### Goals

Measure how the electrical resistance of a cylinder depends on its dimensions.

Invent a constant of proportionality incorporated into this relationship, giving a measure of the intrinsic electrical properties of a material called its resistivity.


### Introduction

Your nervous system depends on electrical currents, and every day you use many devices based on electrical currents without even thinking about it. Despite its ordinariness, the phenomenon of electric currents passing through liquids (e.g., cellular fluids) and solids (e.g., copper wires) is a subtle one. For example, we now know that atoms are composed of smaller, subatomic particles called electrons and nuclei, and that the electrons and nuclei are electrically charged, i.e., matter is electrical. Thus, we now have a picture of these electrically charged particles sitting around in matter, ready to create an electric current by moving in response to an externally applied voltage. Electricity had been used for practical purposes for a hundred years, however, before the electrical nature of matter was proved at the turn of the 20th century.

We observe that some materials, such as metals, conduct electricity more easily than others, such as wood. This suggests that measurements of the resistance of an object made out of a certain substance can be used to learn things about the atomic-

level structure of the substance (e.g., how many free charge carriers it has per unit volume), or as a test to identify an unknown substance.

color	meaning
black	0
brown	1
red	2
orange	3
yellow	4
green	5
blue	6
violet	7
gray	8
white	9
silver	$\pm 10\%$
gold	$\pm 5\%$

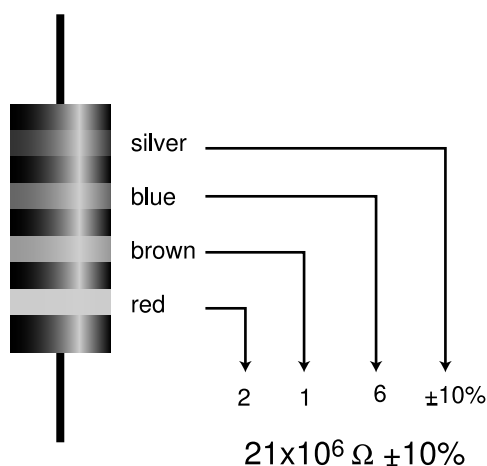

### Observations

#### A A known resistor

In lab 1, you found the resistance of a lightbulb. To keep things simple, that setup used only a sin-

gle multimeter, but that required rearranging how everything was hooked up every time you wanted to switch between measuring voltage and measuring current. In today's lab, you will use a similar setup, but with two meters simultaneously connected, one for each purpose. Build your setup, and test it by measuring data on a known resistor. To avoid creating large currents or making the resistor so hot that it will burn your fingers, use one with a fairly high value, at least in the kilohm range.

Resistors are usually too small to make it convenient to print numerical resistance values on them, so they are labeled with a color code, as shown in the table and example above.

As a test of whether your resistor is actually ohmic, take data at two different voltages and determine whether the current varies in proportion.

## B Dependence of resistance on dimensions

For a fixed substance such as sand or plywood, the resistance is not a fixed number of ohms. The resistance of a sample depends on its size, its shape, and the way in which the electrical contacts are made to it. In this part of the lab you will take measurements with cylindrical samples of play-doh, with the electrical contacts made using the thin disk-shaped 1 g aluminum weights.

Two different numbers are required in order to define the dimensions of a cylinder. (There are various ways to define these; you can make your own choice.) Measure the resistance for different values of these two variables. Practice control of variables by keeping one constant while varying the other. You will learn the most by taking extreme values of the variables.

Leaving the power supply on for a long time causes corrosion to the electrical contacts and chemical changes in the play-doh, so don't do that.

To get good results:

Quite a bit of the resistance comes from the surface of contact between the aluminum disk and the play-dough. What seems to matter the most is that you need a very smooth surface, because on a rough surface, most of the metal isn't even making contact.

It seems to help if you cut into the play-dough and then press the play-dough against the contact with your fingers to make better contact.

I got better results by subtracting the resistance found by placing the two contacts very

close, where the resistance was almost entirely due to the contacts.

## Analysis

For each variable in part B, use the graphing technique given in appendix 5 to see if you can find a power law that relates the resistance to that variable. If the results are at least approximately consistent with power laws, then you can combine them to create a relationship of the form  $R = \rho(\dots)^p(\dots)^q$ , where  $(\dots)$  represents the two variables,  $p$  and  $q$  are constants found from your graph, and the constant of proportionality  $\rho$  (Greek letter rho) is a property of the substance, called its resistivity. If you find that  $p$  and  $q$  have values that are close to nice simple numbers, carry out your analysis under the assumption that those simple numbers are right. What units does  $\rho$  have? Find the resistivity of play-doh.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** Check that you understand the interpretations of the following color-coded resistor labels:

blue gray orange silver	= 68 k $\Omega$ $\pm$ 10%
blue gray orange gold	= 68 k $\Omega$ $\pm$ 5%
blue gray red silver	= 6.8 k $\Omega$ $\pm$ 10%
black brown blue silver	= 1 M $\Omega$ $\pm$ 10%

Now interpret the following color code:

green orange yellow silver = ?

# 3 The Loop and Junction Rules

## Apparatus

DC power supply (Thornton) ..... 1/group  
multimeter (PRO-100, in lab benches) .... 1/group  
resistors

## Goal

Test the loop and junction rules in two electrical circuits.

## Introduction

If you ask physicists what are the most fundamentally important principles of their science, almost all of them will start talking to you about conservation laws. A conservation law is a statement that a certain measurable quantity cannot be changed. A conservation law that is easy to understand is the conservation of mass. No matter what you do, you cannot create or destroy mass.

The two conservation laws with which we will be concerned in this lab are conservation of energy and conservation of charge. Energy is related to voltage, because voltage is defined as  $V = PE/q$ . Charge is related to current, because current is defined as  $I = \Delta q/\Delta t$ .

Conservation of charge has an important consequence for electrical circuits:

When two or more wires come together at a point in a DC circuit, the total current entering that point equals the total current leaving it.

Such a coming-together of wires in a circuit is called a junction. If the current leaving a junction was, say, greater than the current entering, then the junction would have to be creating electric charge out of nowhere. (Of course, charge could have been stored up at that point and released later, but then it wouldn't be a DC circuit — the flow of current would change over time as the stored charge was used up.)

Conservation of energy can also be applied to an electrical circuit. The charge carriers are typically electrons in copper wires, and an electron has a potential energy equal to  $-eV$ . Suppose the electron sets off on a journey through a circuit made of re-

sistors. Passing through the first resistor, our subatomic protagonist passes through a voltage difference of  $\Delta V_1$ , so its potential energy changes by  $-e\Delta V_1$ . To use a human analogy, this would be like going up a hill of a certain height and gaining some gravitational potential energy. Continuing on, it passes through more voltage differences,  $-e\Delta V_2$ ,  $-e\Delta V_3$ , and so on. Finally, in a moment of religious transcendence, the electron realizes that life is one big circuit — you always end up coming back where you started from. If it passed through  $N$  resistors before getting back to its starting point, then the total change in its potential energy was

$$-e(\Delta V_1 + \dots + \Delta V_N) \quad .$$

But just as there is no such thing as a round-trip hike that is all downhill, it is not possible for the electron to have any net change in potential energy after passing through this loop — if so, we would have created some energy out of nothing. Since the total change in the electron's potential energy must be zero, it must be true that  $\Delta V_1 + \dots + \Delta V_N = 0$ . This is the loop rule:

The sum of the voltage differences around any closed loop in a circuit must equal zero.


When you are hiking, there is an important distinction between uphill and downhill, which depends entirely on which direction you happen to be traveling on the trail. Similarly, it is important when applying the loop rule to be consistent about the signs you give to the voltage differences, say positive if the electron sees an increase in voltage and negative if it sees a decrease along its direction of motion.

## Observations

### A The junction rule

Construct a circuit like the one in the figure, using the Thornton power supply as your voltage source. To make things more interesting, don't use equal resistors. Use resistors with values in the range of about 1 k $\Omega$  to 10 M $\Omega$ . If they're much higher than that, the currents will be too low for the PRO-100 meters to measure accurately. If they're much smaller than that, you could burn up the resistors, and the multimeter's internal resistance when used as an ammeter might not be negligible in comparison. Insert your multimeter in the circuit to measure all three

currents that you need in order to test the junction rule.

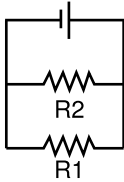

### B The loop rule

Now come up with a circuit to test the loop rule. Since the loop rule is always supposed to be true, it's hard to go wrong here! Make sure that (1) you have at least three resistors in a loop, (2) the whole circuit is not just a single loop, and (3) you hook in the power supply in a way that creates non-zero voltage differences across all the resistors. Measure the voltage differences you need to measure to test the loop rule. Here it is best to use fairly small resistances, so that the multimeter's large internal resistance when used in parallel as a voltmeter will not significantly reduce the resistance of the circuit. Do not use resistances of less than about  $100\ \Omega$ , however, or you may blow a fuse or burn up a resistor.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** Draw a schematic showing where you will insert the multimeter in the circuit to measure the currents in part A.

**P2** Invent a circuit for part B, and draw a schematic. You need not indicate actual resistor values, since you will have to choose from among the values actually available in lab.

**P3** Pick a loop from your circuit, and draw a schematic showing how you will attach the multimeter in the circuit to measure the voltage differences in part B.

**P4** Explain why the following statement is incorrect: "We found that the loop rule was not quite true, but the small error could have been because the resistor's value was off by a few percent compared to the color-code value."

## Self-Check

Do the analysis in lab.

## Analysis

Discuss whether you think your observations agree with the loop and junction rules, taking into account systematic and random errors.

## Programmed Introduction to Practical Electrical Circuits

The following practical skills are developed in this lab:

(1) Use a multimeter without being given an explicit schematic showing how to connect it to your circuit. This means connecting it in parallel in order to measure voltages and in series in order to measure currents.

(2) Use your understanding of the loop and junction rules to simplify electrical measurements. These rules often guarantee that you can get the same current or voltage reading by measuring in more than one place in a circuit. In real life, it is often much easier to connect a meter to one place than another, and you can therefore save yourself a lot of trouble using the rules.

# 4 Electric Fields

## Apparatus

board and U-shaped probe  
ruler  
DC power supply (Thornton)  
multimeter  
scissors  
stencils for drawing electrode shapes on paper

## Goals

To be better able to visualize electric fields and understand their meaning.

To examine the electric fields around certain charge distributions.

## Introduction

By definition, the electric field,  $E$ , at a particular point equals the force on a test charge at that point divided by the amount of charge,  $E = F/q$ . We can plot the electric field around any charge distribution by placing a test charge at different locations and making note of the direction and magnitude of the force on it. The direction of the electric field at any point P is the same as the direction of the force on a positive test charge at P. The result would be a page covered with arrows of various lengths and directions, known as a “sea of arrows” diagram..

In practice, Radio Shack does not sell equipment for preparing a known test charge and measuring the force on it, so there is no easy way to measure electric fields. What really is practical to measure at any given point is the voltage,  $V$ , defined as the electrical energy (potential energy) that a test charge would have at that point, divided by the amount of charge ( $E/Q$ ). This quantity would have units of J/C (Joules per Coulomb), but for convenience we normally abbreviate this combination of units as volts. Just as many mechanical phenomena can be described using either the language of force or the language of energy, it may be equally useful to describe electrical phenomena either by their electric fields or by the voltages involved.

Since it is only ever the difference in potential energy (interaction energy) between two points that

can be defined unambiguously, the same is true for voltages. Every voltmeter has two probes, and the meter tells you the difference in voltage between the two places at which you connect them. Two points have a nonzero voltage difference between them if it takes work (either positive or negative) to move a charge from one place to another. If there is a voltage difference between two points in a conducting substance, charges will move between them just like water will flow if there is a difference in levels. The charge will always flow in the direction of lower potential energy (just like water flows downhill).


All of this can be visualized most easily in terms of maps of constant-voltage curves (also known as equipotentials); you may be familiar with topographical maps, which are very similar. On a topographical map, curves are drawn to connect points having the same height above sea level. For instance, a cone-shaped volcano would be represented by concentric circles. The outermost circle might connect all the points at an altitude of 500 m, and inside it you might have concentric circles showing higher levels such as 600, 700, 800, and 900 m. Now imagine a similar representation of the voltage surrounding an isolated point charge. There is no “sea level” here, so we might just imagine connecting one probe of the voltmeter to a point within the region to be mapped, and the other probe to a fixed reference point very far away. The outermost circle on your map might connect all the points having a voltage of 0.3 V relative to the distant reference point, and within that would lie a 0.4-V circle, a 0.5-V circle, and so on. These curves are referred to as constant-voltage curves, because they connect points of equal voltage. In this lab, you are going to map out constant-voltage curves, but not just for an isolated point charge, which is just a simple example like the idealized example of a conical volcano.

You could move a charge along a constant-voltage curve in either direction without doing any work, because you are not moving it to a place of higher potential energy. If you do not do any work when moving along a constant-voltage curve, there must not be a component of electric force along the surface (or you would be doing work). A metal wire is a constant-voltage curve. We know that electrons in a metal are free to move. If there were a force along the wire, electrons would move because of it. In fact the electrons would move until they were distributed




in such a way that there is no longer any force on them. At that point they would all stay put and then there would be no force along the wire and it would be a constant-voltage curve. (More generally, any flat piece of conductor or any three-dimensional volume consisting of conducting material will be a constant-voltage region.)

There are geometrical and numerical relationships between the electric field and the voltage, so even though the voltage is what you'll measure directly in this lab, you can also relate your data to electric fields. Since there is not any component of electric force parallel to a constant-voltage curve, electric field lines always pass through constant-voltage curves at right angles. (Analogously, a stream flowing straight downhill will cross the lines on a topographical map at right angles.) Also, if you divide the work equation  $(\Delta \text{energy}) = Fd$  by  $q$ , you get  $(\Delta \text{energy})/q = (F/q)d$ , which translates into  $\Delta V = -Ed$ . (The minus sign is because  $V$  goes down when some other form of energy is released.) This means that you can find the electric field strength at a point  $P$  by dividing the voltage difference between the two constant-voltage curves on either side of  $P$  by the distance between them. You can see that units of  $\text{V/m}$  can be used for the  $E$  field as an alternative to the units of  $\text{N/C}$  suggested by its definition — the units are completely equivalent.

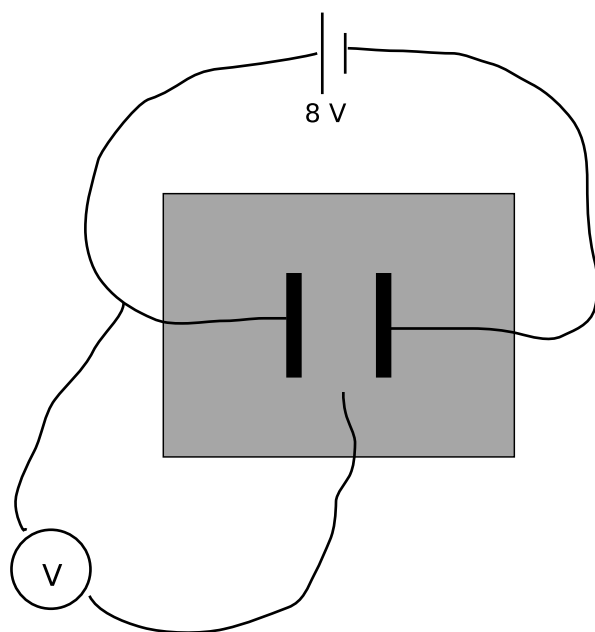

A simplified schematic of the apparatus, being used with pattern 1 on page 18.

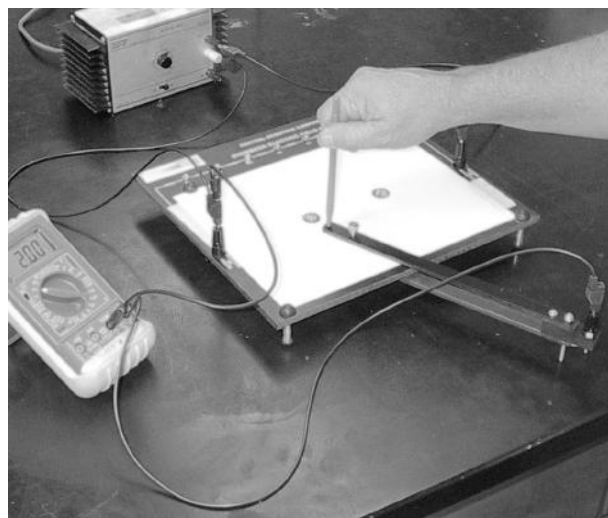

A photo of the apparatus, being used with pattern 3 on page 18.


## Method

The first figure shows a simplified schematic of the apparatus. The power supply provides an 8 V voltage difference between the two metal electrodes, drawn in black. A voltmeter measures the voltage difference between an arbitrary reference voltage and a point of interest in the gray area around the electrodes. The result will be somewhere between 0 and 8 V. A voltmeter won't actually work if it's not part of a complete circuit, but the gray area is intentionally made from a material that isn't a very good insulator, so enough current flows to allow the voltmeter to operate.

The photo shows the actual apparatus. The electrodes are painted with silver paint on a detachable board, which goes underneath the big board. What you actually see on top is just a piece of paper on which you'll trace the equipotentials with a pen. The voltmeter is connected to a U-shaped probe with a metal contact that slides underneath the board, and a hole in the top piece for your pen.

Turn your large board upside down. Find the small detachable board with the parallel-plate capacitor pattern (pattern 1 on page 18) on it, and screw it to the underside of the equipotential board, with the silver-painted side facing down toward the tabletop. Use the washers to protect the silver paint so that it doesn't get scraped off when you tighten the screws. Now connect the voltage source (using the provided wires) to the two large screws on either side of the board. Connect the multimeter so that you can mea-

sure the voltage difference across the terminals of the voltage source. Adjust the voltage source to give 8 volts.

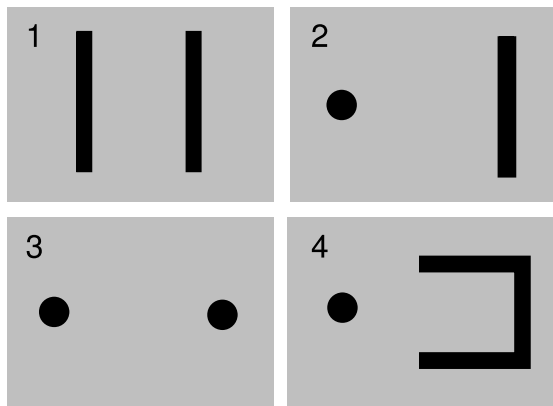

If you press down on the board, you can slip the paper between the board and the four buttons you see at the corners of the board. Tape the paper to your board, because the buttons aren't very dependable. There are plastic stencils in some of the envelopes, and you can use these to draw the electrodes accurately onto your paper so you know where they are. The photo, for example, shows pattern 3 traced onto the paper.

Now put the U-probe in place so that the top is above the equipotential board and the bottom of it is below the board. You will first be looking for places on the pattern board where the voltage is one volt — look for places where the meter reads 1.0 and mark them through the hole on the top of your U-probe with a pencil or pen. You should find a whole bunch of places there the voltage equals one volt, so that you can draw a nice constant-voltage curve connecting them. (If the line goes very far or curves strangely, you may have to do more.) You can then repeat the procedure for 2 V, 3 V, and so on. Label each constant-voltage curve. Once you've finished tracing the equipotentials, everyone in your group will need one copy of each of the two patterns you do, so you will need to photocopy them or simply trace them by hand.

If you're using the PRO-100 meters, they will try to outsmart you by automatically choosing a range. Most people find this annoying. To defeat this misfeature, press the RANGE button, and you'll see the AUTO indicator on the screen turn off.

Repeat this procedure with another pattern. Groups 1 and 4 should do patterns 1 and 2; groups 2 and 5 patterns 1 and 3; groups 3, 6, and 7 patterns 1 and 4.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** Looking at a plot of constant-voltage curves, how could you tell where the strongest electric fields would be? (Don't just say that the field is strongest when you're close to "the charge," because you may have a complex charge distribution, and we don't have any way to see or measure the charge distribution.)

**P2** What would the constant-voltage curves look like in a region of uniform electric field (i.e., one in which the  $\mathbf{E}$  vectors are all the same strength, and all in the same direction)?

## Self-Check

Calculate at least one numerical electric field value to make sure you understand how to do it.

You have probably found some constant-voltage curves that form closed loops. Do the electric field patterns ever seem to close back on themselves? Make sure you understand why or why not.

Make sure the people in your group all have a copy of each pattern.

## Analysis

On each plot, find the strongest and weakest electric fields, and calculate them.

On top of your plots, draw in electric field vectors. You will then have two different representations of the field superimposed on one another.

As always when drawing vectors, the lengths of the arrows should represent the magnitudes of the vectors, although you don't need to calculate them all numerically or use an actual scale. Remember that electric field vectors are always perpendicular to constant-voltage curves. The electric field lines point from high voltage to low voltage, just as the force on a rolling ball points downhill.



# 5 Relativity

## Apparatus

magnetic balance ..... 1/group  
meter stick ..... 1/group  
multimeter (BK or PRO-100, not HP) .... 1/group  
laser ..... 1/group  
vernier calipers ..... 1/group  
photocopy paper, for use as a weight  
DC power supply (Mastech, 30 A)  
box of special cables  
scissors

## Goal


Measure the speed of light.

## Introduction

Oersted discovered that magnetism is an interaction of moving charges with moving charges, but it wasn't until almost a hundred years later that Einstein showed why such an interaction must exist: magnetism occurs as a direct result of his theory of relativity. Since magnetism is a purely relativistic effect, and relativistic effects depend on the speed of light, any measurement of a magnetic effect can be used to determine the speed of light.

## Setup

The idea is to set up opposite currents in two wires, A and B, one under the other, and use the repulsion between the currents to create an upward force on the top wire, A. The top wire is on the arm of a balance, which has a stable equilibrium because of the weight C hanging below it. You initially set up the balance with no current through the wires, adjusting the counterweight D so that the distance between the wires is as small as possible. What we care about is really the center-to-center distance (which we'll call  $R$ ), so even if the wires are almost touching, there's still a millimeter or two worth of distance between them. By shining a laser at the mirror, E, and observing the spot it makes on the wall, you can very accurately determine this particular position of the balance, and tell later on when you've reproduced it. If you put a current through the wires, it will raise

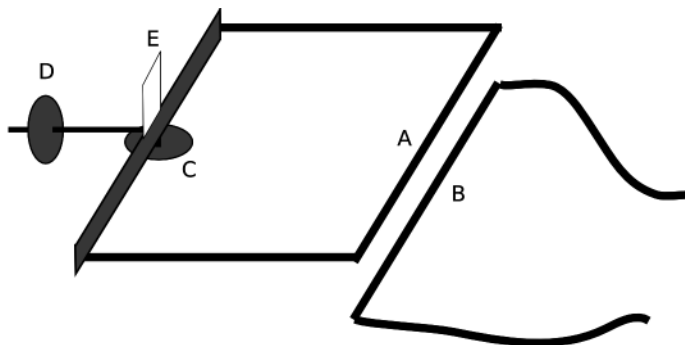

wire A. The torque made by the magnetic repulsion is now canceling the torque made by gravity directly on all the hardware, such as the masses C and D. This gravitational torque was zero before, but now you don't know what it is. The trick is to put a tiny weight on top of wire A, and adjust the current so that the balance returns to the position it originally had, as determined by the laser dot on the wall. You now know that the gravitational torque acting on the original apparatus (everything except for the weight) is back to zero, so the only torques acting are the torque of gravity on the staple and the magnetic torque. Since both these torques are applied at the same distance from the axis, the forces creating these torques must be equal as well. You can therefore infer the magnetic force that was acting.

For a weight, you can carefully and accurately cut a small rectangular piece out of a sheet of photocopy paper. In fall 2013, my students found that 500 sheets of SolCopy 20 lb paper were 2307.0 g. About 1/100 of a sheet seemed to be a good weight to use.

It's very important to get the wires A and B perfectly parallel. The result depends strongly on the small distance  $R$  between their centers, and if the wires aren't straight and parallel, you won't even have a well defined value of  $R$ .

The following technique allows  $R$  to be measured accurately. The idea is to compare the position of the laser spot on the wall when the balance is in its normal position, versus the position where the wires are touching. Using a small-angle approximation, you can then find the angle  $\theta_r$  by which the reflected beam moved. This is twice the angle  $\theta_m = \theta_r/2$  by which the mirror moved.<sup>1</sup> Once you

<sup>1</sup>To see this, imagine the following example that is unreal-

know the angle by which the moving arm of the apparatus moved, you can accurately find the air gap between the wires, and then add in twice the radius of the wires, which can be measured accurately with vernier calipers. For comparison, try to do as good a job as you can of measuring  $R$  directly by positioning the edges of the vernier calipers at the centers of the wires. If the two values of  $R$  don't agree, go back and figure out what went wrong; one possibility is that your wire is slightly bent and needs to be straightened.

You need to minimize the resistance of the apparatus, or else you won't be able to get enough current through it to cancel the weight of the staple. Most of the resistance is at the polished metal knife-edges that the moving part of the balance rests on. It may be necessary to clean the surfaces, or even to freshen them a little with a file to remove any layer of oxidation. Use the separate BK meter to measure the current — not the meter built into the power supply.


The power supply has some strange behavior that makes it not work unless you power it up in exactly the right way. It has four knobs, going from left to right: (1) current regulation, (2) over-voltage protection, (3) fine voltage control, (4) coarse voltage control. Before turning the power supply on, turn knobs 1 and 2 all the way up, and knobs 3 and 4 all the way down. Turn the power supply on. Now use knobs 3 and 4 to control how much current flows.

## Analysis


The first figure below shows a model that explains the repulsion felt by one of the charges in wire A due to all the charges in wire B. This is represented in the frame of the lab. For convenience of analysis, we give the model some unrealistic features: rather than having positively charged nuclei at rest and negatively charged electrons moving, we pretend that both are moving, in opposite directions. Since wire B has zero net density of charge everywhere, it creates no electric fields. (If you like, you can verify this during lab by putting tiny pieces of paper near the wires and verifying that they do not feel any static-electrical attraction.) Since there is no electric field, the force on the charge in wire A must be purely magnetic.

The second figure shows the same scene from the

istic but easy to figure out. Suppose that the incident beam is horizontal, and the mirror is initially vertical, so that the reflected beam is also horizontal. If the mirror is then tilted backward by 45 degrees, the reflected beam will be straight up,  $\theta_r = 90$  degrees.

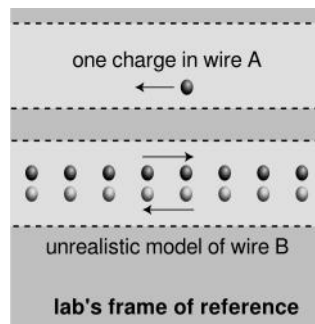

point of view of the charge in wire A. This charge considers itself to be at rest, and it also sees the light-colored charges in B as being at rest. In this frame the dark-colored charges in B are the only ones moving, and they move with twice the speed they had in the lab frame. In this frame, the particle in A is at rest, so it can't feel any magnetic force. The force is now considered to be purely electric. This electric force exists because the dark charges are relativistically contracted, which makes them more dense than their light-colored neighbors, causing a nonzero net density of charge in wire B.

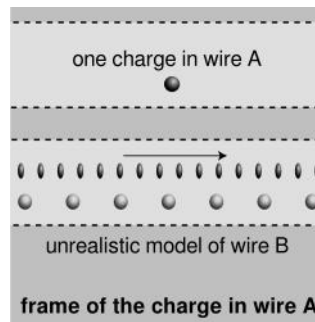

We've considered the force acting on a single charge in wire A. The actual force we observe in the experiment is the sum of all the forces acting on all such charges (of both signs). As in the slightly different example analyzed in section 23.3 of *Light and Matter*, this effect is proportional to the product of the speeds of the charges in the two wires, divided by  $c^2$ . Therefore the effect must be proportional to the product of the currents over  $c^2$ . In this experiment, the same current flows through wire A and then comes back through B in the opposite direction, so we conclude that the force must be proportional to  $I^2/c^2$ .

In the second frame, the force is purely electrical, and as shown in example 4 in section 22.3 of *Light and Matter*, the electric field of a charged wire falls off in proportion to  $1/R$ , where  $R$  is the distance from the wire. Electrical forces are also proportional

to the Coulomb constant  $k$ .

The longer the wires, the more charges interact, so we must also have a proportionality to the length  $\ell$ .

Putting all these factors together, we find that the force is proportional to  $kI^2\ell/c^2R$ . We can easily verify that the units of this expression are newtons, so the only possible missing factor is something unitless. This unitless factor turns out to be 2 — essentially the same 2 found in example 14 in section 22.7. The result for the repulsive force between the two wires is

$$F = \frac{k}{c^2} \cdot \frac{2I^2\ell}{R} \quad .$$

By solving this equation you can find  $c$ . Your final result is the speed of light, with error bars. Compare with the previously measured value of  $c$  and give a probabilistic interpretation, as in the examples in appendix 2.

In your writeup, give both the values of  $R$  (laser and eyeball). The laser technique is inherently better, so that's the value you should use in extracting  $c$ , but I want to see both values of  $R$  because some groups in the past have had a bigger discrepancy than I would have expected. If you have a large discrepancy, get my attention during lab and we can see whether it might be due to a bent wire, or some other cause.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

Do the laser safety checklist, Appendix 7, tear it out, and turn it in at the beginning of lab. If you don't understand something, don't initial that point, and ask your instructor for clarification before you start the lab.

**P1** Show that the equation for the force between the wires has units of newtons.

**P2** Do the algebra to solve for  $c$  in terms of the measured quantities.



# 6 Magnetism

## Apparatus

bar magnet (stack of 6 Nd)  
compass  
Hall effect magnetic field probes  
LabPro interfaces, DC power supplies, and USB cables  
2-meter stick  
Heath solenoids .....2/group  
Mastech power supply .....1/group  
wood blocks .....2/group  
PRO-100 multimeter (in lab bench .....1/group  
another multimeter .....1/group  
D-cell batteries and holders  
Cenco decade resistor box .....1/group

## Goal

Find how the magnetic field of a magnet changes with distance along one of the magnet's lines of symmetry.

## Introduction

### A Variation of Field With Distance: Deflection of a Magnetic Compass

You can infer the strength of the bar magnet's field at a given point by putting the compass there and seeing how much it is deflected from north.

The task can be simplified quite a bit if you restrict yourself to measuring the magnetic field at points along one of the magnet's two lines of symmetry, shown in the top figure on the page three pages after this one.

If the magnet is flipped across the vertical axis, the north and south poles remain just where they were, and the field is unchanged. That means the entire magnetic field is also unchanged, and the field at a point such as point b, along the line of symmetry, must therefore point straight up.

If the magnet is flipped across the horizontal axis, then the north and south poles are swapped, and the field everywhere has to reverse its direction. Thus, the field at points along this axis, e.g., point a, must point straight up or down.

Line up your magnet so it is pointing east-west.

Choose one of the two symmetry axes of your magnet, and measure the deflection of the compass at two points along that axis, as shown in the figure at the end of the lab. As part of your prelab, you will use vector addition to find an equation for  $B_m/B_e$ , the magnet's field in units of the Earth's, in terms of the deflection angle  $\theta$ . For your first point, find the distance  $r$  at which the deflection is 70 degrees; this angle is chosen because it's about as big as it can be without giving very poor relative precision in the determination of the magnetic field. For your second data-point, use twice that distance. By what factor does the field decrease when you double  $r$ ?

The lab benches contain iron or steel parts that distort the magnetic field. You can easily observe this simply by putting a compass on the top of the bench and sliding it around to different places. To work around this problem, lay a 2-meter stick across the space between two lab benches, and carry out the experiment along the line formed by the stick. Even in the air between the lab benches, the magnetic field due to the building materials in the building is significant, and this field varies from place to place. Therefore you should move the magnet while keeping the compass in one place. Then the field from the building becomes a fixed part of the background experienced by the compass, just like the earth's field.

Note that the measurements are very sensitive to the relative position and orientation of the bar magnet and compass.


Based on your two data-points, form a hypothesis about the variation of the magnet's field with distance according to a power law  $B \propto r^p$ .

### B Variation of Field With Distance: Hall Effect Magnetometer


In this part of the lab, you will test your hypothesis about the power law relationship  $B \propto r^p$ ; you will find out whether the field really does obey such a law, and if it does, you will determine  $p$  accurately.

This part of the lab uses a device called a Hall effect magnetometer for measuring magnetic fields. It works by sending an electric current through a substance, and measuring the force exerted on those moving charges by the surrounding magnetic field. The probe only measures the component of the magnetic field vector that is parallel to its own axis. Plug the probe into CH 1 of the LabPro interface, connect




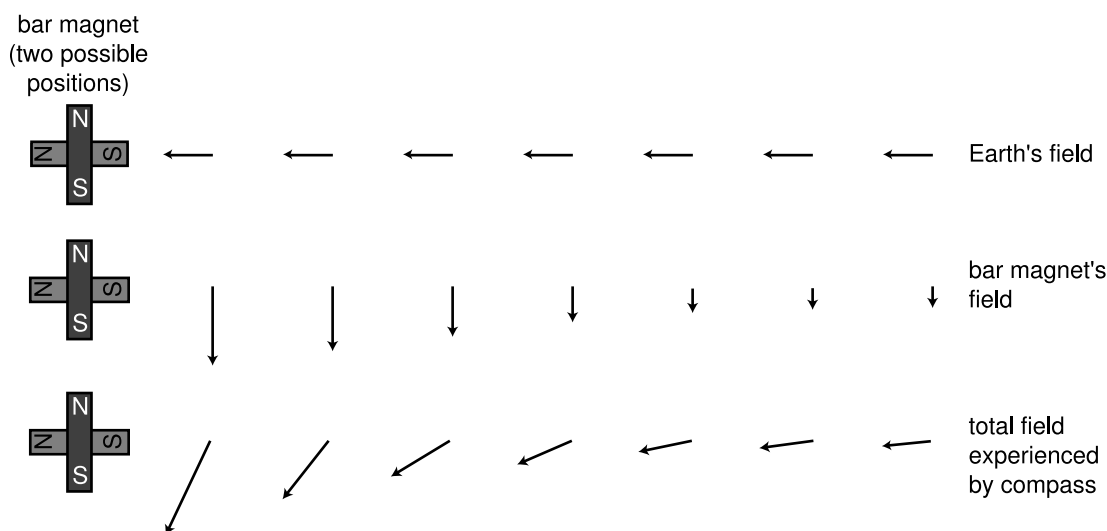

Part A, measuring the variation of the bar magnet's field with respect to distance.

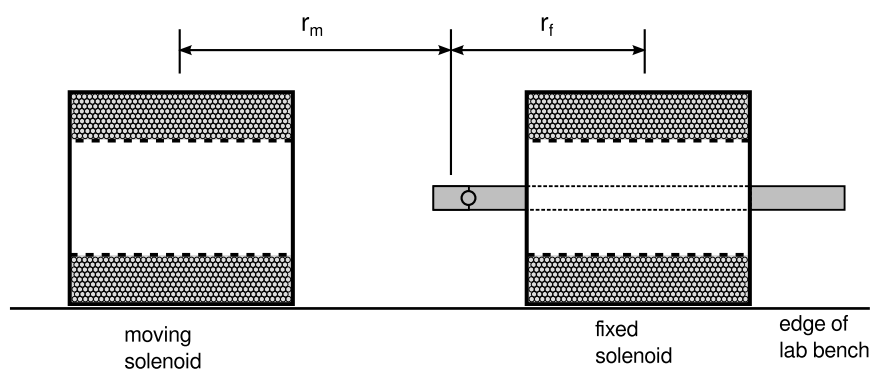

Part B, a different method of measuring the variation of field with distance. The solenoids are shown in cross-section, with empty space on their interiors and their axes running right-left.

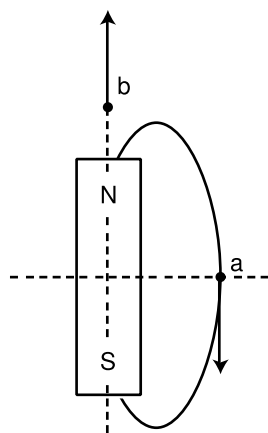

the interface to the computer's USB port, and plug the interface's DC power supply in to it. Start up version 3 of Logger Pro, and it will automatically recognize the probe and start displaying magnetic fields on the screen, in units of mT (millitesla). The

probe has two ranges, one that can read fields up to 0.3 mT, and one that goes up to 6.4 mT. Select the more sensitive 0.3 mT scale using the switch on the probe.

The technique is shown in the bottom figure on the last page of the lab. Identical solenoids (cylindrical coils of wire) are positioned with their axes coinciding, by lining up their edges with the edge of the lab bench. When an electrical current passes through a coil, it creates a magnetic field. At distances that are large compared to the size of the solenoid, we expect that this field will have the same universal pattern as with any magnetic dipole. The sensor is positioned on the axis, with wood blocks (not shown) to hold it up. One solenoid is fixed, while the other is moved to different positions along the axis, including positions (more distant than the one shown) at which we expect its contribution to the field at the sensor to be of the universal dipole form.

The key to the high precision of the measurement is that in this configuration, the fields of the two solenoids can be made to cancel at the position of the probe. Because of the solenoids' unequal distances from the probe, this requires unequal currents. Because the fields cancel, the probe can be used on its most sensitive and accurate scale; it can also be zeroed when the circuits are open, so that the effect of any ambient field is removed. For example, suppose that at a certain distance  $r_m$ , the current  $I_m$  through the moving coil has to be five times greater than the current  $I_f$  through the fixed coil at the constant distance  $r_f$ . Then we have determined that the field pattern of these coils is such that increasing the distance along the axis from  $r_f$  to  $r_m$  causes the field to fall off by a factor of five.

It's a good idea to take data all the way down to  $r_m = 0$ , since this makes it possible to see on a graph where the field does and doesn't behave like a dipole. Note that the distances  $r_f$  and  $r_m$  can't be measured directly with good precision.

The Mastech power supply is capable of delivering a large amount of current, so it can be used to provide  $I_m$ , which needs to be high when  $r_m$  is large. The power supply has some strange behavior that makes it not work unless you power it up in exactly the right way. It has four knobs, going from left to right: (1) current regulation, (2) over-voltage protection, (3) fine voltage control, (4) coarse voltage control. Before turning the power supply on, turn knobs 1 and 2 all the way up, and knobs 3 and 4 all the way down. Turn the power supply on. Now use knobs 3 and 4 to control how much current flows.

At large values of  $r_m$ , it can be difficult to get a power supply to give a *small* enough  $I_f$ . Try using a battery, and further reducing the current by placing another resistance in series with the coil. The Cenco decade resistance boxes can be used for this purpose; they are variable resistors whose resistance can be dialed up as desired using decimal knobs. Use the plugs on the resistance box labeled H and L.

For every current measurement, make sure to use the most sensitive possible scale on the meter to get as many sig figs as possible. This is why the ammeter built into the Mastech power supply is not useful here. I found it to be a hassle to measure  $I_m$  with an ammeter, because the currents required were often quite large, and I kept inadvertently blowing the fuse on the milliamp scale. For this reason, you may actually want to measure  $V_m$ , the voltage difference across the moving solenoid. Conceptually, magnetic fields are caused by moving charges, cur-

rent is a measure of moving charge, and therefore current is what is relevant here. But if the DC resistance of the coil is fixed, the current and voltage are proportional to one another, assuming that the voltage is measured directly across the coil and the resistance of the banana-plug connections is either negligible or constant.

As shown in a lecture demonstration, deactivating the electromagnet requires getting rid of the energy stored in the magnetic field, and this can be done in more than one way. If you use your hand to break the circuit by pulling out a banana plug, the energy is dissipated in a spark, and a large value of  $I_m$  is being used the result can be an unpleasant shock. To avoid this, deactivate the moving coil by turning down the knob on the power supply rather than by breaking the circuit.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** In part A, suppose that when the compass is 11.0 cm from the magnet, it is 45 degrees away from north. What is the strength of the bar magnet's field at this location in space, in units of the Earth's field?

**P2** Find  $B_m/B_e$  in terms of the deflection angle  $\theta$  measured in part A. As a special case, you should be able to recover your answer to P1.

## Analysis

Determine the variation of the solenoid's magnetic field with distance. Look for a power-law relationship using the log-log graphing technique described in appendix 5. Does the power law hold for all the distances you investigated, or only at large distances? No error analysis is required.



# 7 Electromagnetism

## Apparatus

oscilloscope (Tektronix TDS 1001B) ..... 1/group  
microphone (RS 33-1067) ..... for 6 groups  
microphone (Shure C606) ..... for 1 group  
various tuning forks, mounted on wooden boxes  
solenoid (Heath) ..... 1/group  
2-meter wire with banana plugs ..... 1/group  
magnet (stack of 6 Nd) ..... 1/group  
masking tape  
string

## Goals

- Learn to use an oscilloscope.
- Observe electric fields induced by changing magnetic fields.
- Build a generator.
- Discover Lenz's law.


## Introduction

Physicists hate complication, and when physicist Michael Faraday was first learning physics in the early 19th century, an embarrassingly complex aspect of the science was the multiplicity of types of forces. Friction, normal forces, gravity, electric forces, magnetic forces, surface tension — the list went on and on. Today, 200 years later, ask a physicist to enumerate the fundamental forces of nature and the most likely response will be “four: gravity, electromagnetism, the strong nuclear force and the weak nuclear force.” Part of the simplification came from the study of matter at the atomic level, which showed that apparently unrelated forces such as friction, normal forces, and surface tension were all manifestations of electrical forces among atoms. The other big simplification came from Faraday's experimental work showing that electric and magnetic forces were intimately related in previously unexpected ways, so intimately related in fact that we now refer to the two sets of force-phenomena under a single term, “electromagnetism.”

Even before Faraday, Oersted had shown that there was at least some relationship between electric and

magnetic forces. An electrical current creates a magnetic field, and magnetic fields exert forces on an electrical current. In other words, electric forces are forces of charges acting on charges, and magnetic forces are forces of moving charges on moving charges. (Even the magnetic field of a bar magnet is due to currents, the currents created by the orbiting electrons in its atoms.)

Faraday took Oersted's work a step further, and showed that the relationship between electricity and magnetism was even deeper. He showed that a changing electric field produces a magnetic field, and a changing magnetic field produces an electric field. Faraday's work forms the basis for such technologies as the transformer, the electric guitar, the transformer, and generator, and the electric motor. It also led to the understanding of light as an electromagnetic wave.

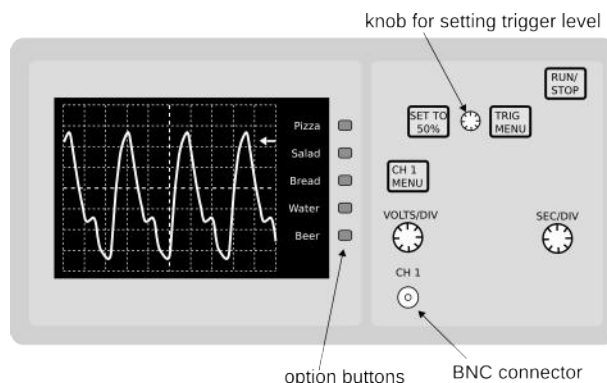


### *The Oscilloscope*


An oscilloscope graphs an electrical signal that varies as a function of time. The graph is drawn from left to right across the screen, being painted in real time as the input signal varies. The purpose of the oscilloscope in this lab is to measure electromagnetic induction, but to get familiar with the oscilloscope, we'll start out by using the signal from a microphone as an input, allowing you to see sound waves.

The input signal is supplied in the form of a voltage. The input connector on the front of the oscilloscope accepts a type of cable known as a BNC cable. A BNC cable is a specific example of coaxial cable (“coax”), which is also used in cable TV, radio, and computer networks. The electric current flows

in one direction through the central conductor, and returns in the opposite direction through the outside conductor, completing the circuit. The outside conductor is normally kept at ground, and also serves as shielding against radio interference. The advantage of coaxial cable is that it is capable of transmitting rapidly varying signals without distortion.

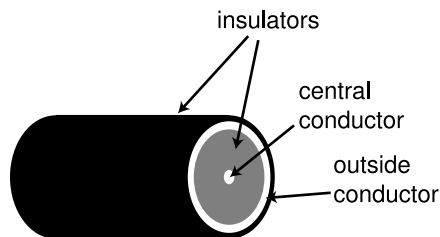

Most of the voltages we wish to measure are not big enough to use directly for the vertical deflection voltage, so the oscilloscope actually amplifies the input voltage, i.e., the small input voltage is used to control a much larger voltage generated internally. The amount of amplification is controlled with a knob on the front of the scope. For instance, setting the knob on 1 mV selects an amplification such that 1 mV at the input deflects the electron beam by one square of the 1-cm grid. Each 1-cm division is referred to as a “division.”

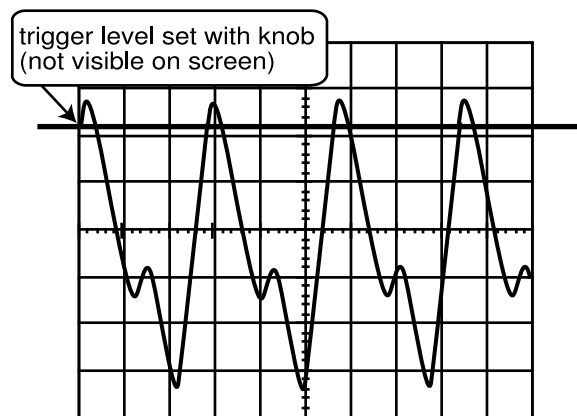

### *The Time Base and Triggering*

Since the  $X$  axis represents time, there also has to be a way to control the time scale, i.e., how fast the imaginary “penpoint” sweeps across the screen. For instance, setting the knob on 10 ms causes it to sweep across one square in 10 ms. This is known as the time base.

In the figure, suppose the time base is 10 ms. The scope has 10 divisions, so the total time required for the beam to sweep from left to right would be 100

ms. This is far too short a time to allow the user to examine the graph. The oscilloscope has a built-in method of overcoming this problem, which works well for periodic (repeating) signals. The amount of time required for a periodic signal to perform its pattern once is called the period. With a periodic signal, all you really care about seeing what one period or a few periods in a row look like — once you’ve seen one, you’ve seen them all. The scope displays one screenful of the signal, and then keeps on overlaying more and more copies of the wave on top of the original one. Each trace is erased when the next one starts, but is being overwritten continually by later, identical copies of the wave form. You simply see one persistent trace.

How does the scope know when to start a new trace? If the time for one sweep across the screen just happened to be exactly equal to, say, four periods of the signal, there would be no problem. But this is unlikely to happen in real life — normally the second trace would start from a different point in the waveform, producing an offset copy of the wave. Thousands of traces per second would be superimposed on the screen, each shifted horizontally by a different amount, and you would only see a blurry band of light.

To make sure that each trace starts from the same point in the waveform, the scope has a triggering circuit. You use a knob to set a certain voltage level, the trigger level, at which you want to start each trace. The scope waits for the input to move across the trigger level, and then begins a trace. Once that trace is complete, it pauses until the input crosses the trigger level again. To make extra sure that it is really starting over again from the same point in the waveform, you can also specify whether you want to start on an increasing voltage or a decreasing voltage — otherwise there would always be at least two points in a period where the voltage crossed your trigger level.

## **Setup**

To start with, we’ll use a sine wave generator, which makes a voltage that varies sinusoidally with time. This gives you a convenient signal to work with while you get the scope working.

The figure on the preceding page is a simplified drawing of the front panel of a digital oscilloscope, showing only the most important controls you’ll need for this lab. When you turn on the oscilloscope, it will take a while to start up.

Press DEFAULT SETUP.

Use the SEC/DIV knob to put the time base on something reasonable compared to the period of the signal you're looking at. The time base is displayed on the screen, e.g., 10 ms/div, or 1 s/div.

Use the VOLTS/DIV knob to put the voltage scale (Y axis) on a reasonable scale compared to the amplitude of the signal you're looking at.

The scope has two channels, i.e., it can accept input through two BNC connectors and display both or either. You'll only be using channel 1, which is the only one represented in the simplified drawing. By default, the oscilloscope draws graphs of both channels' inputs; to get rid of ch. 2, hold down the CH 2 MENU button (not shown in the diagram) for a couple of seconds. You also want to make sure that the scope is triggering on CH 1, rather than CH 2. To do that, press the TRIG MENU button, and use an option button to select CH 1 as the source. Set the triggering mode to normal, which is the mode in which the triggering works as I've described above. If the trigger level is set to a level that the signal never actually reaches, you can play with the knob that sets the trigger level until you get something. A quick and easy way to do this without trial and error is to use the SET TO 50% button, which automatically sets the trigger level to midway between the top and bottom peaks of the signal.

You want to select AC, not DC or GND, on the channel you're using. You are looking at a voltage that is alternating, creating an alternating current, "AC." The "DC" setting is only necessary when dealing with constant or very slowly varying voltages. The "GND" simply draws a graph using  $y = 0$ , which is only useful in certain situations, such as when you can't find the trace. To select AC, press the CH 1 MENU button, and select AC coupling.

Observe the effect of changing the voltage scale and time base on the scope. Try changing the frequency and amplitude on the sine wave generator.

You can freeze the display by pressing RUN/STOP, and then unfreeze it by pressing the button again.

## Preliminary Observations

Now try observing signals from the microphone.

Once you have your setup working, try measuring the period and frequency of the sound from a tuning fork, and make sure your result for the frequency is the same as what's written on the tuning fork.

## Qualitative Observations

In this lab you will use a permanent magnet to produce changing magnetic fields. This causes an electric field to be induced, which you will detect using a solenoid (spool of wire) connected to an oscilloscope. The electric field drives electrons around the solenoid, producing a current which is detected by the oscilloscope.

Note that although I've described the standard way of triggering a scope, when the time base is very long, triggering becomes unnecessary. These scopes are programmed so that when the time base is very long, they simply continuously display traces.

### A A constant magnetic field

Do you detect any signal on the oscilloscope when the magnet is simply placed at rest inside the solenoid? Try the most sensitive voltage scale.

### B A changing magnetic field

Do you detect any signal when you move the magnet or wiggle it inside the solenoid or near it? What happens if you change the speed at which you move the magnet?

### C Moving the solenoid

What happens if you hold the magnet still and move the solenoid?

The poles of the magnet are its flat faces. In later parts of the lab you will need to know which is north. Determine this now by hanging it from a string and seeing how it aligns itself with the Earth's field. The pole that points north is called the north pole of the magnet. The field pattern funnels into the body of the magnet through its south pole, and reemerges at its north pole.

### D A generator

Tape the magnet securely to the eraser end of a pencil so that its flat face (one of its two poles) is like the head of a hammer, and mark the north and south poles of the magnet for later reference. Spin the pencil near the solenoid and observe the induced signal.

You have built a generator. (I have unfortunately not had any luck lighting a lightbulb with the setup, due to the relatively high internal resistance of the solenoid.)

## Trying Out Your Understanding

### E Changing the speed of the generator

If you change the speed at which you spin the pencil, you will of course cause the induced signal to have a longer or shorter period. Does it also have any effect on the *amplitude* of the wave?

### F A solenoid with fewer loops

Use the two-meter cable to make a second solenoid with the same diameter but fewer loops. Compare the strength of the induced signals.

### G Dependence on distance

How does the signal picked up by your generator change with distance?

Try to explain what you have observed, and discuss your interpretations with your instructor.

## Lenz's Law

Lenz's law describes how the clockwise or counterclockwise direction of the induced electric field's whirlpool pattern relates to the changing magnetic field. The main result of this lab is a determination of how Lenz's law works. To focus your reasoning, here are four possible forms for Lenz's law:

1. The electric field forms a pattern that is clockwise when viewed along the direction of the  $B$  vector of the changing magnetic field.
2. The electric field forms a pattern that is counterclockwise when viewed along the direction of the  $B$  vector of the changing magnetic field.
3. The electric field forms a pattern that is clockwise when viewed along the direction of the  $\Delta B$  vector of the changing magnetic field.
4. The electric field forms a pattern that is counterclockwise when viewed along the direction of the  $\Delta B$  vector of the changing magnetic field.


Your job is to figure out which is correct.

The most direct way to figure out Lenz's law is to make a tomahawk-chopping motion that ends up with the magnet in the solenoid, observing whether the pulse induced is positive or negative. What hap-

pens when you reverse the chopping motion, or when you reverse the north and south poles of the magnet? Try all four possible combinations and record your results.

To set up the scope, press DEFAULT SETUP. This should have the effect of setting the scope on DC coupling, which is what you want. (If it's on AC coupling, it tries to filter out any DC part of the input signals, which distorts the results.) To check that you're on DC coupling, you can do CH 1 MENU, and check that Coupling says DC. Set the triggering mode ("Mode") to Auto.

Make sure the scope is on DC coupling, not AC coupling, or your pulses will be distorted.

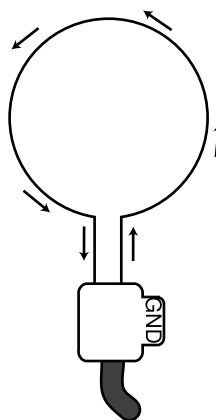

It can be tricky to make the connection between the polarity of the signal on the screen of the oscilloscope and the direction of the electric field pattern. The figure shows an example of how to interpret a positive pulse: the current must have flowed through the scope from the center conductor of the coax cable to its outer conductor (marked GND on the coax-to-banana converter).

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** In the sample oscilloscope trace shown on page 29, what is the period of the waveform? What is its frequency? The time base is 10 ms.

**P2** In the same example, again assume the time base is 10 ms/division. The voltage scale is 2 mV/division. Assume the zero voltage level is at the middle

of the vertical scale. (The whole graph can actually be shifted up and down using a knob called “position.”) What is the trigger level currently set to? If the trigger level was changed to 2 mV, what would happen to the trace?

**P3** Referring to the chapter of your textbook on sound, which of the following would be a reasonable time base to use for an audio-frequency signal? 10 ns,  $1\mu$  s, 1 ms, 1 s

**P4** Does the oscilloscope show you the signal’s period, or its wavelength? Explain. [Skip this question if you’re in Physics 222.]

**P5** The time-scale for all the signals is determined by the fact that you’re wiggling and waving the magnet by hand, so what’s a reasonable order of magnitude to choose for the time base on the oscilloscope? [Skip this question if you’re in Physics 222.]

## Self-Check

Determine which version of Lenz’s law is correct.





# 8 The Charge to Mass Ratio of the Electron

## Apparatus

vacuum tube with Helmholtz	
coils (Leybold )	1
Cenco 33034 HV supply	1
12-V DC power supplies (Thornton)	1
multimeters (Fluke or HP)	2
compass	1
ruler	1
banana-plug cables	

## Goal


Measure the charge-to-mass ratio of the electron.

## Introduction

Why should you believe electrons exist? By the turn of the twentieth century, not all scientists believed in the literal reality of atoms, and few could imagine smaller objects from which the atoms themselves were constructed. Over two thousand years had elapsed since the Greeks first speculated that atoms existed based on philosophical arguments without experimental evidence. During the Middle Ages in Europe, “atomism” had been considered highly suspect, and possibly heretical. Finally by the Victorian era, enough evidence had accumulated from chemical experiments to make a persuasive case for atoms, but subatomic particles were not even discussed.

If it had taken two millennia to settle the question of atoms, it is remarkable that another, subatomic level of structure was brought to light over a period of only about five years, from 1895 to 1900. Most of the crucial work was carried out in a series of experiments by J.J. Thomson, who is therefore often considered the discoverer of the electron.

In this lab, you will carry out a variation on a crucial experiment by Thomson, in which he measured the ratio of the charge of the electron to its mass,  $q/m$ . The basic idea is to observe a beam of electrons in a region of space where there is an approximately uniform magnetic field,  $B$ . The electrons are emitted perpendicular to the field, and, it turns out, travel in a circle in a plane perpendicular to it. The force

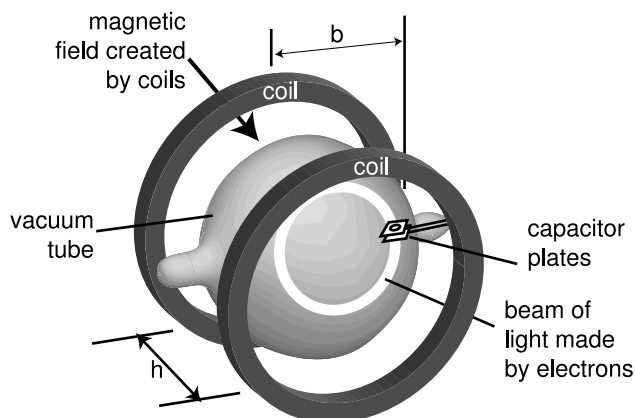

of the magnetic field on the electrons is

$$F = qvB \quad , \quad (1)$$

directed towards the center of the circle. Their acceleration is

$$a = \frac{v^2}{r} \quad , \quad (2)$$

so using  $F = ma$ , we can write

$$qvB = \frac{mv^2}{r} \quad . \quad (3)$$

If the initial velocity of the electrons is provided by accelerating them through a voltage difference  $V$ , they have a kinetic energy equal to  $qV$ , so

$$\frac{1}{2}mv^2 = qV \quad . \quad (4)$$

From equations 3 and 4, you can determine  $q/m$ . Note that since the force of a magnetic field on a moving charged particle is always perpendicular to the direction of the particle's motion, the magnetic field can never do any work on it, and the particle's KE and speed are therefore constant.

You will be able to see where the electrons are going, because the vacuum tube is filled with a hydrogen gas at a low pressure. Most electrons travel large distances through the gas without ever colliding with a hydrogen atom, but a few do collide, and the atoms then give off blue light, which you can see. Although I will loosely refer to “seeing the beam,” you are really seeing the light from the collisions, not the beam of electrons itself. The manufacturer of the

tube has put in just enough gas to make the beam visible; more gas would make a brighter beam, but would cause it to spread out and become too broad to measure it precisely.

The field is supplied by an electromagnet consisting of two circular coils, each with 130 turns of wire (the same on all the tubes we have). The coils are placed on the same axis, with the vacuum tube at the center. A pair of coils arranged in this type of geometry are called Helmholtz coils. Such a setup provides a nearly uniform field in a large volume of space between the coils, and that space is more accessible than the inside of a solenoid.

### *Safety*

You will use the Cenco high-voltage supply to make a DC voltage of about 300 V. Two things automatically keep this from being very dangerous:

Several hundred DC volts are far less dangerous than a similar AC voltage. The household AC voltages of 110 and 220 V are more dangerous because AC is more readily conducted by body tissues.

The HV supply will blow a fuse if too much current flows.

Do the high voltage safety checklist, Appendix 6, tear it out, and turn it in at the beginning of lab. If you don't understand something, don't initial that point, and ask your instructor for clarification before you start the lab.

## Setup

Before beginning, make sure you do not have any computer disks near the apparatus, because the magnetic field could erase them.

Heater circuit: As with all vacuum tubes, the cathode is heated to make it release electrons more easily. There is a separate low-voltage power supply built into the high-voltage supply. It has a set of green plugs that, in different combinations, allow you to get various low voltage values. Use it to supply 6 V to the terminals marked "heater" on the vacuum tube. The tube should start to glow.

Electromagnet circuit: Connect the other Thornton power supply, in series with an ammeter, to the terminals marked "coil." The current from this power supply goes through both coils to make the magnetic

field. Verify that the magnet is working by using it to deflect a nearby compass.

High-voltage circuit: Leave the Cenco HV supply unplugged. It is really three HV circuits in one box. You'll be using the circuit that goes up to 500 V. Connect it to the terminals marked "anode." Ask your instructor to check your circuit. Now plug in the HV supply and turn up the voltage to 300 V. You should see the electron beam. If you don't see anything, try it with the lights dimmed.

## Observations

Make the necessary observations in order to find  $q/m$ , carrying out your plan to deal with the effects of the Earth's field. The high voltage is supposed to be 300 V, but to get an accurate measurement of what it really is you'll need to use a multimeter rather than the poorly calibrated meter on the front of the high voltage supply.

The beam can be measured accurately by using the glass rod inside the tube, which has a centimeter scale marked on it.

Be sure to compute  $q/m$  before you leave the lab. That way you'll know you didn't forget to measure something important, and that your result is reasonable compared to the currently accepted value.

There is a glass rod inside the vacuum tube with a centimeter scale on it, so you can measure the diameter  $d$  of the beam circle simply by looking at the place where the glowing beam hits the scale. This is much more accurate than holding a ruler up to the tube, because it eliminates the parallax error that would be caused by viewing the beam and the ruler along a line that wasn't perpendicular to the plane of the beam. However, the manufacturing process used in making these tubes (they're probably hand-blown by a glass blower) isn't very precise, and on many of the tubes you can easily tell by comparison with the a ruler that, e.g., the 10.0 cm point on the glass rod is not really 10.0 cm away from the hole from which the beam emerges. Past students have painstakingly determined the appropriate corrections,  $a$ , to add to the observed diameters by the following electrical method. If you look at your answer to prelab question P1, you'll see that the product  $Br$  is always a fixed quantity in this experiment. It therefore follows that  $Id$  is also supposed to be constant. They measured  $I$  and  $d$  at two different values of  $I$ , and determined the correction  $a$  that had to be added to their  $d$  values in order to make the two values of  $Id$  equal. The results are as follows:

<i>serial number</i>	<i>a</i> (cm)
98-16	0.0
9849	0.0
99-08	+0.15
99-10	-0.2
99-17	+0.2
99-56	+0.3
031427	-0.3

If your apparatus is one that hasn't already had its  $a$  determined, then you should do the necessary measurements to calibrate it.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

The week before you are to do the lab, briefly familiarize yourself visually with the apparatus.

Do the high voltage safety checklist, Appendix 6, tear it out, and turn it in at the beginning of lab. If you don't understand something, don't initial that point, and ask your instructor for clarification before you start the lab.

**P1** Derive an equation for  $q/m$  in terms of  $V$ ,  $r$  and  $B$ .

**P2** For an electromagnet consisting of a single circular loop of wire of radius  $b$ , the field at a point on its axis, at a distance  $z$  from the plane of the loop, is given by

$$B = \frac{2\pi k I b^2}{c^2(b^2 + z^2)^{3/2}} \quad .$$

Starting from this equation, derive an equation for the magnetic field at the center of a pair of Helmholtz coils. Let the number of turns in each coil be  $N$  (in our case,  $N = 130$ ), let their radius be  $b$ , and let the distance between them be  $h$ . (In the actual experiment, the electrons are never exactly on the axis of the Helmholtz coils. In practice, the equation you will derive is sufficiently accurate as an approximation to the actual field experienced by the electrons.) If you have trouble with this derivation, see your instructor in his/her office hours.

**P3** Find the currently accepted value of  $q/m$  for the electron.

**P4** The electrons will be affected by the Earth's magnetic field, as well as the (larger) field of the

coils. Devise a plan to eliminate, correct for, or at least estimate the effect of the Earth's magnetic field on your final  $q/m$  value.

**P5** Of the three circuits involved in this experiment, which ones need to be hooked up with the right polarity, and for which ones is the polarity irrelevant?

**P6** What would you infer if you found the beam of electrons formed a helix rather than a circle?

## Analysis

Determine  $q/m$ , with error bars.

Answer the following questions:

**Q1.** Thomson started to become convinced during his experiments that the "cathode rays" observed coming from the cathodes of vacuum tubes were building blocks of atoms — what we now call electrons. He then carried out observations with cathodes made of a variety of metals, and found that  $q/m$  was the same in every case. How would that observation serve to test his hypothesis?

**Q2.** Why is it not possible to determine  $q$  and  $m$  themselves, rather than just their ratio, by observing electrons' motion in electric or magnetic fields?

**Q3.** Thomson found that the  $q/m$  of an electron was thousands of times larger than that of ions in electrolysis. Would this imply that the electrons had more charge? Less mass? Would there be no way to tell? Explain.



# 9 Radioactivity

Note to the lab technician: The isotope generator kits came with 250 mL bottles of eluting solution (0.9% NaCl in 0.04M HCl, made with deionized water). If we ever run out of the solution, we can make more from materials in the chem stockroom. The GM counters have 9 V batteries, which should be checked before lab.

## Apparatus

isotope generator kit  
Geiger-Müller (GM) counter  
computer with Logger Pro software and LabPro interface  
“grabber” clamp and stand  
wood blocks, 25 mm thick  
pieces of steel, 17 mm thick

## Goal

Determine the properties of an unknown radioactive source.

## Introduction

You’re a science major, but even if you weren’t, it would be important for you as a citizen and a voter to understand the properties of radiation. As an example of an important social issue, many environmentalists who had previously opposed nuclear power now believe that its benefits, due to reduction of global warming, outweigh its problems, such as disposal of waste. To understand such issues, you need to learn to reason about radioactivity quantitatively.

A radioactive substance contains atoms whose nuclei spontaneously decay into nuclei of a different type. Nobody has ever succeeded in finding a physical law that would predict when a particular nucleus will “choose” to decay. The process is random, but we can make quantitative statements about how quickly the process tends to happen. A radioactive substance has a certain *half-life*, defined as the time required before (on the average) 50% of its nuclei will have decayed.

## Safety

The radioactive source used in this lab is very weak. It is so weak that it is exempt from government regulation, it can be sent in the mail, and when the college buys new equipment, it is legal to throw out the old sources in the trash. The following table compares some radiation doses, including an estimate of the typical dose you might receive in this lab. These are in units of microSieverts ( $\mu\text{Sv}$ ).

CT scan	$\sim 10,000 \mu\text{Sv}$
natural background per year	2,000-7,000 $\mu\text{Sv}$
health guidelines for exposure to a fetus	1,000 $\mu\text{Sv}$
flying from New York to Tokyo	150 $\mu\text{Sv}$
this lab	$\sim 80 \mu\text{Sv}$
chest x-ray	50 $\mu\text{Sv}$

A variety of experiments seem to show cases in which low levels of radiation activate cellular damage control mechanisms, increasing the health of the organism. For example, there is evidence that exposure to radiation up to a certain level makes mice grow faster; makes guinea pigs’ immune systems function better against diphtheria; increases fertility in female humans, trout, and mice; improves fetal mice’s resistance to disease; reduces genetic abnormalities in humans; increases the life-spans of flour beetles and mice; and reduces mortality from cancer in mice and humans. This type of effect is called radiation hormesis. Nobody knows for sure, but it’s possible that you will receive a very tiny improvement in your health from the radiation exposure you experience in lab today.

Although low doses of radiation may be beneficial, governments, employers, and schools generally practice a philosophy called ALARA, which means to make radiation doses As Low As Reasonably Achievable. You should adhere to this approach in this lab. In general, internal exposure to radiation produces more of an effect than external exposure, so you should not eat or drink during this lab, and you should avoid getting any of the radioactive substances in an open cut. You should also reduce your exposure by not spending an unnecessarily large amount of time with your body very close to the source, e.g., you should not hold it in the palm of your hand for the entire lab period.

## The source and the GM counter

You are supplied with a radioactive source packaged inside a small plastic disk about the size of the spindle that fits inside a roll of scotch tape.

Our radiation detector for this lab is called a Geiger-Müller (GM) counter. It is a cylinder full of gas, with the outside of the cylinder at a certain voltage and a wire running down its axis at another voltage. The voltage difference creates a strong electric field. When ionizing radiation enters the cylinder, it can ionize the gas, separating negatively charged ions (electrons) and positively charged ones (atoms lacking some electrons). The electric field accelerates the ions, making them hit other air molecules, and causing a cascade of ions strong enough to be measured as an electric current.

If you look at the top side of the GM counter (behind the top of the front panel), you'll see a small window. Non-penetrating radiation can only get in through this thin layer of mica. (Gammas can go right in through the plastic housing.)

Put the bottom switch on Audio. The top switch doesn't have any effect on the data collection we'll be doing with the computer.

## Poisson Statistics

Although we will not be formally estimating error bars in this lab, the following information will be helpful in understanding what's going on when you're taking data. When we have a large number of things that may happen with some small probability, the total number of them that do happen is called a Poisson random variable (accent on the second syllable). For example, the number of houses burglarized in Fullerton this year is a Poisson random variable. When you count the number of nuclear decays in a certain time interval, the result is Poisson. The helpful thing to know is that when a Poisson variable has an average value  $N$ , its statistical uncertainty is  $\sqrt{N}$ . So for example if your GM counter counts 100 clicks in one minute, this is  $100 \pm 10$ . Knowing this will help you to have some idea whether, for example, an apparent change in the count rate is actually too small to be statistically meaningful, or whether you need to count for longer in order to get reliable results.

## Observations

### A Background

Use the GM counter to observe the background radiation in the room. This radiation is probably a combination of gamma rays from naturally occurring minerals in the ground plus betas and gammas from building materials such as concrete. If you like, you can walk around the room and see if you can detect any variations in the intensity of the background.

Estimate the rate at which the GM detector counts when it is exposed only to background. Starting at this point, it is more convenient to interface the GM counter to the computer. Plug the cable into DIG/SONIC 1 on the LabPro interface. Start LoggerPro 3 on the computer, and open the file Probes and Sensors : Radiation monitor : Counts versus time. The interface is not able to automatically identify this particular sensor, so the software will ask you to confirm that you really do have this type of sensor hooked up; confirm this by clicking on Connect.

Now when you hit the Collect button, the sensor will start graphing the number of counts it receives during successive 5-second intervals. The y axis of the graph is counts per 5 seconds, and the x axis is time.

Once you've made some preliminary observations, try to get a good measurement of the background, counting for several minutes so as to reduce the statistical errors. It is important to get a good measurement of the background rate, because in later parts of the lab you'll need to subtract it from all the other count rates you measure. For longer runs like this, it is convenient to let the software collect data for about a minute at a time, rather than 5 seconds. To get it to do this, do Experiment : Data collection : 60 seconds/sample.

### B Type of radiation

Your first task in the lab is to figure out whether the source emits alpha, beta, or gamma radiation — or perhaps some mixture of these. It is up to you to decide how to do this, but essentially you want to use the fact that they are absorbed differently in matter; referring to your textbook, you'll see that the historical labels  $\alpha$ ,  $\beta$ , and  $\gamma$  were assigned purely on the basis of the differences in absorption, before anyone even knew what they really were. Note that the fact that you are able to detect the radiation at all implies that at least some of the radiation is able to penetrate the thin plastic walls of the source.

Also keep in mind that the descriptions of absorption in a textbook are generalizations that do not take into account the energy of the particles. For example, low-energy betas could be absorbed by a kleenex, whereas high-energy betas could penetrate cardboard.

There are six things you could try to prove through your measurements:

1. The source inside the plastic container emits alphas.
2. The source does not emit alphas.
3. The source emits betas.
4. The source does not emit betas.
5. The source emits gammas.
6. The source does not emit gammas.

Try to figure out which of these six statements you can either definitively prove or definitively disprove. Because you are not allowed to extract the source from the packaging, there will be some cases in which you cannot draw any definitive conclusion one way or the other.

### C Distance dependence

Measure the count rate at several different distances from the source. The goal is to find the mathematical form of this function (see Analysis, below). Distances of less than about 10 cm do not work well, because the size of the GM tube is comparable to 10 cm.


### D Absorption

Measure the reduction in count rate when a 25 mm thick wood block is interposed between the source and the detector, and likewise for 17 mm of steel. Keep the distance from the source to the detector constant throughout. Let the counter run for at least five minutes each time. Based on these observations, predict the count rate you would get with two 17-mm thicknesses of steel instead of one. Test your prediction.

### E Decay curve

The source consists of a particular isotope of cesium; we'll refer to it as  $^{NNN}\text{Cs}$ , since the main goal of this lab is to determine what the unknown isotope actually is. It decays to an isotope of barium, and rather than decaying to the ground state of the barium nucleus, it nearly always decays to an excited state,

which then emits the radiation you characterized in part B. Although the half-life of the parent cesium isotope is many years, the half-life of the excited isotope in the barium daughter is short enough that it can be observed during the lab period. However, if the cesium and barium are not separated, then no time variation will be observed, because the supply of barium nuclei is being continuously replenished by decay of the parent cesium.

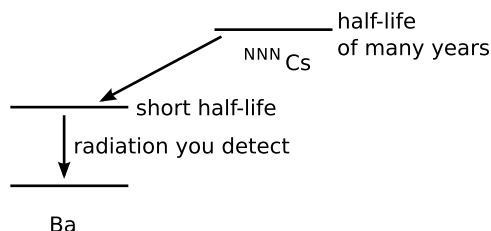

To get around this, the source is packaged so that when a weak acid solution is forced through it, a small amount of barium is washed out. Note that the yellow tape around the circumference of the source has an arrow on it. This arrow points in the direction that you're supposed to make the solution flow. The isotope generator kit has coin-sized steel trays on which to collect a few drops of the radioactive solution.

Use the syringe to draw 1 mL of the acid solution from the 250-mL bottle (labeled "eluting solution"). Take one of the tiny coin-sized steel trays out of the isotope generator kit and lay it on the lab bench. Remove the little plugs from the top and bottom of the radioactive source. Stick the syringe in the in-flow hole, and use the plunger to force seven drops of liquid out onto the tray. Note that the amount of liquid that flows from the syringe into the source is quite a bit more than the amount that comes out into the tray. If you have solution left in the syringe at the end of lab, squeeze it back out into the 250-mL bottle.

Use the computer to collect data on the rate of decay as a function of time. About 5 or 10 seconds per sample works well.

When you're done, make sure to shut off the GM counter so that its battery doesn't get drained.

## Waste disposal

To get an idea of what a non-issue radioactive waste disposal is in this lab, recall that it would be legal to throw the *entire* source in the trash — although we won't actually do that. The amount of radioactive material that you wash out in part E is a tiny



fraction of this. Furthermore, essentially all the radioactivity is gone by the end of lab. It is therefore not a problem to dump your seven drops of material down the drain at the end of class.

There is also no chemical disposal issue with this tiny amount of solution. It's a few drops of very dilute acid, equivalent to a little spritz of lemon juice.

## Analysis

In part C, you should first subtract the background rate from each datum. Then make a log-log plot as described in appendix 5, and see if you can successfully describe the data using a power law. Note that, just like a human, the GM counter cannot count faster than a certain rate. This is because every time it gets a count it completely discharges its voltage, and then it has to recharge itself again. For this reason, it is possible that your data from very small distances will not agree with the behavior of the data at larger distances. The documentation for these GM counters says that they can count at up to about 3500 counts per second; this is only a very rough guide, but it gives you some idea what count rates should be expected to start departing from ideal behavior.

Estimate the half-lives of any isotopes present in the data. If you find that only one half-life is present, you can simply determine the amount of time required for the count rate to fall off by a factor of two. If the natural background count rate measured in part A is significant, you will need to subtract it from the raw data. If more than one half-life is present, try plotting the logarithm of the count rate as a function of time, and seeing if there are linear sections on the graph. Note that this is all referring to the half-lives of any decay chain that occurs *after* the cesium decays to barium. The half-life of the cesium parent nucleus is much longer, and is not measured in this experiment.

Consult the Wikipedia article "Isotopes of cesium." It has an extremely lengthy table of all the known isotopes from very light ones (with far too few neutrons to be stable) to very heavy ones (with far too many). Since the source was shipped to us through the mail, and sits on the shelf in the physics stockroom for semester after semester, you can tell that the half-life of the cesium isotope must be fairly long — at least on the order of years, not months. From this information, you should be able to narrow down the range of possibilities. (Half-lives in units of years are listed with "a," for "annum," as the unit of time.

The notations  $m1$  and  $m2$  mean energy states that are above the lowest-energy state.) The radioactive isotopes from this remaining list of possibilities all have their own Wikipedia articles, and these articles give the properties of the daughter nuclei (isotopes of barium and xenon), including the half-lives of any gamma-emitting states. Look for one that has a half-life that seems to match the one you measured in lab. Having tentatively identified the unknown isotope as this isotope of cesium, check against the results of part B, where you determined the type of radiation emitted.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** Suppose that in part C you obtain the following data:

r (cm)	count rate (counts/2 min)
10	707
20	207
30	95

Suppose that the background rate you measure in part A is 30 counts per 2 min. Use the technique described in appendix 5 to see if the data can be described by a power law, and if so, determine the exponent.

**P2** If a source emits gamma or beta radiation, then the radiation spreads out in all directions, like an expanding sphere. Based on the scaling of a sphere's surface area with increasing radius, how would you theoretically expect the intensity of the radiation to fall off with distance? Would it be a power law? If so, what power? Why would you *not* expect the same behavior for a source emitting alpha particles into air?

# 10 Two-Source Interference

## Apparatus


ripple tank	1/group
yellow foam pads	4/group
lamp and unfrosted straight-filament bulb	
1/group wave generator	1/group
big metal L-shaped arms for hanging the wave generator	
the wave generator	1/group
little metal L-shaped arms with yellow plastic balls	2/group
rubber bands	2/group
Thornton DC voltage source (in lab bench)	1/group
small rubber stopper	1/group
power strip	1/group
bucket	1/group
mop	1
flathead screwdriver	1
rulers and protractors	
kimwipes and alcohol for cleaning	
butcher paper	

## Goals

Observe how a 2-source interference pattern of water waves depends on the distance between the sources.


## Observations

Light is really made of waves, not rays, so when we treated it as rays, we were making an approximation. You might think that when the time came to treat light as a wave, things would get very difficult, and it would be hard to predict or understand anything without doing complicated calculations.

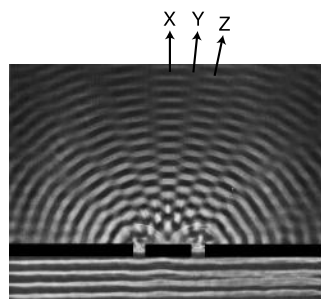

Life isn't that bad. It turns out that all of the most important ideas about light as a wave can be seen

in one simple experiment, shown in the first figure.<sup>1</sup> A wave comes up from the bottom of the page, and encounters a wall with two slits chopped out of it. The result is a fan pattern, with strong wave motion coming out along directions like X and Z, but no vibration of the water at all along lines like Y. The reason for this pattern is shown in the second figure. The two parts of the wave that get through the slits create an overlapping pattern of ripples. To get to a point on line X, both waves have to go the same distance, so they're in step with each other, and reinforce. But at a point on line Y, due to the unequal distances involved, one wave is going up while the other wave is going down, so there is cancellation. The angular spacing of the fan pattern depends on both the wavelength of the waves,  $\lambda$ , and the distance between the slits,  $d$ .

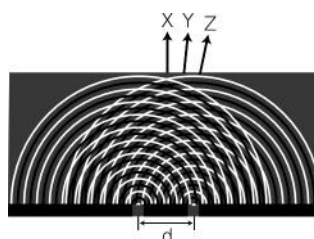

The ripple tank is tank that sits about 30 cm above the floor. You put a little water in the tank, and produce waves. There is a lamp above it that makes a point-like source of light, and the waves cast patterns of light on a screen placed on the floor. The patterns of light on the screen are easier to see and measure than the ripples themselves.

In reality, it's not very convenient to produce a double-slit diffraction pattern exactly as depicted in the first figure, because the waves beyond the slits are so weak that they are difficult to observe clearly. Instead, you'll simply produce synchronized circular ripples from two sources driven by a motor.

Put the tank on the floor. Plug the hole in the side of the tank with the black rubber stopper. If the plastic is dirty, clean it off with alcohol and kimwipes. Wet the four yellow foam pads, and place them around the sides of the tank. Pour in water to a depth of about 5-7 mm. Adjust the metal feet to level the tank, so that the water is of equal depth throughout the tank. (Do not rotate the wooden legs them-

<sup>1</sup>The photo is from the textbook PSSC Physics, which has a blanket permission for free use after 1970.

selves, just the feet.) If too many bubbles form on the plastic, wipe them off with a ruler.

Make sure the straight-filament bulb in the light source is rotated so that when you look in through the hole, you are looking along the length of the filament. This way the lamp acts like a point source of light above the tank. To test that it's oriented correctly, check that you can cast a perfectly sharp image of the tip of a pen.

The light source is intended to be clamped to the wooden post, but I've found that that works very poorly, since the clamp doesn't hold it firmly enough. Instead, clamp the light source to the lab bench's lip or its leg. Turn it on. Put the butcher paper on the floor under the tank. If you make ripples in the water, you should be able to see the wave pattern on the screen.

The wave generator consists of a piece of wood that hangs by rubber bands from the two L-shaped metal hangers. There is a DC motor attached, which spins an intentionally unbalanced wheel, resulting in vibration of the wood. The wood itself can be used to make straight waves directly in the water, but in this experiment you'll be using the two little L-shaped pieces of metal with the yellow balls on the end to make two sources of circular ripples. The DC motor runs off of the DC voltage source, and the more voltage you supply, the faster the motor runs.


Start just by sticking one little L-shaped arm in the piece of wood, and observing the circular wave pattern it makes. Now try two sources at once, in neighboring holes. Pick a speed (frequency) for the motor that you'll use throughout the experiment — a fairly low speed works well. Measure the angular spacing of the resulting diffraction pattern for several values of the spacing,  $d$ , between the two sources of ripples.

Use the methods explained in Appendix 5 and look for any kind of a power law relationship for the dependence of the angular spacing on  $d$ .

# 11 Refraction and Images

## Apparatus

plastic box  
propanol (1 liter/group, to be reused)  
laser  
spiral plastic tube and fiber optic cable for demonstrating total internal reflection  
ruler  
protractor  
butcher paper  
funnel

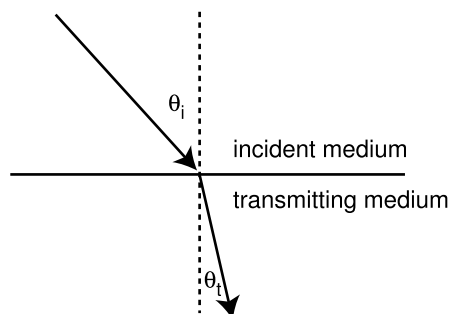

## Goals

Test whether the index of refraction of a liquid is proportional to its density.

Observe the phenomena of refraction and total internal reflection.

Locate a virtual image in a plastic block by ray tracing, and compare with the theoretically predicted position of the image.

## Introduction

Without the phenomenon of refraction, the lens of your eye could not focus light on your retina, and you would not be able to see. Refraction is the bending of rays of light that occurs when they pass through the boundary between two media in which the speed of light is different.

Refraction occurs for the following reason. Imagine, for example, a beam of light entering a swimming pool at an angle. Because of the angle, one side of the beam hits the water first, and is slowed down. The other side of the beam, however, gets to travel in air, at its faster speed, for longer, because it enters the water later — by the time it enters the water, the other side of the beam has been limping along through the water for a little while, and has not gotten as far. The wavefront is therefore twisted around a little, in the same way that a marching band turns by having the people on one side take smaller steps.

Quantitatively, the amount of bending is given by Snell's law:

$$n_i \sin \theta_i = n_t \sin \theta_t,$$

where the subscript  $i$  refers to the incident light and

incident medium, and  $t$  refers to the transmitted light and the transmitting medium. This relation can be taken as defining the quantities  $n_i$  and  $n_t$ , which are known as the indices of refraction of the two media. Note that the angles are defined with respect to the normal, i.e., the imaginary line perpendicular to the boundary.

Also, not all of the light is transmitted. Some is reflected — the amount depends on the angles. In fact, for certain values of  $n_i$ ,  $n_t$ , and  $\theta_i$ , there is no value of  $\theta_t$  that will obey Snell's law ( $\sin \theta_t$  would have to be greater than one). In such a situation, 100% of the light must be reflected. This phenomenon is known as total internal reflection. The word internal is used because the phenomenon only occurs for  $n_i > n_t$ . If one medium is air and the other is plastic or glass, then this can only happen when the incident light is in the plastic or glass, i.e., the light is trying to escape but can't. Total internal reflection is used to good advantage in fiber-optic cables used to transmit long-distance phone calls or data on the internet — light traveling down the cable cannot leak out, assuming it is initially aimed at an angle close enough to the axis of the cable.

Although most of the practical applications of the phenomenon of refraction involve lenses, which have curved shapes, in this lab you will be dealing almost exclusively with flat surfaces.

### Preliminaries

Check whether your laser's beam seems to be roughly parallel.

## Observations

### A Index of refraction of alcohol

The index of refraction is sometimes referred to as the optical density. This usage makes sense, because when a substance is compressed, its index of refraction goes up. In this part of the lab, you will test whether the indices of refraction of different liquids are proportional to their mass densities. Water has a density of  $1.00 \text{ g/cm}^3$  and an index of refraction of 1.33. Propanol has a density of  $0.79 \text{ g/cm}^3$ . You will find out whether its index of refraction is lower than water's in the same proportion. The idea is to pour some alcohol into a transparent plastic box and measure the amount of refraction at the interface between air and alcohol.

Make the measurements you have planned in order to determine the index of refraction of the alcohol. The laser and the box can simply be laid flat on the table. Make sure that the laser is pointing towards the wall.

### B Total internal reflection

Try shining the laser into one end of the spiral-shaped plastic rod. If you aim it nearly along the axis of the cable, none will leak out, and if you put your hand in front of the other end of the rod, you will see the light coming out the other end. (It will not be a well-collimated beam any more because the beam is spread out and distorted when it undergoes the many reflections on the rough and curved inside the rod.)



There's no data to take. The point of having this as part of the lab is simply that it's hard to demonstrate to a whole class all at once.

### C A virtual image


Pour the alcohol back into the container for reuse, and pour water into the box to replace it.

Pick up the block, and have your partner look sideways through it at your finger, touching the surface of the block. Have your partner hold her own finger next to the block, and move it around until it appears to be as far away as your own finger. Her brain achieves a perception of depth by subconsciously comparing the images it receives from her two eyes. Your partner doesn't actually need to be able to see her own finger, because her brain knows how to position her arm at a certain point in space. Measure the distance  $d_i$ , which is the depth of the image of your finger relative to the front of the block.

Next we will use the laser to simulate the rays com-

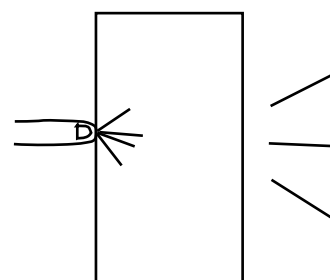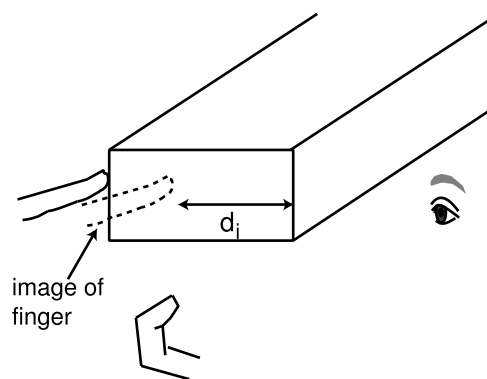

Light rays spreading out from the finger by diffuse reflection. The emerging rays all appear to have come from a point inside the block.

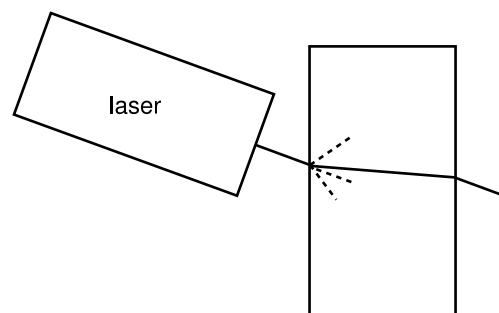

Simulating one of the rays using the laser.

ing from the finger, as shown in the figure. Shine the laser at the point where your finger was originally touching the block, observe the refracted beam, and draw it in. Repeat this whole procedure several times, with the laser at a variety of angles. Finally, extrapolate the rays leaving the block back into the block. They should all appear to have come from the same point, where you saw the virtual image. You'll need to photocopy the tracing so that each person can turn in a copy with his or her writeup.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

Print out, complete, and turn in the laser safety checklist, appendix 7.

**P1** Laser beams are supposed to be very nearly parallel (not spreading out or contracting to a focal point). Think of a way to test, roughly, whether this is true for your laser.

**P2** Plan how you will determine the index of refraction in part A.

**P3** You have complete freedom to choose any incident angle you like in part A. Discuss what choice would give the highest possible precision for the measurement of the index of refraction.

## Analysis

Using your data for part A, extract the index of refraction of propanol, with error bars. Test the hypothesis that the index of refraction is proportional to the density in the case of water and propanol.

Using trigonometry and Snell's law, make a theoretical calculation of  $d_i$ . You'll need to use the small-angle approximation  $\sin \theta \approx \tan \theta \approx \theta$ , for  $\theta$  measured in units of radians. (For large angles, i.e. viewing the finger from way off to one side, the rays will not converge very closely to form a clear virtual image.)

Explain your results in part C and their meaning.

Compare your three values for  $d_i$ : the experimental value based on depth perception, the experimental value found by ray-tracing with the laser, and the theoretical value found by trigonometry.



# 12 Geometric Optics

## Apparatus

optical bench ..... 1/group  
converging lens,  $f = 50$  cm, 5 cm diam  
1/group converging lens,  $f = 5$  cm ..... 1/group  
lamp and arrow-shaped mask ..... 1/group  
frosted glass screen ..... 1/group

## Goals

Observe a real image formed by a convex lens, and determine its focal length.

Construct a telescope and measure its angular magnification.

## Introduction

The credit for invention of the telescope is disputed, but Galileo was probably the first person to use one for astronomy. He first heard of the new invention when a foreigner visited the court of his royal patrons and attempted to sell it for an exorbitant price. Hearing through second-hand reports that it consisted of two lenses, Galileo sent an urgent message to his benefactors not to buy it, and proceeded to reproduce the device himself. An early advocate of simple scientific terminology, he wanted the instrument to be called the “occhialini,” Italian for “eyething,” rather than the Greek “telescope.”

His astronomical observations soon poked some gaping holes in the accepted Aristotelian view of the heavens. Contrary to Aristotle’s assertion that the heavenly bodies were perfect and without blemishes, he found that the moon had mountains and the sun had spots (the marks on the moon visible to the naked eye had been explained as optical illusions or atmospheric phenomena). This put the heavens on an equal footing with earthly objects, paving the way for physical theories that would apply to the whole universe, and specifically for Newton’s law of gravity. He also discovered the four largest moons of Jupiter, and demonstrated his political savvy by naming them the “Medicean satellites” after the powerful Medici family. The fact that they revolved around Jupiter rather than the earth helped make more plausible Copernicus’ theory that the planets did not revolve around the earth but around the sun.

Galileo’s ideas were considered subversive, and many people refused to look through his telescope, either because they thought it was an illusion or simply because it was supposed to show things that were contrary to Aristotle.

The figure on the next page shows the simplest refracting telescope. The object is assumed to be at infinity, so a real image is formed at a distance from the objective lens equal to its focal length,  $f_o$ . By setting up the eyepiece at a distance from the image equal to its own focal length,  $f_E$ , light rays that were parallel are again made parallel.

The point of the whole arrangement is angular magnification. The small angle  $\alpha_1$  is converted to a large  $\alpha_2$ . It is the small angular size of distant objects that makes them hard to see, not their distance. There is no way to tell visually whether an object is a thirty meters away or thirty billion. (For objects within a few meters, your brain-eye system gives you a sense of depth based on parallax.) The Pleiades star cluster can be seen more easily across many light years than Mick Jagger’s aging lips across a stadium. People who say the flying saucer “looked as big as an aircraft carrier” or that the moon “looks as big as a house” don’t know what they’re talking about. The telescope does not make things “seem closer” — since the rays coming at your eye are parallel, the final virtual image you see is at infinity. The angular magnification is given by

$$M_A = \alpha_2 / \alpha_1$$

(to be measured directly in this lab)

$$M_A = f_o / f_E$$


(theory).

## Observations

### A Focal length of the lenses

In this part of the lab, you’ll accurately determine the focal lengths of the two lenses being used for the telescope. They are poorly quality-controlled, and I’ve found that the labeled values are off by as much as 10%. If you’re only doing part A of the lab, then you will only determine the focal length of one lens, which is given to you as an unknown.



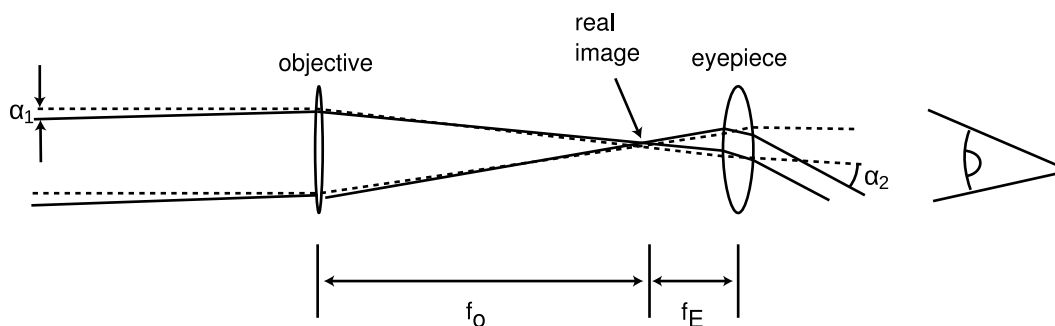

A refracting telescope

Start with the short-focal-length lens you're going to use as your eyepiece. Use the lens to project a real image on the frosted glass screen. For your object, use the lamp with the arrow-shaped aperture in front of it. Make sure to lock down the parts on the optical bench, or else they may tip over and break the optics!

For the long-focal-length lens you're going to use as your objective, you will probably be unable to do a similar determination on a one-meter optical bench. Improvising a similar setup without the bench will still give you a much more accurate value than the one written on the label.

A careful measurement here pays off later by making the focus in part B much easier to find. This is especially true for the longer-focal-length lens. To improve the quality of your result, do the kind of thing they do at the optometrist — “which is better, 1 or 2?” Have several people do independent determinations of the best focus.

## B The telescope

Use your optical bench and your two known lenses to build a telescope. Since the telescope is a device for viewing objects at infinity, you'll want to take it outside.

The best method for determining the angular magnification is to observe the same object with both eyes open, with one eye looking through the telescope and one seeing the object without the telescope. Good precision can be obtained, for example, by looking at a large object like a coke machine, and determining that a small part of it, whose size you can measure with a ruler, appears, when magnified, to cover some larger part of it, which you can also measure. The figures on p. 51 show a simulation of what the superimposed images should look like and of how it would look if the telescope is not yet adjusted quite correctly.

Your brain is not capable of focusing one eye at one

distance, and the other at another distance. Therefore it's important to get your telescope adjusted precisely so that the image is at infinity. You can do this by focusing your naked eye on a distant object, and then moving the objective until the image pops into focus in the other eye. Theoretically this would be accomplished simply by setting the lenses at the distance shown in the diagram, but in reality, a small amount of further adjustment is necessary, because of the uncertainty in the measured focal lengths.

A good quick test of the focus is to pick someone who's nearsighted and see if they can focus on the image without their glasses on. If they can, then the image is not at infinity, because nearsighted people can't focus on an image that's at infinity.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** Skip this question if you were assigned problem 31-22 from *Light and Matter*. In part A, do you want the object to be closer to the lens than the lens' focal length, exactly at a distance of one focal length, or farther than the focal length? What about the screen?


**P2** Plan what measurements you will make in part A and how you will use them to determine the lenses' focal length.

**P3** Skip this question if you are only doing part A of the lab. It's disappointing to construct a telescope with a very small magnification. Given a selection of lenses, plan how you can make a telescope with the greatest possible angular magnification.


## Analysis

Determine the focal length of the two lenses, with error bars.

Find the angular magnification of your telescope from your data, with error bars, and compare with theory. Do they agree to within the accuracy of the measurement? Give a probabilistic interpretation, as in the example in appendix 2. See the example at the end of appendix 3 of how to test for equality between two numbers that both have error bars.

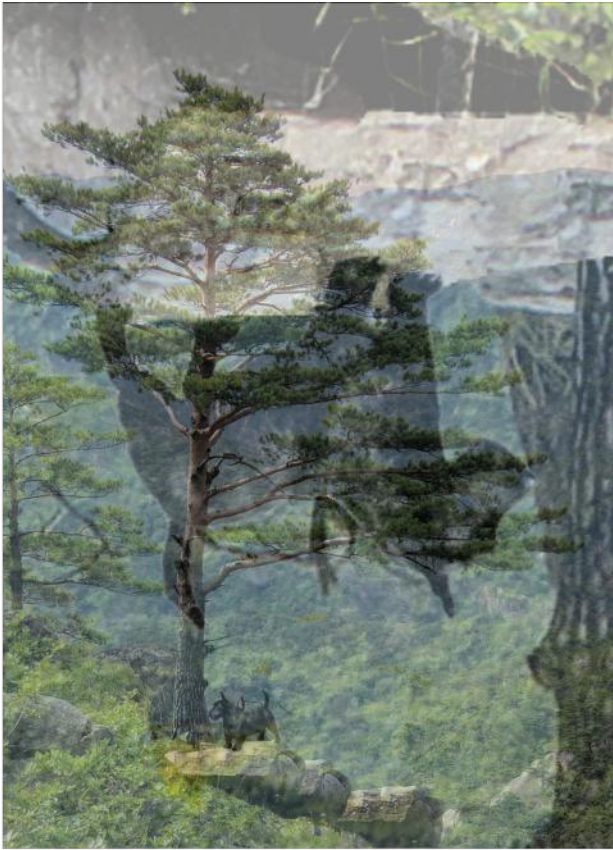

The goal: magnified and unmagnified images superimposed and in focus.

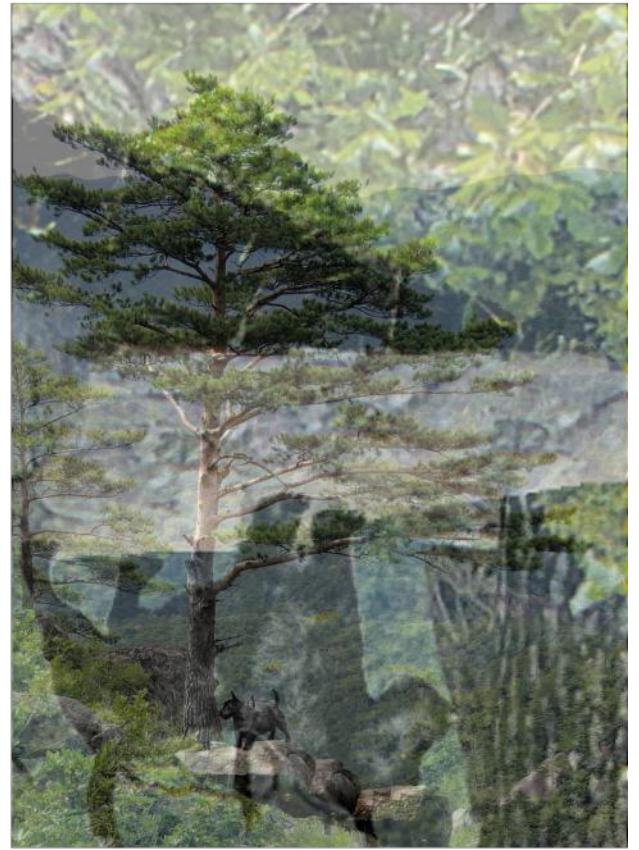

Your aim is wrong.



The telescope's image is not at infinity. Your brain focuses your eye on the dog.



The telescope's image is not at infinity. Your brain focuses your eye on the tree.

These images can be viewed in a web browser at [lightandmatter.com/lec/telescope](http://lightandmatter.com/lec/telescope) . See the file NOTES in the source code for photo credits.

# 13 Wave Optics

## Apparatus

helium-neon laser  
1/group optical bench with posts & holders 1/group  
high-precision double slits ..... 1/group  
rulers  
meter sticks  
tape measures  
butcher paper  
black cloths for covering light sources

## Goals

Observe evidence for the wave nature of light.

Determine the wavelength of the red light emitted by your laser, by measuring a double-slit diffraction pattern. (The part of the spectrum that appears red to the human eye covers quite a large range of wavelengths. A given type of laser, e.g., He-Ne or solid-state, will produce one very specific wavelength.)

Determine the approximate diameter of a human hair, using its diffraction pattern.

## Introduction

Isaac Newton's epitaph, written by Alexander Pope, reads:

Nature and Nature's laws lay hid in night.

God said let Newton be, and all was light.

Notwithstanding Newton's stature as the greatest physical scientist who ever lived, it's a little ironic that Pope chose light as a metaphor, because it was in the study of light that Newton made some of his worst mistakes. Newton was a firm believer in the dogma, then unsupported by observation, that matter was composed of atoms, and it seemed logical to him that light as well should be composed of tiny particles, or "corpuscles." His opinions on the subject were so strong that he influenced generations of his successors to discount the arguments of Huygens and Grimaldi for the wave nature of light. It was not until 150 years later that Thomas Young demonstrated conclusively that light was a wave.

Young's experiment was incredibly simple, and could probably have been done in ancient times if some savvy Greek or Chinese philosopher had only thought of it. He simply let sunlight through a pinhole in a window shade, forming what we would now call a coherent beam of light (that is, a beam consisting of plane waves marching in step). Then he held a thin card edge-on to the beam, observed a diffraction pattern on a wall, and correctly inferred the wave nature and wavelength of light. Since Roemer had already measured the speed of light, Young was also able to determine the frequency of oscillation of the light.

Today, with the advent of the laser, the production of a bright and coherent beam of light has become as simple as flipping a switch, and the wave nature of light can be demonstrated very easily. In this lab, you will carry out observations similar to Young's, but with the benefit of hindsight and modern equipment.

## Observations

### A Determination of the wavelength of red light

Set up your laser on your optical bench. You will want as much space as possible between the laser and the wall, in order to let the diffraction pattern spread out as much as possible and reveal its fine details.

Tear off two small scraps of paper with straight edges. Hold them close together so they form a single slit. Hold this improvised single-slit grating in the laser beam and try to get a single-slit diffraction pattern. You may have to play around with different widths for the slit. No quantitative data are required. This is just to familiarize you with single-slit diffraction.

Make a diffraction pattern with the double-slit grating. See what happens when you hold it in your hand and rotate it around the axis of the beam.

The diffraction pattern of the double-slit grating consists of a rapidly varying pattern of bright and dark bars, with a more slowly varying pattern superimposed on top. (See the figure two pages after this page.) The rapidly varying pattern is the one that is numerically related to the wavelength,  $\lambda$ , and the distance between the slits,  $d$ , by the equation

$$\Delta\theta = \lambda/d,$$

This spacing,  $\Delta\theta$ , is the one you want to measure. It is related to  $d$ , the center-to-center distance between the slits, by  $\Delta\theta = \lambda/d$ .



If your screen is too close to the slits, the finer pattern may be invisible, and this spacing may be all you see. This larger angular spacing is related to the width of the slits, not to  $d$ .

A double-slit diffraction pattern.


where  $\theta$  is measured in radians. To make sure you can see the fine spacing, put your slits across the room from the wall. To make it less likely that someone will walk through the beam and get the beam in their eye, put some of the small desks under the beam. The slit patterns we're using actually have three sets of slits, with the following dimensions:

	$w$ (mm)	$d$ (mm)
A	.12	.6
B	.24	.6
C	.24	1.2

The small value of  $d$  is typically better, for two reasons: (1) it produces a wider diffraction pattern, which is easier to see; (2) it's easy to get the beam of the laser to cover both slits.

If your diffraction pattern doesn't look like the one in the figure on the following page, typically the reason is that you're only hitting one slit with the beam (in which case you get a single-slit diffraction pattern), or you're not illuminating the two slits equally (giving a funny-looking pattern with little dog-bones and things in it).

As shown in the figure below, it is also possible to have the beam illuminate only *part* of each slit, so that the slits act effectively as if they had a smaller value of  $d$ . The beam spreads as it comes out of the laser, so you can avoid this problem by putting it fairly far away from the laser (at the far end of the optical bench).

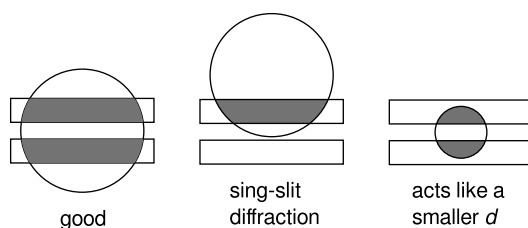

Think about the best way to measure the spacing of the pattern accurately. Is it best to measure from a

bright part to another bright part, or from dark to dark? Is it best to measure a single spacing, or take several spacings and divide by the number to find what one spacing is?

Determine the wavelength of the light, in units of nanometers. Make sure it is in the right range for red light.

Check that the  $\Delta\theta$  you obtain is in the range predicted in prelab question P1. In the past, I've seen cases where groups got goofy data, and I suspect that it was because they were hitting a place on the slits where there was a scratch, bump, or speck of dust.

## B Diameter of a human hair

Pull out one of your own hairs, hold it in the laser beam, and observe a diffraction pattern. It turns out that the diffraction pattern caused by a narrow obstruction, such as your hair, has the same spacing as the pattern that would be created by a single slit whose width was the same as the diameter of your hair. (This is an example of a general theorem called Babinet's principle.) Measure the spacing of the diffraction pattern. (Since the hair's diameter is the only dimension involved, there is only one diffraction pattern with one spacing, not superimposed fine and coarse patterns as in part A.) Determine the diameter of your hair. Make sure the value you get is reasonable, and compare with the order-of-magnitude guess you made in your prelab writeup.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise

you're just setting yourself up for failure in lab.

Read the safety checklist.

**P1** Look up the approximate range of wavelengths that the human eye perceives as red. With  $d = 0.6$  mm, predict the lowest and highest possible values of  $\Delta\theta$  that could occur with red light.

**P2** It is not practical to measure  $\Delta\theta$  directly using a protractor. Suppose that a lab group finds that 27 fringes extend over 29.7 cm on their butcher paper, which is on a wall 389 cm away from the slits. They calculate  $\Delta\theta = \tan^{-1}(29.7 \text{ cm}/(27 \times 389 \text{ cm})) = 2.83 \times 10^{-3}$  rad. Simplify this calculation using a small-angle approximation, and show that the resulting error is negligible.

**P3** Make a rough order-of-magnitude guess of the diameter of a human hair.

## Analysis

Determine the wavelength of the light and the diameter of the hair, with error bars.



# 14 The Photoelectric Effect

Note to the lab technician: The Pasco SE-5509 Hg gas discharge tubes have too many ways of being destroyed by students,<sup>1</sup> so let's use the OS-9286 tubes instead (the big black ones with the fins). That means we don't need the separate power supply for the discharge tube (Pasco BEM-5007) or the track that comes with the photodiode. We need wood blocks to raise the photodiode to the same height as the discharge tube.

## Apparatus

hand-held diffraction gratings  
Hg gas discharge tube (Pasco OS-9286)  
photodiode (Pasco SE-5509)  
power supply (Pasco BEM-5001)  
high-sensitivity ammeter (Pasco BEM-5004)  
wood blocks ..... 3/group

## Goals

Use the photoelectric effect to test predictions of the wave and wave-particle models of light.

## Introduction

The photoelectric effect, a phenomenon in which light shakes an electron loose from an object, provided the first evidence for wave-particle duality: the idea that the basic building blocks of light and matter show a strange mixture of particle and wave behaviors. At the turn of the twentieth century, physicists assumed that particle and wave phenomena were completely distinct. Young had shown that light could undergo interference effects such as diffraction, so it had to be a wave. Since light was a wave composed of oscillating electric and magnetic fields, it made sense that when light encountered matter, it would tend to shake the electrons. It was only to be expected that something like the photoelectric effect could happen, with the light shaking the electrons vigorously enough to knock them out of the atom.

<sup>1</sup>The two main modes of destruction seem to be: (1) While the discharge tube is connected to the power supply, students monkey with the red 110V/220V switch on the power supply; and (2) they don't use the power supply and connect two AC connectors together in order to connect the discharge tube directly to the wall socket.

But once the effect was observed, physicists began running into trouble interpreting how it behaved. There were various variables they could adjust, such as:

- the light's frequency (color), and
- the light's intensity (brightness).

Given these input conditions, there were outputs they could look at, including

- any time delay before electrons began to pop out,
- the rate at which the electrons then flowed (measured as a current on an ammeter), and
- the amount of kinetic energy they had.

At the time this was considered an obscure technical topic, but experimentalists began generating data, which theorists then had zero success in interpreting. Albert Einstein, better known today for the theory of relativity, was the first to come up with the radical, and correct, explanation, which involved a fundamental rewriting of the laws of physics.

## Setup

The Hg gas discharge tube emits light with several different wavelengths. Turn on the discharge tube immediately, because it takes a long time to warm up.

The photodiode is a vacuum tube housed inside another box, with a small hole to allow light to come in and hit one of the electrodes (the cathode) inside the vacuum tube. On the front of the box, covering the hole, are two rotating wheels. The wheel that you can see has five colored filters. Each of these filters lets through only one of the five wavelengths of light emitted by the Hg tube, so that you can control the frequency.

The hole through which the light enters is actually a hole in a second wheel, located behind the filter wheel. By clicking that wheel into different position we can select holes of different sizes, which lets in different amounts of light. Although this is convenient, it doesn't actually control the intensity of the light



(watts per square meter) but only the total amount (watts). Therefore we'll just leave this set for the 2 mm hole.

There is an easier way to control the actual intensity of the beam, which is that you can simply change the distance between the discharge tube and the photodiode. The light from the discharge tube spreads out in a cone, so that as you vary the distance to the photodiode, the intensity falls off as the inverse square of the distance.

There is a power supply used for applying an external voltage to the photodiode. Although you don't actually want to apply that external voltage during this part of the lab, it seems that the ammeter won't work until you hook up the power supply, so you need to do that now. On the bottom right side of the power supply are two banana plugs. Connect these to the two plugs near the bottom of the photodiode. To get the polarity right, connect the positive (red) output of the power supply to the anode (marked A).


The photodiode has an output that can be connected to a extremely sensitive ammeter to measure the rate at which electrons are ejected from the cathode and absorbed at the anode. Before connecting the ammeter (labeled "DC Current Amplifier" on the front panel), set it to its most sensitive scale ( $10^{-13}$  A) and depress the button labeled "Calibration." Since the meter is not connected to anything, the current is truly zero. Use the knob to force the meter to read zero, and then let the calibration button out. Now use a BNC cable to connect the photodiode to the ammeter.

Once the discharge tube has warmed up, arrange things so that you can see a spot of its light, e.g., by letting it fall on a white piece of paper. Hold the diffraction grating up to your eye and look at the light. In the first-order ( $m = \pm 1$ ) fringes, you should be able to see that the light contains a mixture of four discrete wavelengths of visible light. There is also an ultraviolet wavelength, which you may be able to see as well if the paper fluoresces.

Now suppose we were considering the following possible models of light: (1) pure particle model, (2) pure wave model, (3) a hybrid model in which light has both particle and wave properties. You have just observed diffraction of the light from your source. Which of the models are consistent with this observation, and which are immediately ruled out?

The wavelengths are as follows:

color	wavelength (nm)
ultraviolet	365
violet	405
blue	436
green	546
yellow	578

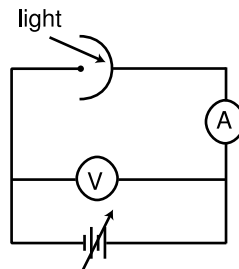

The circuit. The light creates an electric current, which trickles back through the large but finite resistance of the voltmeter.

### Circuit

The figure above shows a circuit diagram of the setup. Light comes in and knocks electrons out of the curved cathode. If the voltage is turned off, there is no electric field, so the electrons travel in straight lines; some will hit the anode, creating a current. If the voltage is turned on, the electric field repels the electrons from the wire electrode, and the current is reduced or perhaps even completely eliminated.

## Observations

Although the photodiode box has filters on the front, no filter is perfect, and therefore these will all let in some stray light of wavelengths other than the intended one. Therefore the room should be very dark when you do your measurements.

### A Time delay

As explored in the prelab, there may be some delay between the time when the light is allowed to hit the cathode and the time when electrons begin to be ejected. If so, then we lack even a rough a priori estimate of this time. Investigate this. If the time seems to be extremely short, do what you planned in the prelab to try to make it long enough to detect. If the time is much too long, do the opposite so that you can actually observe the photoelectric effect. If you're able to get the time delay into a range where it's measurable using eyeballs and a clock, do so. If not, then try to set an upper or lower limit on it.

## B Energy of the electrons

Until you do the lab, it's not obvious how much energy the electrons would have as they pop out of the cathode. It could be some fixed number, or it could depend on the conditions you choose, and it could also have multiple values for the electrons produced under a given set of conditions. In fact, we may expect a range of values for two reasons.

(1) As explored in a prelab question, the electrons will have random kinetic energies to start with, due to their thermal motion

(2) The light penetrates to some depth in the cathode, and an electron that starts at some depth will lose some amount of energy as it then comes out to the surface. The electron's direction of motion is random. If it happens to be toward the surface, then the thickness of material that it traverses will depend both on its initial depth and on its random direction of motion.

This range of energies may have an upper limit for a given set of experimental conditions. If so, then by applying a high enough voltage you should be able to eliminate the current completely. The minimum voltage required to do this would be called the stopping voltage.

Now you want to apply a voltage to the photodiode using the DC power supply, which you previously connected but didn't use. There is a button between the two LED readouts. Let this button out in order to select a range of voltages from 0 to -4.5 V.

Find out whether there is a stopping voltage, and if so, measure it for the conditions you've chosen. If there is never a sharp cutoff, you should still be able to determine some quantitative measure of the voltage that corresponds to a *typical* energy for the electrons, e.g., the voltage at which some fixed fraction of the current is eliminated. From now on in the lab manual, I'll just refer to this as "the voltage" for a given set of experimental conditions.

Determine the voltage, and compare with the estimate in the prelab of what voltage would correspond to the thermal motion of the electrons. Based on this comparison, is the amount of energy involved in the photoelectric effect much less than, much more than, or on the same order of magnitude as the thermal energy?

## C Dependence of voltage on intensity

Vary the intensity of the light and determine whether and how the voltage depends on intensity.

We have observed strange results sometimes, which seem to occur because when the discharge tube is very close to the photodiode, light gets in and hits the *anode* (not just the cathode) and causes the photoelectric effect in the wrong direction. This seems to occur with the shorter wavelengths of light. You should be able to tell if this is happening because when you turn up the voltage high enough, you get a current in the opposite direction. Check for this behavior and avoid taking data under conditions that produce it.

## D Dependence of voltage on frequency

Do the same thing for the frequency.

## Analysis

The point of the analysis is to try to compare the observed results with our expectations based on two models of light: the pure wave model, and a model in which light is both a particle and a wave — we've seen above that a pure particle model is untenable. For conceptual simplicity, however, we may find it helpful to visualize the wave-particle model as if the light was purely particle-like, so that a beam of light would be like a stream of machine-gun bullets. This is essentially what Einstein does in his 1905 paper. He admits that this is not literally possible, but doesn't attempt a more detailed reconciliation of the particle and wave pictures, which he doesn't know how to achieve. For this reason, the title of the paper refers to the particle picture as a "heuristic," which means a kind of non-rigorous way of getting an answer without using correct fundamental principles.

**Time delay:** You investigated the possible time delay in the wave model in some detail in the prelab and/or homework. In a particle model, what would be your expectations about a time delay? Compare with experiment.

**Energy:** Based on the particle model, we would expect one "bullet" of light to give its fixed amount of energy to one electron, so that under a given set of experimental conditions, there would be a maximum possible kinetic energy for the electrons, achieved when the electron originated very close to the surface of the cathode. In the wave model it is more difficult to make a definite prediction. Did you observe that there was a definite stopping voltage, or that there was no definite cut-off in the current? Does this allow you to test either model?

**Dependence of energy on intensity:** In the particle model, a more intense beam of light would be one

containing a larger number of particles (per unit of cross-sectional area). In the wave model, a more intense beam would be a higher-amplitude wave. Does either model lead you to predict anything specific about the dependence of voltage on intensity? Test against experiment.

Dependence of energy on frequency: Does the (typical or maximum) energy of the electrons  $eV$  depend on the light's frequency  $f$ ? If so, then in the particle model this probably means that we're observing a change in the amount of energy per particle of light. We can't just equate the electron energy  $eV$  to the energy of the particles of light, because the electrons lose a fixed amount of energy (called the work function) as they emerge from the surface of the metal. We can get around the issue of this constant offset by finding the slope of  $eV$  versus  $f$ ; a constant offset on the  $y$  axis doesn't affect the slope of a graph. Let's call this constant  $h$ . Estimate its numerical value. What units does it have? There is no such universal constant in any of the classical laws of physics, so if it pops up here, it indicates that some entirely new physical theory is being probed.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** Under the hypothesis that light is purely a wave, the energy of a beam of light would be smoothly and continuously accumulated by whatever it hit. Therefore it should take some amount of time before an electron could accumulate enough energy to pop out of the metal. You may have estimated this time scale in a homework problem (*Light and Matter* problem 34-11), but that calculation depends on many crude assumptions and rough estimates, so it's hard to know whether to trust it even as an order of magnitude estimate. Therefore if we want to try to observe this time delay in this experiment, it's impossible to know in advance what to look for, and it may be either too short (in which case we can't measure it) or too long (in which case it will look like the apparatus simply isn't working). Suppose that the time delay is too short to detect in the conditions you initially pick. How could you change the conditions in order to make it longer, and possibly detectable?

**P2** Heat is the random motion of the microscopic particles of which matter is composed. Depending on the level of detail in the treatment of thermodynamics in your first-semester course, you may know that the typical energy per particle can be calculated as  $kT$ , where  $T$  is the absolute temperature and  $k$  is a constant called Boltzmann's constant (not to be confused with the Coulomb constant). The result at room temperature is on the order of  $10^{-20}$  joules per particle. In the photoelectric effect, an electron will absorb some additional amount of energy from the light, which is enough to pop it out of the cathode. Until doing the experiment, we do not know how much this additional amount is, but during the lab you will be able to probe this question by using a voltage  $V$  to try to stop the electrons from making it across the gap. If this additional energy was on the *same* order of magnitude as the thermal energy (which it may not be), estimate the voltage required.

**P3** In this experiment, the light comes out of the discharge tube in a spreading cone. Geometrically, how should the intensity  $I$  of the light depend on the distance  $r$ ? State a proportionality.

**P4** Look up the wavelength of visible light and the typical distance between atoms in a solid. How do they compare? In the wave model, should we expect a particular wave-train to hit one atom, or many?

# 15 Electron Diffraction

## Apparatus

cathode ray tube (Leybold 555 626)  
high-voltage power supply (new Leybold)  
100-k $\Omega$  resistor with banana-plug connectors  
Vernier calipers


## Goals

Observe wave interference patterns (diffraction patterns) of electrons, demonstrating that electrons exhibit wave behavior as well as particle behavior.

Learn what it is that determines the wavelength of an electron.

## Introduction

The most momentous discovery of 20th-century physics has been that light and matter are not simply made of waves or particles — the basic building blocks of light and matter are strange entities which display both wave and particle properties at the same time. In our course, we have already learned about the experimental evidence from the photoelectric effect showing that light is made of units called photons, which are both particles and waves. That probably disturbed you less than it might have, since you most likely had no preconceived ideas about whether light was a particle or a wave. In this lab, however, you will see direct evidence that electrons, which you had been completely convinced were particles, also display the wave-like property of interference. Your schooling had probably ingrained the particle interpretation of electrons in you so strongly that you used particle concepts without realizing it. When you wrote symbols for chemical ions such as  $\text{Cl}^-$  and  $\text{Ca}^{2+}$ , you understood them to mean a chlorine atom with one excess electron and a calcium atom with two electrons stripped off. By teaching you to count electrons, your teachers were luring you into the assumption that electrons were particles. If this lab's evidence for the wave properties of electrons disturbs you, then you are on your way to a deeper understanding of what an electron really is — both a particle and a wave.

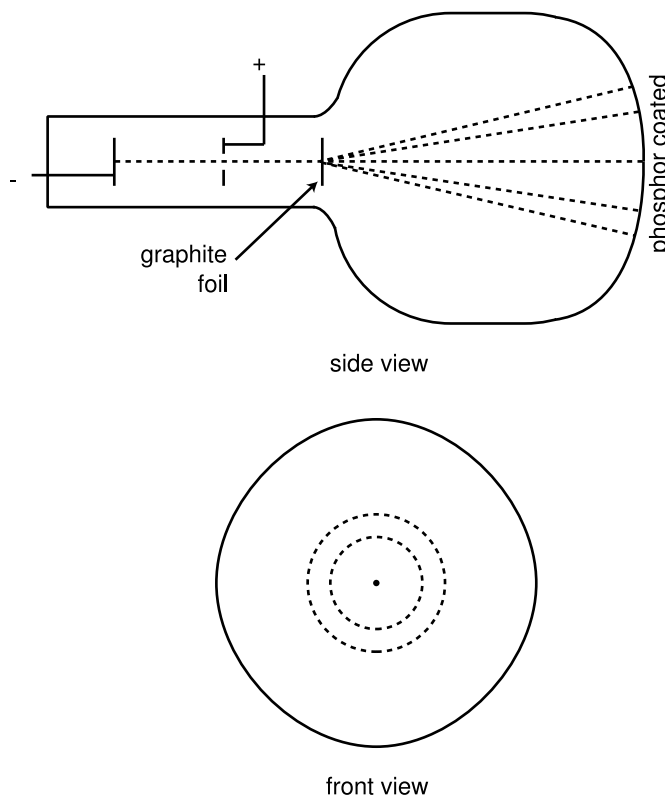

The electron diffraction tube. The distance labeled as 13.5 cm in the figure actually varies from about 12.8 cm to 13.8 cm, even for tubes that otherwise appear identical. This doesn't affect your results, since you're only searching for a proportionality.

## Method

What you are working with is basically the same kind of vacuum tube as the picture tube in your television. As in a TV, electrons are accelerated through a voltage and shot in a beam to the front (big end) of the tube, where they hit a phosphorescent coating and produce a glow. You cannot see the electron beam itself. There is a very thin carbon foil (it looks like a tiny piece of soap bubble) near where the neck joins the spherical part of the tube, and the electrons must pass through the foil before crossing over to the phosphorescent screen.


The purpose of the carbon foil is to provide an ultra-fine diffraction grating — the “grating” consists of the crystal lattice of the carbon atoms themselves! As you will see in this lab, the wavelengths of the

electrons are very short (a fraction of a nanometer), which makes a conventional ruled diffraction grating useless — the closest spacing that can be achieved on a conventional grating is on the order of one micrometer. The carbon atoms in graphite are arranged in sheets, each of which consists of a hexagonal pattern of atoms like chicken wire. That means they are not lined up in straight rows, so the diffraction pattern is slightly different from the pattern produced by a ruled grating.

Also, the carbon foil consists of many tiny graphite crystals, each with a random orientation of its crystal lattice. The net result is that you will see a bright spot surrounded by two faint circles. The two circles represent cones of electrons that intersect the phosphor. Each cone makes an angle  $\theta$  with respect to the central axis of the tube, and just as with a ruled grating, the angle is given by

$$\sin \theta = \lambda/d,$$

where  $\lambda$  is the wavelength of the wave. For a ruled grating,  $d$  would be the spacing between the lines. In this case, we will have two different cones with two different  $\theta$ 's,  $\theta_1$  and  $\theta_2$ , corresponding to two different  $d$ 's,  $d_1$  and  $d_2$ . Their geometrical meaning is shown below.<sup>1</sup>

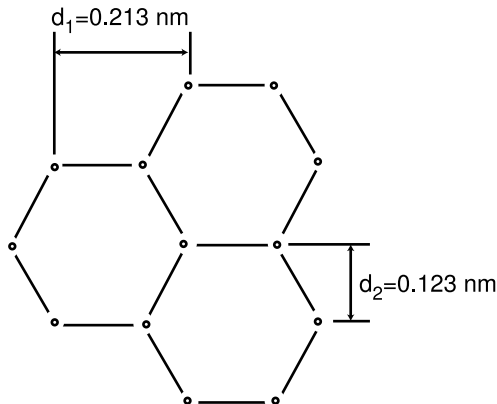

The carbon atoms in the graphite crystal are arranged hexagonally.

### Safety

This lab involves the use of voltages of up to 6000 V. Do not be afraid of the equipment, however; there is a fuse in the high-voltage supply that limits the amount of current that it can produce, so it is not particularly dangerous. Read the safety checklist on

<sup>1</sup>See <http://bit.ly/XxoEYr> for more information.

high voltage in Appendix 6. Before beginning the lab, make sure you understand the safety rules, initial them, and show your safety checklist to your instructor. If you don't understand something, ask your instructor for clarification.

In addition to the high-voltage safety precautions, please observe the following rules to avoid damaging the apparatus:

----- The tubes cost \$1000. Please treat them with respect! Don't drop them! Dropping them would also be a safety hazard, since they're vacuum tubes, so they'll implode violently if they break.

----- Do not turn on anything until your instructor has checked your circuit.

----- Don't operate the tube continuously at the highest voltage values (5000-6000 V). It produces x-rays when used at these voltages, and the strong beam also decreases the life of the tube. You can use the circuit on the right side of the HV supply's panel, which limits its own voltage to 5000 V. Don't leave the tube's heater on when you're not actually taking data, because it will decrease the life of the tube.

## Setup

Your setup will consist of two circuits, a heater circuit and the high-voltage circuit.

The heater circuit is to heat the cathode, increasing the velocity with which the electrons move in the metal and making it easier for some of them to escape from the cathode. This will produce the friendly and nostalgia-producing yellow glow which is characteristic of all vacuum-tube equipment. The heater is simply a thin piece of wire, which acts as a resistor when a small voltage is placed across it, producing heat. Connect the heater connections, labeled F1 and F2, to the 6-V AC outlet at the back of the HV supply.

The high-voltage circuit's job is to accelerate the electrons up to the desired speed. An electron that happens to jump out of the cathode will head "down-hill" to the anode. (The anode is at a *higher* voltage than the cathode, which would make it seem like it would be uphill from the cathode to the anode. However, electrons have negative charge, so they're like negative-mass water that flows uphill.) The high voltage power supply is actually two different power supplies in one housing, with a left-hand panel for one and a right-hand panel for the other. Connect the anode (A) and cathode (C) to the right-hand

panel of the HV supply, and switch the switch on the HV supply to the right, so it knows you're using the right-hand panel.

The following connections are specified in the documentation, although I don't entirely understand what they're for. First, connect the electrode X to the same plug as the cathode.<sup>2</sup> Also, connect F1 to C with the wire that has the 100-k $\Omega$  resistor spliced into it. The circuit diagram on page 63 summarizes all this.

Check your circuit with your instructor before turning it on!

## Observations

You are now ready to see for yourself the evidence of the wave nature of electrons, observe the diffraction pattern for various values of the high voltage, and figure out what determines the wavelength of the electrons. You will need to do your measurements in the dark.

Important: As of 2018, some of our tubes are starting to die, and we will not be able to buy replacements until 2020. For this reason, please take the following steps to extend the remaining lifetimes of the working tubes. (1) Don't take too many data points. Change the voltage in steps of 1.0 kV, not smaller steps. (2) Turn the knob on the high voltage power supply all the way down to zero except when you're actually measuring a diffraction fringe. (3) Try to get all your data-taking done without leaving the heater circuit on for more than about 30 minutes.

You will measure the  $\theta$ 's, and thus determine the wavelength,  $\lambda$ , for several different voltages. Each voltage will produce electrons with a different velocity, momentum, and energy.

Hints:

While measuring the diffraction pattern, don't touch the vacuum tube — the static electric fields of one's body seem to be able to perturb the pattern.

It is easiest to take measurements at the highest voltages, where the electrons pack a wallop and make nice bright rings on the phosphor. Start with the highest voltages and take data at lower and lower voltages until you can't see the rings well enough to take precise data. To

<sup>2</sup>If you look inside the tube, you can see that X is an extra electrode sandwiched in between the anode and the cathode. I think it's meant to help produce a focused beam.

get unambiguous results, you'll need to take data with the widest possible range of voltages.

## Analysis

Once you have your data, the idea is to plot  $\lambda$  as a function of quantities such as  $KE$ ,  $p$ ,  $1/KE$ , or  $1/p$ . If the graph is a straight line through the origin, then the experiment supports the hypothesis that the wavelength is proportional to that quantity. You can simplify your analysis by leaving out constant factors, and P5 asks you to consider how you can rule out some of these possibilities without having to make all the graphs.

What does  $\lambda$  seem to be proportional to? Your data may cover a small enough range of voltage that more than one graph may look linear. However, only one will be consistent with a line that passes *through the origin*, as it must for a proportionality. This is why it is important to have your graph include the origin.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.


The week before you are to do the lab, briefly familiarize yourself visually with the apparatus.

Read the high voltage safety checklist.

**P1** The figure shows the vacuum tube as having a particular shape, which is a sphere with the foil and phosphor at opposite ends of a diameter. In reality, the tubes we're using now are not quite that shape. To me, they look like they may have been shaped so that the phosphor surface is a piece of a sphere centered on the foil. Therefore the arc lengths across the phosphor can be connected to diffraction angles very simply via the definition of radian measure. Plan how you will do this.

**P2** If the voltage difference across which the electrons are accelerated is  $V$ , and the known mass and charge of the electron are  $m$  and  $e$ , what are the electrons' kinetic energy and momentum, in terms of  $V$ ,  $m$ , and  $e$ ? (As a numerical check on your results, you should find that  $V = 5700$  V gives  $KE = 9.1 \times 10^{-16}$  J and  $p = 4.1 \times 10^{-23}$  kg·m/s.)

**P3** All you're trying to do based on your graphs is

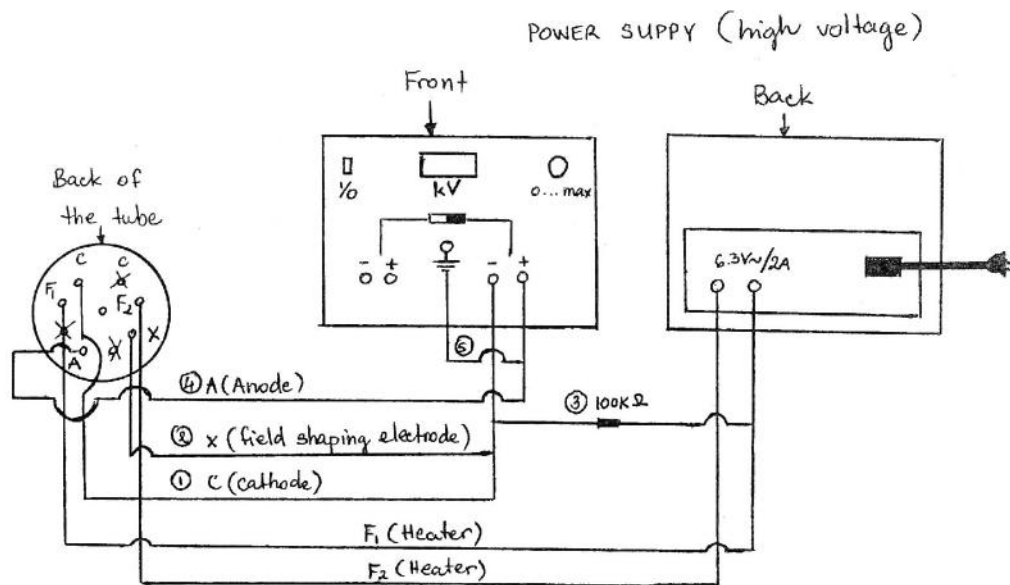

- : ① C (cathode) , ② X (field shaping electrode) , ③ 100K $\Omega$  (resistor) - goes to the back of the power supply
- + : ④ A (Anode) , ⑤ goes to Ground  $\underline{\underline{\text{G}}}$

The circuit for the new setup.

judge which one could be a graph of a proportionality, i.e., a line passing through the origin. Because of this, you can omit any constant factors from the equations you found in P2. When you do this, what do your expressions turn out to be?

**P4** Why is it not logically possible for the wavelength to be proportional to both  $p$  and  $KE$ ? To both  $1/p$  and  $1/KE$ ?

**P5** I have suggested plotting  $\lambda$  as a function of  $p$ ,  $KE$ ,  $1/p$  and  $1/KE$  to see if  $\lambda$  is directly proportional to any of them. Once you have your raw data, how can you immediately rule out two of these four possibilities and avoid drawing the graphs?

**P6** On each graph, you will have two data-points for each voltage, corresponding to two different measurements of the same wavelength. The two wavelengths will be almost the same, but not exactly the same because of random errors in measuring the rings. Should you get the wavelengths by combining the smaller angle with  $d_1$  and the larger angle with  $d_2$ , or vice versa?

# 16a Setup of the Spectrometer

## Apparatus

Hg gas discharge tube (PASCO OS-9286)	.....3
spectrometer	..... 1/group
diffraction grating, 600 lines/mm	..... 1/group
small screwdriver	.....1
black cloth	..... 1
piece of plywood	..... 1
block of wood	..... 1
penlight	..... 1/group
light block	

## Goals

The lab has three parts. This one, part a, is about setting up the optics of the spectrometer. This is to be done once by the instructor or lab technician. It never needs to be redone unless something gets messed up.

## Introduction

### Method

The apparatus is shown in the first figure below. For a given wavelength, the grating produces diffracted light at many different angles: a central zeroth-order line at  $\theta = 0$ , first-order lines on both the left and right, and so on through higher-order lines at larger angles. The line of order  $m$  occurs at an angle satisfying the equation  $m\lambda = d \sin \theta$ .

To measure a wavelength, students will move the telescope until the diffracted first-order image of the slit is lined up with the telescope's cross-hairs and then read off the angle. Note that the angular scale on the table of the spectroscopy actually gives the angle labeled  $\alpha$  in the figure, not  $\theta$ .

### Sources of systematic errors

There are three sources of systematic error:

*angular scale out of alignment:* If the angular scale is out of alignment, then all the angles will be off by a constant amount.

*factory's calibration of  $d$ :* The factory that made the grating labeled it with a certain spacing (in lines per millimeter) which can be con-

verted to  $d$  (center-to-center distance between lines). But their manufacturing process is not all that accurate, so the actual spacing of the lines is a little different from what the label says.

*orientation of the grating:* Errors will be caused if the grating is not perpendicular to the beam from the collimator, or if the lines on the grating are not vertical (perpendicular to the plane of the circle).

### Eliminating systematic errors

A trick to eliminate the error due to misalignment of the angular scale is to observe the same line on both the right and the left, and take  $\theta$  to be half the difference between the two angles, i.e.,  $\theta = (\alpha_R - \alpha_L)/2$ . Because you are subtracting two angles, any source of error that adds a constant offset onto the angles is eliminated. A few of the spectrometers have their angular scales out of alignment with the collimators by as much as a full degree, but that's of absolutely no consequence if this technique is used.

Regarding the calibration of  $d$ , the first person who ever did this type of experiment simply had to make a diffraction grating whose  $d$  was very precisely constructed. But once someone has accurately measured at least one wavelength of one emission line of one element, one can simply determine the spacing,  $d$ , of any grating using a line whose wavelength is known.

You might think that these two tricks would be enough to get rid of any error due to misorientation of the grating, but they're not. They will get rid of any error of the form  $\theta \rightarrow \theta + c$  or  $\sin \theta \rightarrow c \sin \theta$ , but misorientation of the grating produces errors of the form  $\sin \theta \rightarrow \sin \theta + c$ . Part A below describes some additional adjustments that help to get rid of additional sources of error.

## Theory of Operation

The second figure below shows the optics from the side, with the telescope simply looking down the throat of the collimator at  $\theta = 0$ . You are actually using the optics to let you see an image of the slit, not the tube itself. The point of using a telescope is that it provides angular magnification, so that a small change in angle can be seen visually.



A lens is used inside the collimator to make the light from the slit into a parallel beam. This is important, because we are using  $m\lambda = d \sin \theta$  to determine the wavelength, but this equation was derived under the assumption that the light was coming in as a parallel beam. To make a parallel beam, the slit must be located accurately at the focal point of the lens. This adjustment should have already been done, but you will check later and make sure. A further advantage of using a lens in the collimator is that a telescope only works for objects far away, not nearby objects from which the reflected light is diverging strongly. The lens in the collimator forms a virtual image at infinity, on which the telescope can work.

The objective lens of the telescope focuses the light, forming a real image inside the tube. The eyepiece then acts like a magnifying glass to let you see the image. In order to see the cross-hairs and the image of the slit both in focus at the same time, the cross-hairs must be located accurately at the focal point of the objective, right on top of the image.

## Adjustments


First you must check that the cross-hairs are at the focal point of the objective. If they are, then the image of the slits formed by the objective will be at the same point in space as the crosshairs. You'll be able to focus your eye on both simultaneously, and there will be no parallax error depending on the exact position of your eye. The easiest way to check this is to look through the telescope at something far away ( $\gtrsim 50$  m), and move your head left and right to see if the crosshairs move relative to the image. Slide the eyepiece in and out to achieve a comfortable focus. If this adjustment is not correct, you may need to move the crosshairs in or out; this is done by sliding the tube that is just outside the eyepiece tube. (You need to use the small screwdriver to loosen the screw on the side, which is recessed inside a hole. The hole may have a dime-sized cover over it.)

The white plastic pedestal should have already been adjusted properly to get the diffraction grating oriented correctly in three dimensions, but you should check it carefully. There are some clever features built into the apparatus to help in accomplishing this. As shown in the third figure, there are three axes about which the grating could be rotated. Rotation about axis 1 is like opening a door, and this is accomplished by rotating the entire pedestal like a lazy Susan. Rotation about axes 2 (like folding down a tailgate) and 3 are accomplished using the tripod of screws underneath the pedestal. The eye-


piece of the telescope is of a type called a Gauss eyepiece, with a diagonal piece of glass in it. When the grating is oriented correctly about axes 1 and 2 and the telescope is at  $\theta = 0$ , a beam of light that enters through the side of the eyepiece is partially reflected to the grating, and then reflected from the grating back to the eye. If these two axes are correctly adjusted, the reflected image of the crosshairs is superimposed on the crosshairs.

First get a rough initial adjustment of the pedestal by moving the telescope to 90 degrees and sighting along it like a gun to line up the grating. Now loosen the screw (not shown in the diagram) that frees the rotation of the pedestal. Put a desk lamp behind the slits of the collimator, line up the telescope with the  $m = 0$  image (which may not be exactly at  $\alpha = 180$  degrees), remove the desk lamp, cover the whole apparatus with the black cloth, and position a penlight so that it shines in through the hole in the side of the eyepiece. Adjust axes 1 and 2. If you're far out of adjustment, you may see part of a circle of light, which is the reflection of the penlight; start by bringing the circle of light into your field of view. When you're done, tighten the screw that keeps the pedestal from rotating. The pedestal is locked down to the tripod screws by the tension in a spring, which keeps the tips of two of the screws secure in dimples underneath the platform. Don't lower the screws too much, or the pedestal will no longer stay locked; make a habit of gently wiggling the pedestal after each adjustment to make sure it's not floating loose. Two of the spectrometers have the diagonal missing from their eyepieces, so if you have one of those, you'll have to borrow an eyepiece from another group to do this adjustment.


For the adjustment of axis 3, place a piece of masking tape so that it covers exactly half of the slits of the collimator. Put the Hg discharge tube behind the slits. The crosshairs should be near the edge of the tape in the  $m = 0$  image. Move the telescope out to a large angle where you see one of the high- $m$  Hg lines, and adjust the tripod screws so that the crosshairs are at the same height relative to the edge of the tape.

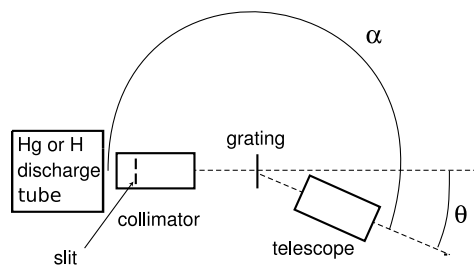

The spectrometer

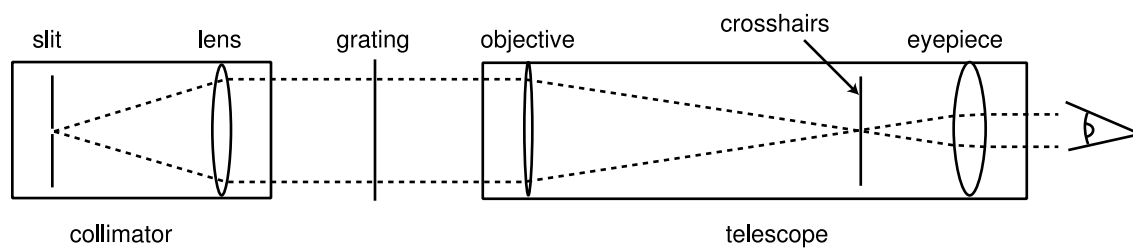

Optics.

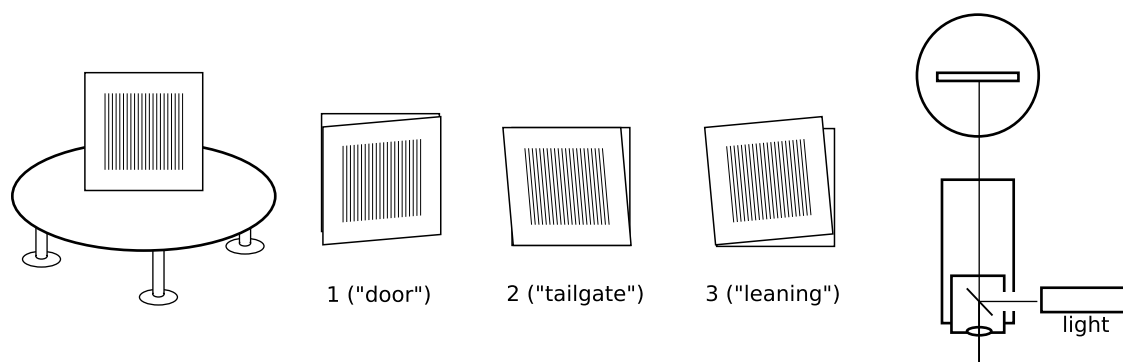

Orienting the grating.



# 16b The Mass of the Electron

## Apparatus

H gas discharge tube	3
Hg gas discharge tube (PASCO OS-9286)	3
spectrometer	1/group
diffraction grating, 600 lines/mm	1/group
small screwdriver	1
black cloth	1
piece of plywood	1
block of wood	1
light block	1

## Goals

The lab is split up into parts a, b, and c. Student lab groups will do either part b or part c. This part, b, is a measurement of the mass of the electron.

## Introduction

What's going on inside an atom? The question would have seemed nonsensical to physicists before the 20th century — the word “atom” is Greek for “unsplit-table,” and there was no evidence for subatomic particles. Only after Thomson and Rutherford had demonstrated the existence of electrons and the nucleus did the atom begin to be imagined as a tiny solar system, with the electrons moving in elliptical orbits around the nucleus under the influence of its electric field. The problem was that physicists knew very well that accelerating charges emit electromagnetic radiation, as for example in a radio antenna, so the acceleration of the electrons should have caused them to emit light, steadily lose energy, and spiral into the nucleus, all within a microsecond.

Luckily for us, atoms do not spontaneously shrink down to nothing, but there was indeed evidence that atoms could emit light. The spectra emitted by very hot gases were observed to consist of patterns of discrete lines, each with a specific wavelength. The process of emitting light always seemed to stop short of finally annihilating the atom — why? Also, why were only those specific wavelengths emitted?

The first step toward understanding the structure of the atom was Einstein's theory that light consisted of particles (photons), whose energy was related to their frequency by the equation  $E_{\text{photon}} = hf$ , or

substituting  $f = c/\lambda$ ,  $E_{\text{photon}} = hc/\lambda$ .

According to this theory, the discrete wavelengths that had been observed came from photons with specific energies. It seemed that the atom could exist only in specific states of specific energies. To get from an initial state with energy  $E_i$  to a final state with a lower energy  $E_f$ , conservation of energy required the atom to release a photon with an energy of  $E_{\text{photon}} = E_i - E_f$ .

Not only could the discrete line spectra be explained, but if the atom possessed a state of least energy (called a “ground state”), then it would always end up in that state, and it could not collapse entirely. Knowing the differences between the energy levels of the atom, it was then possible to work backwards and figure out the atomic energy levels themselves.

## Method

The apparatus you will use to observe the spectrum of hydrogen or nitrogen is shown in the first figure below. For a given wavelength, the grating produces diffracted light at many different angles: a central zeroth-order line at  $\theta = 0$ , first-order lines on both the left and right, and so on through higher-order lines at larger angles. The line of order  $m$  occurs at an angle satisfying the equation  $m\lambda = d \sin \theta$ .

To measure a wavelength, you will move the telescope until the diffracted first-order image of the slit is lined up with the telescope's cross-hairs and then read off the angle. Note that the angular scale on the table of the spectroscope actually gives the angle labeled  $\alpha$  in the figure, not  $\theta$ .

### Eliminating systematic errors

A trick to eliminate the error due to misalignment of the angular scale is to observe the same line on both the right and the left, and take  $\theta$  to be half the difference between the two angles, i.e.,  $\theta = (\alpha_R - \alpha_L)/2$ . Because you are subtracting two angles, any source of error that adds a constant offset onto the angles is eliminated. A few of the spectrometers have their angular scales out of alignment with the collimators by as much as a full degree, but that's of absolutely no consequence if this technique is used.

Regarding the calibration of  $d$ , the first person who ever did this type of experiment simply had to make

a diffraction grating whose  $d$  was very precisely constructed. But once someone has accurately measured at least one wavelength of one emission line of one element, one can simply determine the spacing,  $d$ , of any grating using a line whose wavelength is known.

## Observations

Turn on the mercury discharge tube right away, to let it get warmed up.

The second figure below shows the optics from the side, with the telescope simply looking down the throat of the collimator at  $\theta = 0$ . You are actually using the optics to let you see an image of the slit, not the tube itself. The point of using a telescope is that it provides angular magnification, so that a small change in angle can be seen visually.

A lens is used inside the collimator to make the light from the slit into a parallel beam. This is important, because we are using  $m\lambda = d\sin\theta$  to determine the wavelength, but this equation was derived under the assumption that the light was coming in as a parallel beam. To make a parallel beam, the slit must be located accurately at the focal point of the lens. This adjustment should have already been done, but you will check later and make sure. A further advantage of using a lens in the collimator is that a telescope only works for objects far away, not nearby objects from which the reflected light is diverging strongly. The lens in the collimator forms a virtual image at infinity, on which the telescope can work.

The objective lens of the telescope focuses the light, forming a real image inside the tube. The eyepiece then acts like a magnifying glass to let you see the image. In order to see the cross-hairs and the image of the slit both in focus at the same time, the cross-hairs must be located accurately at the focal point of the objective, right on top of the image.

## Setup

Skim lab 16a so you have some idea of the way the apparatus has been carefully aligned in advance by the instructor or lab technician.

### A Calibration

You will use the blue line from mercury as a calibration. In theory it shouldn't matter what known line we use for calibration, but in practice there may be small aberrations in the spectrometer, and their effect is minimized by using calibration lines of nearly

the same wavelength as the unknown lines to be measured.

Put the mercury tube behind the collimator. Make sure the hottest part of the tube is directly in front of the slits. You will need to use pieces of wood to get the height right. You want the tube as close to the slits as possible, and lined up with the slits as well as possible; you can adjust this while looking through the telescope at an  $m = 1$  line, so as to make the line as bright as possible.

If your optics are adjusted correctly, you should be able to see the microscopic bumps and scratches on the knife edges of the collimator, and there should be no parallax of the crosshairs relative to the image of the slits.

Here is a list of the wavelengths of the most prominent visible Hg lines, in nm, to high precision.<sup>1</sup>

Mercury:

404.656	violet	There is a dimmer violet line nearby at 407.781 nm.
435.833	blue	
491.604	blue-green	Dim. You may also see another blue-green line that is even dimmer.
546.074	green yellow	This is actually a complex set of lines, so it's not useful for calibration.

Start by making sure that you can find all of the lines in the correct sequence — if not, then you have probably found some first-order lines and some second-order ones. If you can find some lines but not others, use your head and search for them in the right area based on where you found the lines you did see. You may see various dim, fuzzy lights through the telescope — don't waste time chasing these, which could be coming from other tubes or from reflections. The real lines will be bright, clear and well-defined. By draping the black cloth over the discharge tube and the collimator, you can get rid of stray light that could cause problems for you or others. The discharge tubes also have holes in the back; to block the stray light from these holes, either put the two discharge tubes back to back or use one of the small "light blocks" that slide over the hole.

We will use the wavelength  $\lambda_c$  of the blue Hg line as a calibration. Measure its two angles  $\alpha_L$  and  $\alpha_R$ , and

<sup>1</sup>The table gives the wavelengths in vacuum. Although we're doing the lab in air, our goal is to find what the hydrogen or nitrogen wavelengths would have been in vacuum; by calibrating using vacuum wavelengths for mercury, we end up getting vacuum wavelengths for our unknowns as well.

check that the resulting value of  $\theta_c$  is close to the approximate ones predicted in prelab question P1. The nominal value of the spacing of the grating given in that prelab question is not very accurate. Having measured  $\theta_c$ , then we can sidestep the determination of the grating's spacing entirely and determine an unknown wavelength  $\lambda$  by using the relation

$$\lambda = \frac{\sin \theta}{\sin \theta_c} \lambda_c \quad .$$

The angles are measured using a vernier scale, which is similar to the one on the vernier calipers you have already used in the first-semester lab course. Your final reading for an angle will consist of degrees plus minutes. (One minute of arc, abbreviated 1', is 1/60 of a degree.) The main scale is marked every 30 minutes. Your initial, rough reading is obtained by noting where the zero of the vernier scale falls on the main scale, and is of the form "xxx °0' plus a little more" or "xxx °30' plus a little more." Next, you should note which line on the vernier scale lines up most closely with one of the lines on the main scale. The corresponding number on the vernier scale tells you how many minutes of arc to add for the "plus a little more."

As a check on your results, everybody in your group should take independent readings of every angle you measure in the lab, nudging the telescope to the side after each reading. Once you have independent results for a particular angle, compare them. If they're consistent to within one or two minutes of arc, average them. If they're not consistent, figure out what went wrong.

## B Spectroscopy of Hydrogen

You will study the spectrum of light emitted by the hydrogen atom, the simplest of all atoms, with just one proton and one electron. In 1885, before electrons and protons had even been imagined, a Swiss schoolteacher named Johann Balmer discovered that the wavelengths emitted by hydrogen were related by mysterious ratios of small integers. For instance, the wavelengths of the red line and the blue-green line form a ratio of exactly 20/27. Balmer even found a mathematical rule that gave all the wavelengths of the hydrogen spectrum (both the visible ones and the invisible ones that lay in the infrared and ultraviolet). The formula was completely empirical, with no theoretical basis, but clearly there were patterns lurking in the seemingly mysterious atomic spectra.

Niels Bohr showed that the energy levels of hydrogen obey a relatively simple equation,

$$E_n = -\frac{mk^2e^4}{2\hbar^2} \cdot \frac{1}{n^2}$$

where  $n$  is an integer labeling the level,  $k$  is the Coulomb constant,  $e$  is the fundamental unit of charge,  $\hbar$  is Planck's constant over  $2\pi$ , and  $m$  is the mass of the electron. All the energies of the photons in the emission spectrum could now be explained as differences in energy between specific states of the atom. For instance the four visible wavelengths observed by Balmer all came from cases where the atom ended up in the  $n = 2$  state, dropping down from the  $n = 3, 4, 5$ , and  $6$  states.

Although the equation's sheer size may appear formidable, keep in mind that the quantity  $mk^2e^4/2\hbar^2$  in front is just a numerical constant, and the variation of energy from one level to the next is of the very simple mathematical form  $1/n^2$ . It was because of this basic simplicity that the wavelength ratios like 20/27 occurred. The minus sign occurs because the equation includes both the electron's potential energy and its kinetic energy, and the standard choice of a reference-level for the potential energy results in negative values.

Now try swapping in the hydrogen tube in place of the mercury tube, and go through a similar process of acquainting yourself with the four lines in its visible spectrum, which are as follows:

violet	dim
purple	
blue-green	
red	

Again you'll again have to make sure the hottest part of the tube is in front of the collimator; this requires putting books and/or blocks of wood under the discharge tube.


We will use the purple and blue-green lines to determine the mass of the electron.

## Prelab

The point of the prelab questions is to make sure you understand what you're doing, why you're doing it, and how to avoid some common mistakes. If you don't know the answers, make sure to come to my office hours before lab and get help! Otherwise you're just setting yourself up for failure in lab.

**P1** The nominal (and not very accurate) spacing of the grating is stated as 600 lines per millimeter. From this information, find  $d$ , and predict the angles  $\alpha_L$  and  $\alpha_R$  at which you will observe the blue mercury line.

**P2** Make sure you understand the first three vernier

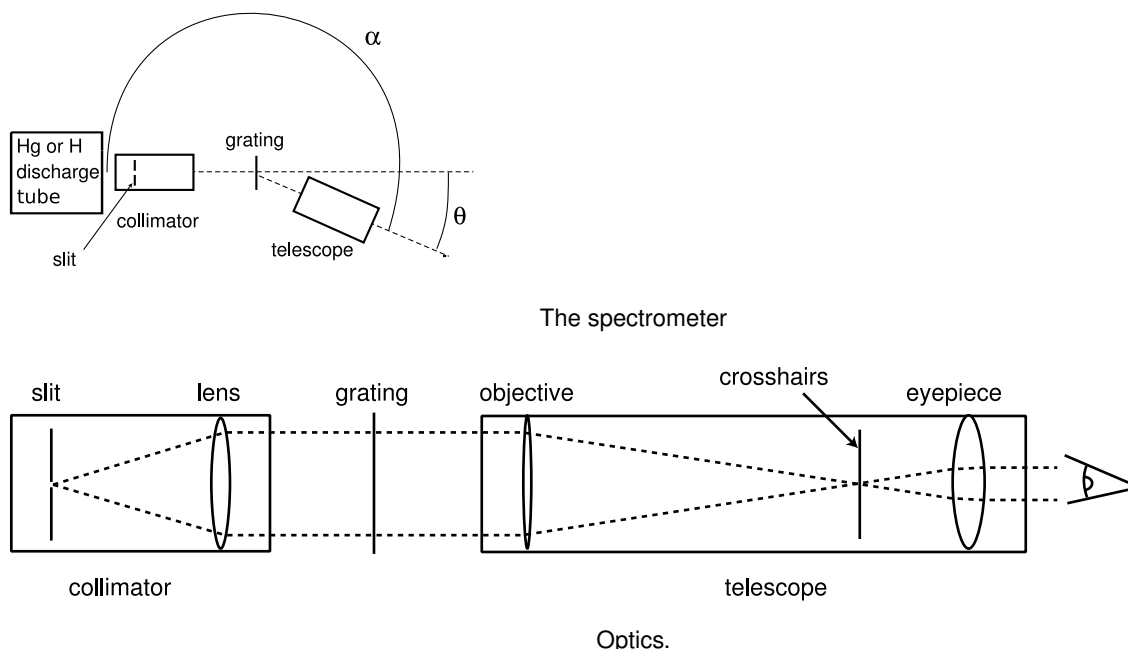

readings in the fourth figure, and then interpret the fourth reading.

**P3** For the calibration with mercury, in what sequence do you expect to see the lines on each side? Make a drawing showing the sequence of the angles as you go out from  $\theta=0$ .

**P4** The visible lines of hydrogen come from the  $3 \rightarrow 2$ ,  $4 \rightarrow 2$ ,  $5 \rightarrow 2$ , and  $6 \rightarrow 2$  transitions. Based on  $E = hf$ , which of these should correspond to which colors?

known:

$$e = 1.6022 \times 10^{-19} \text{ C}$$

$$k = 8.9876 \times 10^9 \text{ N}\cdot\text{m}^2/\text{C}^2$$

$$h = 6.6261 \times 10^{-34} \text{ J}\cdot\text{s}$$

$$c = 2.9979 \times 10^8 \text{ m/s}$$

The energies of the four types of visible photons emitted by a hydrogen atom equal  $E_n - E_2$ , where  $n = 3, 4, 5$ , and  $6$ . Using the Bohr equation, we have

$$E_{\text{photon}} = A \left( \frac{1}{4} - \frac{1}{n^2} \right),$$

where  $A$  is the expression from the Bohr equation that depends on the mass of the electron. From the two lines you've measured, extract a value for  $A$ . If your data passed the self-check above, then you should find that these values for  $A$  agree to no worse than a few parts per thousand at worst. Compute an average value of  $A$ , and extract the mass of the electron, with error bars.

Finally, there is a small correction that should be made to the result for the mass of the electron because actually the proton isn't infinitely massive compared to the electron; in terms of the quantity  $m$  given by the equation on page 70, the mass of the electron,  $m_e$ , would actually be given by  $m_e = m/(1 - m/m_p)$ , where  $m_p$  is the mass of the proton,  $1.6726 \times 10^{-27} \text{ kg}$ .

## Self-Check


In homework problem 36-11 in *Light and Matter*, you calculated the ratio

$\lambda_{\text{blue-green}}/\lambda_{\text{purple}}$ . Before leaving lab, make sure that your wavelengths are consistent with this prediction, to a precision of no worse than about one part per thousand.

## Analysis

Throughout your analysis, remember that this is a high-precision experiment, so you don't want to round off to less than five significant figures.

We assume that the following constants are already

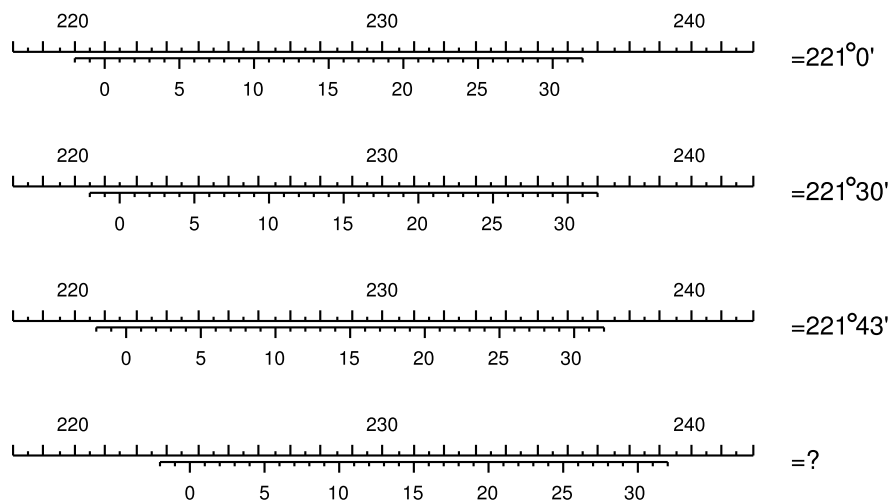

Prelab question 2.





# 16c The Nitrogen Molecule

## Apparatus

He gas discharge tube .....	3
N <sub>2</sub> gas discharge tube (in green carousel) .....	3
spectrometer .....	1/group
diffraction grating, 600 lines/mm .....	1/group
small screwdriver .....	1
black cloth .....	1
piece of plywood .....	1
block of wood .....	1
penlight .....	1/group
light block	

## Goals

The lab has three parts. Each group will only do two parts. In this lab, part c, you will use an energy sum to test a hypothesis about the energy levels of the nitrogen molecule, N<sub>2</sub>.

## Method

The apparatus you will use to observe the spectrum of hydrogen or nitrogen is shown in the first figure below. For a given wavelength, the grating produces diffracted light at many different angles: a central zeroth-order line at  $\theta = 0$ , first-order lines on both the left and right, and so on through higher-order lines at larger angles. The line of order  $m$  occurs at an angle satisfying the equation  $m\lambda = d \sin \theta$ .

To measure a wavelength, you will move the telescope until the diffracted first-order image of the slit is lined up with the telescope's cross-hairs and then read off the angle. Note that the angular scale on the table of the spectroscope actually gives the angle labeled  $\alpha$  in the figure, not  $\theta$ .

### Eliminating systematic errors


A trick to eliminate the error due to misalignment of the angular scale is to observe the same line on both the right and the left, and take  $\theta$  to be half the difference between the two angles, i.e.,  $\theta = (\alpha_R - \alpha_L)/2$ . Because you are subtracting two angles, any source of error that adds a constant offset onto the angles is eliminated. A few of the spectrometers have their angular scales out of alignment with the collimators by as much as a full degree, but that's of absolutely

no consequence if this technique is used.

Regarding the calibration of  $d$ , the first person who ever did this type of experiment simply had to make a diffraction grating whose  $d$  was very precisely constructed. But once someone has accurately measured at least one wavelength of one emission line of one element, one can simply determine the spacing,  $d$ , of any grating using a line whose wavelength is known.

## B Energy Sums

The nitrogen discharge tube is housed in a green plastic carousel. With the power off, rotate the carousel so that the nitrogen tube is the one that is in the active position, and then turn on the power. If you hold a diffraction grating up to your eye and look at the tube, you will see a remarkable spectrum, unlike the visible light spectrum of almost any other gas. This is because the N<sub>2</sub> molecule has an extremely strong bond, requiring twice the energy to break compared to otherwise similar gases such as H<sub>2</sub> or O<sub>2</sub>. Whereas these other gases would break up into individual atoms under the extreme conditions present in a discharge tube, the nitrogen molecule holds together, so that you are seeing the spectrum of the molecule, not the atom. For this reason, the spectrum of nitrogen contains a large number of lines.

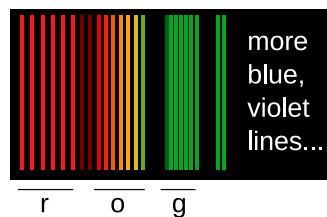


But these lines are not random. They occur in sets, each of which looks like a comb with an approximately equal spacing between the “teeth.” The figure shows a portion of the spectrum, including three sets of lines, which I have labeled r (red), o (orange to green), and g (green).

I have spent some time trying to interpret the origin of these lines, and I believe the interpretation is something like this. Each of these lines is the emission of a photon as the molecule goes from an initial state to a final state that has less energy. The initial state has some energy because the electrons are in an excited state (labeled B by spectroscopists)

and also some energy because the molecule is vibrating, like two masses connected by a spring. The final state has the electrons in a lower-energy state (labeled A), and is also vibrating. The initial and final electronic states B and A are the same in all cases, but the vibrational states differ. An idealized quantum-mechanical vibrator turns out to have a series of energy states like a ladder with nearly evenly spaced rungs. States higher on the “ladder” are vibrating more violently — classically, they vibrate with greater amplitude. The rungs of the vibrational ladder are labeled  $v = 0, 1, 2$ , and so on. (Because of the Heisenberg uncertainty principle, some vibrational energy is present even in the  $v = 0$  state.) I think the states in the red set are from a state  $v$  to a state  $v - 3$ , i.e., a change of 3 units in the vibrational quantum number. The o set would be a change of 4, and g a change of 5.

This hypothesis can be tested as follows. If it is true, we could pick out a set of energy levels like the following example:

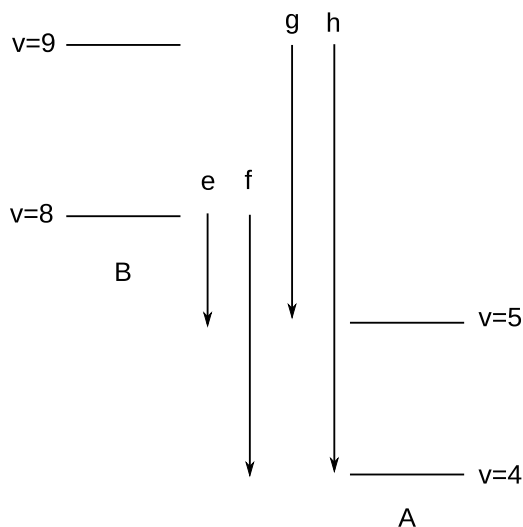

The letters  $e$ ,  $f$ ,  $g$ , and  $h$  are the energy differences that would be observed as the energies of the photons. In a set like this, we would have

$$h - g = f - e,$$

since each side of the equation would be equal to the energy difference between the  $v = 4$  and  $5$  states of ladder A. Try to find a set of lines that would be consistent with this interpretation. This may require some trial and error, but I think it may work if  $e$  is one of the lines near the middle of the r set,  $g$  is the orange line that is fourth from the short-wavelength end of the o set,  $f$  is the fifth in that set, and  $h$  is in the g set.

## C Calibration

You will use the yellow line from helium as a calibration. In theory it shouldn't matter what known line we use for calibration, but in practice there may be small aberrations in the spectrometer, and their effect is minimized by using calibration lines of nearly the same wavelength as the unknown lines to be measured.

Put the helium tube behind the collimator. Make sure the hottest part of the tube is directly in front of the slits. You will need to use pieces of wood to get the height right. You want the tube as close to the slits as possible, and lined up with the slits as well as possible; you can adjust this while looking through the telescope at an  $m = 1$  line, so as to make the line as bright as possible.

If your optics are adjusted correctly, you should be able to see the microscopic bumps and scratches on the knife edges of the collimator, and there should be no parallax of the crosshairs relative to the image of the slits.

Here is a list of the wavelengths of the most prominent visible He lines, in nm, to high precision.<sup>1</sup>

Helium:	
447.148	bright blue-purple
471.314	dim blue
492.193	dim green
501.567	bright green
587.562	yellow
667.815	dim red
706.5	very dim red

Start by making sure that you can find all of the lines in the correct sequence — if not, then you have probably found some first-order lines and some second-order ones. If you can find some lines but not others, use your head and search for them in the right area based on where you found the lines you did see. You may see various dim, fuzzy lights through the telescope — don't waste time chasing these, which could be coming from other tubes or from reflections. The real lines will be bright, clear and well-defined. By draping the black cloth over the discharge tube and the collimator, you can get rid of stray light that could cause problems for you or others. The discharge tubes also have holes in the back; to block the stray light from these holes, either put the two discharge tubes back to back or use one of the small

<sup>1</sup>The table gives the wavelengths in vacuum. Although we're doing the lab in air, our goal is to find what the nitrogen wavelengths would have been in vacuum; by calibrating using vacuum wavelengths for mercury, we end up getting vacuum wavelengths for our unknowns as well.

“light blocks” that slide over the hole.

We will use the wavelength  $\lambda_c$  of the hollow He line as a calibration. Measure its two angles  $\alpha_L$  and  $\alpha_R$ , and check that the resulting value of  $\theta_c$  is close to the approximate ones predicted in prelab question P1. The nominal value of the spacing of the grating given in that prelab question is not very accurate. Having measured  $\theta_c$ , then we can sidestep the determination of the grating’s spacing entirely and determine an unknown wavelength  $\lambda$  by using the relation

$$\lambda = \frac{\sin \theta}{\sin \theta_c} \lambda_c \quad .$$

The angles are measured using a vernier scale, which is similar to the one on the vernier calipers you have already used in the first-semester lab course. Your final reading for an angle will consist of degrees plus minutes. (One minute of arc, abbreviated 1', is 1/60 of a degree.) The main scale is marked every 30 minutes. Your initial, rough reading is obtained by noting where the zero of the vernier scale falls on the main scale, and is of the form “xxx°0' plus a little more” or “xxx°30' plus a little more.” Next, you should note which line on the vernier scale lines up most closely with one of the lines on the main scale. The corresponding number on the vernier scale tells you how many minutes of arc to add for the “plus a little more.”

As a check on your results, everybody in your group should take independent readings of every angle you measure in the lab, nudging the telescope to the side after each reading. Once you have independent results for a particular angle, compare them. If they’re consistent to within one or two minutes of arc, average them. If they’re not consistent, figure out what went wrong.

## Status as of August 2018

In spring of 2018, my students and I worked on measuring and interpreting this spectrum. A summary of our results is in a Google Docs spreadsheet at [goo.gl/akrbcY](https://goo.gl/akrbcY). There is a basic explanation of the physics in Simple Nature section 14.2. More detailed information about my interpretation of the lines is at [physics.stackexchange.com/a/334451/4552](https://physics.stackexchange.com/a/334451/4552).

In the notation used in the material on stackexchange, the states are labeled with a quantum number  $v$ . The energy sum based on our data come out quite nice for  $v = 9$  and 8 going to 5 and 4, accurate to about 0.3%, which is reasonable for this

technique. The energy sum for  $v = 10$  and 9 going to 6 and 5 also works well, but with slightly higher error, maybe partly because 10-5 is a doublet. There are two dim green lines that in my labeling system would be g0 and g-1. These may be the ones we would need in order to get a couple more energy sums.

I experimented with doing the measurements photographically. A student took a photo of a diffraction pattern using a cell phone. I first did a rough calibration against student data, then improved this calibration by doing a linear fit to my own spectrometer data. This worked fairly well, but doesn’t work without actual spectrometer data, which are needed for calibration.

We have digital spectrometers which may be helpful here and are worth trying. Their resolution is supposed to be about 1 or 1.5 nm, which is an order of magnitude worse than the analog spectrometers, but may be adequate for this purpose, and they can display a spectrum as a graph, which may be help enough to make up for the lower resolution.

I couldn’t resolve the green band. I asked a couple of students the next day, and they seemed to think that it was doable to resolve these lines. Possibly my tube was behaving differently than theirs, but I’m not sure I believe their data. The red and orange bands came out nice, all wavelengths being within a few tenths of a nm of Lofthus’s values.

## Prelab


The point of the prelab questions is to make sure you understand what you’re doing, why you’re doing it, and how to avoid some common mistakes. If you don’t know the answers, make sure to come to my office hours before lab and get help! Otherwise you’re just setting yourself up for failure in lab.

The week before you are to do the lab, briefly familiarize yourself visually with the apparatus.


**P1** The nominal (and not very accurate) spacing of the grating is stated as 600 lines per millimeter. From this information, find  $d$ , and predict the angles  $\alpha_L$  and  $\alpha_R$  at which you will observe the yellow helium line.

**P2** Make sure you understand the first three vernier readings in the fourth figure, and then interpret the fourth reading.

**P3** For the calibration with helium, in what sequence do you expect to see the lines on each side?

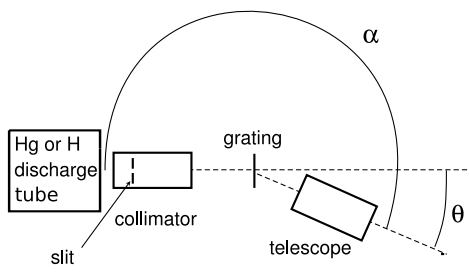


The spectrometer

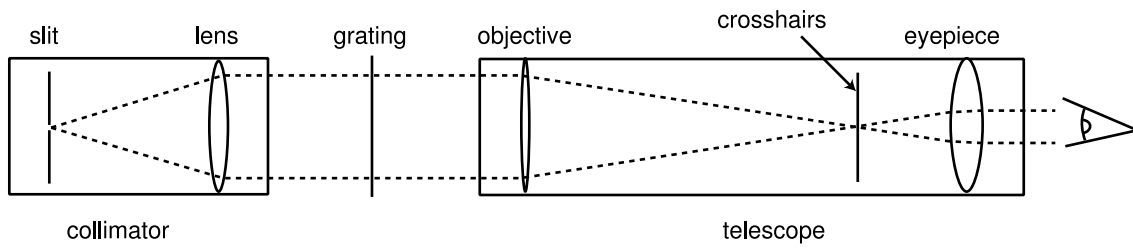


Optics.

Make a drawing showing the sequence of the angles as you go out from  $\theta=0$ .

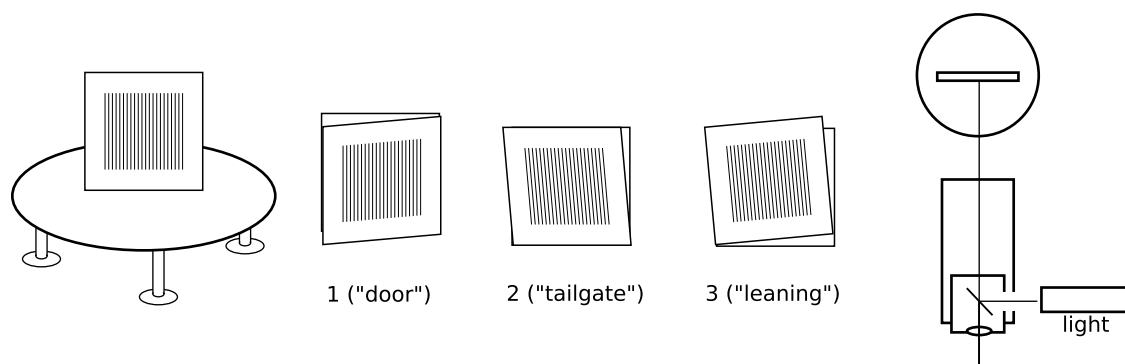

Orienting the grating.

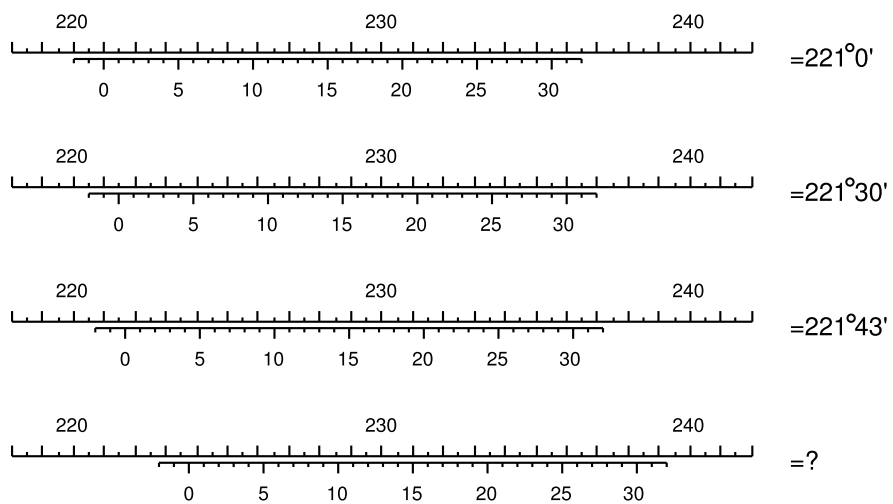

Prelab question 2.



# Appendix 1: Format of Lab Writeups

Lab reports must be three pages or less, not counting your raw data. The format should be as follows:

## Title

**Raw data** — *Keep actual observations separate from what you later did with them.*

These are the results of the measurements you take down during the lab, hence they come first. Write your raw data directly in your lab book; don't write them on scratch paper and recopy them later. Don't use pencil. The point is to separate facts from opinions, observations from inferences.

**Procedure** — *Did you have to create your own methods for getting some of the raw data?*

Do not copy down the procedure from the manual. In this section, you only need to explain any methods you had to come up with on your own, or cases where the methods suggested in the handout didn't work and you had to do something different. Don't write anything here unless you think I will really care and want to change how we do the lab in the future. In most cases this section can be totally blank. Do not discuss how you did your calculations here, just how you got your raw data.

**Abstract** — *What did you find out? Why is it important?*

The "abstract" of a scientific paper is a *short* paragraph at the top that summarizes the experiment's results in a few sentences.

Many of our labs are comparisons of theory and experiment. The abstract for such a lab needs to say whether you think the experiment was consistent with theory, or not consistent with theory. If your results deviated from the ideal equations, don't be afraid to say so. After all, this is real life, and many of the equations we learn are only approximations, or are only valid in certain circumstances. However, (1) if you simply mess up, it is your responsibility to realize it in lab and do it again, right; (2) you will never get exact agreement with theory, because measurements are not perfectly exact — the important issue is whether your results agree with theory to roughly within the error bars.

The abstract is not a statement of what you hoped to find out. It's a statement of what you *did* find out. It's like the brief statement at the beginning of a debate: "The U.S. should have free trade with China." It's not this: "In this debate, we will discuss

whether the U.S. should have free trade with China."

If this is a lab that has just one important numerical result (or maybe two or three of them), put them in your abstract, with error bars where appropriate. There should normally be no more than two to four numbers here. Do not recapitulate your raw data here — this is for your final results.

If you're presenting a final result with error bars, make sure that the number of significant figures is consistent with your error bars. For example, if you write a result as  $323.54 \pm 6$  m/s, that's wrong. Your error bars say that you could be off by 6 in the ones' place, so the 5 in the tenths' place and the four in the hundredths' place are completely meaningless.

If you're presenting a number in scientific notation, with error bars, don't do it like this

$$1.234 \times 10^{-89} \text{ m/s} \pm 3 \times 10^{-92} \text{ m/s} \quad ,$$

do it like this

$$(1.234 \pm 0.003) \times 10^{-89} \text{ m/s} \quad ,$$

so that we can see easily which digit of the result the error bars apply to.

**Calculations and Reasoning** — *Convince me of what you claimed in your abstract.*

Often this section consists of nothing more than the calculations that you started during lab. If those calculations are clear enough to understand, and there is nothing else of interest to explain, then it is not necessary to write up a separate narrative of your analysis here. If you have a long series of similar calculations, you may just show one as a sample. If your prelab involved deriving equations that you will need, repeat them here without the derivation.

In some labs, you will need to go into some detail here by giving logical arguments to convince me that the statements you made in the abstract follow logically from your data. Continuing the debate metaphor, if your abstract said the U.S. should have free trade with China, this is the rest of the debate, where you convince me, based on data and logic, that we should have free trade.





## Appendix 2: Basic Error Analysis

### No measurement is perfectly exact.

One of the most common misconceptions about science is that science is “exact.” It is always a struggle to get beginning science students to believe that no measurement is perfectly correct. They tend to think that if a measurement is a little off from the “true” result, it must be because of a mistake — if a pro had done it, it would have been right on the mark. Not true!

What scientists can do is to estimate just how far off they might be. This type of estimate is called an error bar, and is expressed with the  $\pm$  symbol, read “plus or minus.” For instance, if I measure my dog’s weight to be  $52 \pm 2$  pounds, I am saying that my best estimate of the weight is 52 pounds, and I think I could be off by roughly 2 pounds either way. The term “error bar” comes from the conventional way of representing this range of uncertainty of a measurement on a graph, but the term is also used when no graph is involved.

Some very good scientific work results in measurements that nevertheless have large error bars. For instance, the best measurement of the age of the universe is now  $15 \pm 5$  billion years. That may not seem like wonderful precision, but the people who did the measurement knew what they were doing. It’s just that the only available techniques for determining the age of the universe are inherently poor.

Even when the techniques for measurement are very precise, there are still error bars. For instance, electrons act like little magnets, and the strength of a very weak magnet such as an individual electron is customarily measured in units called Bohr magnetons. Even though the magnetic strength of an electron is one of the most precisely measured quantities ever, the best experimental value still has error bars:  $1.0011596524 \pm 0.0000000002$  Bohr magnetons.

There are several reasons why it is important in scientific work to come up with a numerical estimate of your error bars. If the point of your experiment is to test whether the result comes out as predicted by a theory, you know there will always be some disagreement, even if the theory is absolutely right. You need to know whether the measurement is reasonably consistent with the theory, or whether the discrepancy is too great to be explained by the lim-

itations of the measuring devices.

Another important reason for stating results with error bars is that other people may use your measurement for purposes you could not have anticipated. If they are to use your result intelligently, they need to have some idea of how accurate it was.

### Error bars are not absolute limits.

Error bars are not absolute limits. The true value may lie outside the error bars. If I got a better scale I might find that the dog’s weight is  $51.3 \pm 0.1$  pounds, inside my original error bars, but it’s also possible that the better result would be  $48.7 \pm 0.1$  pounds. Since there’s always some chance of being off by a somewhat more than your error bars, or even a lot more than your error bars, there is no point in being extremely conservative in an effort to make absolutely sure the true value lies within your stated range. When a scientist states a measurement with error bars, she is not saying “If the true value is outside this range, I deserve to be drummed out of the profession.” If that was the case, then every scientist would give ridiculously inflated error bars to avoid having her career ended by one fluke out of hundreds of published results. What scientists are communicating to each other with error bars is a typical amount by which they might be off, not an upper limit.

The important thing is therefore to define error bars in a standard way, so that different people’s statements can be compared on the same footing. By convention, it is usually assumed that people estimate their error bars so that about two times out of three, their range will include the true value (or the results of a later, more accurate measurement with an improved technique).

### Random and systematic errors.

Suppose you measure the length of a sofa with a tape measure as well as you can, reading it off to the nearest millimeter. If you repeat the measurement again, you will get a different answer. (This is assuming that you don’t allow yourself to be psychologically biased to repeat your previous answer, and that 1 mm is about the limit of how well you can see.) If you kept on repeating the measurement,

you might get a list of values that looked like this:

203.1 cm	203.4	202.8	203.3	203.2
203.4	203.1	202.9	202.9	203.1

Variations of this type are called random errors, because the result is different every time you do the measurement.

The effects of random errors can be minimized by averaging together many measurements. Some of the measurements included in the average are too high, and some are too low, so the average tends to be better than any individual measurement. The more measurements you average in, the more precise the average is. The average of the above measurements is 203.1 cm. Averaging together many measurements cannot completely eliminate the random errors, but it can reduce them.


On the other hand, what if the tape measure was a little bit stretched out, so that your measurements always tended to come out too low by 0.3 cm? That would be an example of a systematic error. Since the systematic error is the same every time, averaging didn't help us to get rid of it. You probably had no easy way of finding out exactly the amount of stretching, so you just had to suspect that there might a systematic error due to stretching of the tape measure.

Some scientific writers make a distinction between the terms “accuracy” and “precision.” A precise measurement is one with small random errors, while an accurate measurement is one that is actually close to the true result, having both small random errors and small systematic errors. Personally, I find the distinction is made more clearly with the more memorable terms “random error” and “systematic error.”


The  $\pm$  sign used with error bars normally implies that random errors are being referred to, since random errors could be either positive or negative, whereas systematic errors would always be in the same direction.

## The goal of error analysis


Very seldom does the final result of an experiment come directly off of a clock, ruler, gauge or meter. It is much more common to have raw data consisting of direct measurements, and then calculations based on the raw data that lead to a final result. As an example, if you want to measure your car's gas mileage, your raw data would be the number of gallons of gas consumed and the number of miles you went. You would then do a calculation, dividing

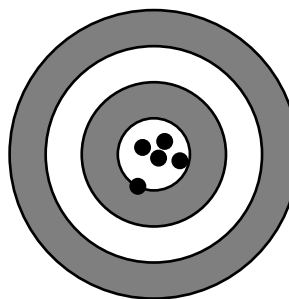

small random errors,  
small systematic error

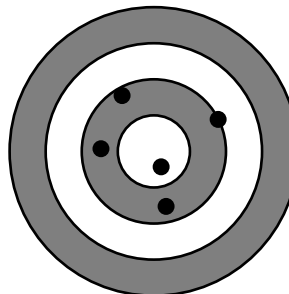

large random errors,  
small systematic error

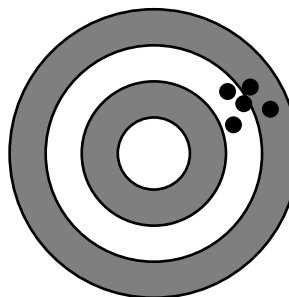

small random errors,  
large systematic error

miles by gallons, to get your final result. When you communicate your result to someone else, they are completely uninterested in how accurately you measured the number of miles and how accurately you measured the gallons. They simply want to know how accurate your final result was. Was it  $22 \pm 2$  mi/gal, or  $22.137 \pm 0.002$  mi/gal?

Of course the accuracy of the final result is ultimately based on and limited by the accuracy of your raw data. If you are off by 0.2 gallons in your measurement of the amount of gasoline, then that amount of error will have an effect on your final result. We say that the errors in the raw data “propagate” through the calculations. When you are requested to do “error analysis” in a lab writeup, that means that you

are to use the techniques explained below to determine the error bars on your final result. There are two sets of techniques you'll need to learn:

techniques for finding the accuracy of your raw data

techniques for using the error bars on your raw data to infer error bars on your final result

## Estimating random errors in raw data

We now examine three possible techniques for estimating random errors in your original measurements, illustrating them with the measurement of the length of the sofa.

### Method #1: Guess

If you're measuring the length of the sofa with a metric tape measure, then you can probably make a reasonable guess as to the precision of your measurements. Since the smallest division on the tape measure is one millimeter, and one millimeter is also near the limit of your ability to see, you know you won't be doing better than  $\pm 1$  mm, or 0.1 cm. Making allowances for errors in getting tape measure straight and so on, we might estimate our random errors to be a couple of millimeters.

Guessing is fine sometimes, but there are at least two ways that it can get you in trouble. One is that students sometimes have too much faith in a measuring device just because it looks fancy. They think that a digital balance must be perfectly accurate, since unlike a low-tech balance with sliding weights on it, it comes up with its result without any involvement by the user. That is incorrect. No measurement is perfectly accurate, and if the digital balance only displays an answer that goes down to tenths of a gram, then there is no way the random errors are any smaller than about a tenth of a gram.

Another way to mess up is to try to guess the error bars on a piece of raw data when you really don't have enough information to make an intelligent estimate. For instance, if you are measuring the range of a rifle, you might shoot it and measure how far the bullet went to the nearest centimeter, concluding that your random errors were only  $\pm 1$  cm. In reality, however, its range might vary randomly by fifty meters, depending on all kinds of random factors you don't know about. In this type of situation, you're better off using some other method of estimating your random errors.

### Method #2: Repeated Measurements and the Two-Thirds Rule

If you take repeated measurements of the same thing, then the amount of variation among the numbers can tell you how big the random errors were. This approach has an advantage over guessing your random errors, since it automatically takes into account all the sources of random error, even ones you didn't know were present.

Roughly speaking, the measurements of the length of the sofa were mostly within a few mm of the average, so that's about how big the random errors were. But let's make sure we are stating our error bars according to the convention that the true result will fall within our range of errors about two times out of three. Of course we don't know the "true" result, but if we sort out our list of measurements in order, we can get a pretty reasonable estimate of our error bars by taking half the range covered by the middle two thirds of the list. Sorting out our list of ten measurements of the sofa, we have

202.8 cm	202.9	202.9	203.1	203.1
203.1	203.2	203.3	203.4	203.4

Two thirds of ten is about 6, and the range covered by the middle six measurements is 203.3 cm - 202.9 cm, or 0.4 cm. Half that is 0.2 cm, so we'd estimate our error bars as  $\pm 0.2$  cm. The average of the measurements is 203.1 cm, so your result would be stated as  $203.1 \pm 0.2$  cm.

One common mistake when estimating random errors by repeated measurements is to round off all your measurements so that they all come out the same, and then conclude that the error bars were zero. For instance, if we'd done some overenthusiastic rounding of our measurements on the sofa, rounding them all off to the nearest cm, every single number on the list would have been 203 cm. That wouldn't mean that our random errors were zero! The same can happen with digital instruments that automatically round off for you. A digital balance might give results rounded off to the nearest tenth of a gram, and you may find that by putting the same object on the balance again and again, you always get the same answer. That doesn't mean it's perfectly precise. Its precision is no better than about  $\pm 0.1$  g.

### Method #3: Repeated Measurements and the Standard Deviation

The most widely accepted method for measuring error bars is called the standard deviation. Here's how the method works, using the sofa example again.

- (1) Take the average of the measurements.

$$\text{average} = 203.1 \text{ cm}$$

- (2) Find the difference, or “deviation,” of each measurement from the average.

-0.3 cm	-0.2	-0.2	0.0	0.0
0.0	0.1	0.1	0.3	0.3

- (3) Take the square of each deviation.

0.09 cm <sup>2</sup>	0.04	0.04	0.00	0.00
0.00	0.01	0.01	0.09	0.09

- (4) Average together all the squared deviations.

$$\text{average} = 0.04 \text{ cm}^2$$

- (5) Take the square root. This is the standard deviation.

$$\text{standard deviation} = 0.2 \text{ cm}$$

If we’re using the symbol  $x$  for the length of the couch, then the result for the length of the couch would be stated as  $x = 203.1 \pm 0.2 \text{ cm}$ , or  $x = 203.1 \text{ cm}$  and  $\sigma_x = 0.2 \text{ cm}$ . Since the Greek letter sigma ( $\sigma$ ) is used as a symbol for the standard deviation, a standard deviation is often referred to as “a sigma.”

Step (3) may seem somewhat mysterious. Why not just skip it? Well, if you just went straight from step (2) to step (4), taking a plain old average of the deviations, you would find that the average is zero! The positive and negative deviations always cancel out exactly. Of course, you could just take absolute values instead of squaring the deviations. The main advantage of doing it the way I’ve outlined above are that it is a standard method, so people will know how you got the answer. (Another advantage is that the standard deviation as I’ve described it has certain nice mathematical properties.)

A common mistake when using the standard deviation technique is to take too few measurements. For instance, someone might take only two measurements of the length of the sofa, and get 203.4 cm and 203.4 cm. They would then infer a standard deviation of zero, which would be unrealistically small because the two measurements happened to come out the same.

In the following material, I’ll use the term “standard deviation” as a synonym for “error bar,” but that does not imply that you must always use the standard deviation method rather than the guessing method or the 2/3 rule.

There is a utility on the class’s web page for calculating standard deviations.

## Probability of deviations

You can see that although 0.2 cm is a good figure for the typical size of the deviations of the measurements of the length of the sofa from the average, some of the deviations are bigger and some are smaller. Experience has shown that the following probability estimates tend to hold true for how frequently deviations of various sizes occur:


> 1 standard deviation about 1 times out of 3

> 2 standard deviations about 1 time out of 20

> 3 standard deviations about 1 in 500

> 4 standard deviations about 1 in 16,000

> 5 standard deviations about 1 in 1,700,000

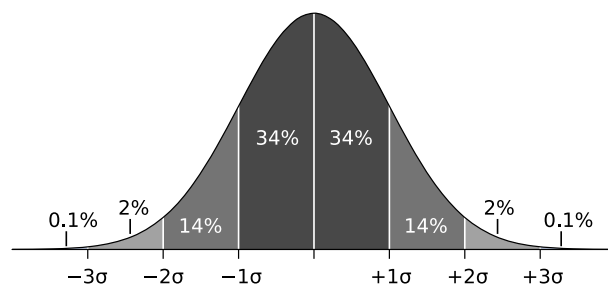

The probability of various sizes of deviations, shown graphically. Areas under the bell curve correspond to probabilities. For example, the probability that the measurement will deviate from the truth by less than one standard deviation ( $\pm 1\sigma$ ) is about  $34 \times 2 = 68\%$ , or about 2 out of 3. (J. Kemp, P. Strandmark, Wikipedia.)

### Example: How significant?

In 1999, astronomers Webb et al. claimed to have found evidence that the strength of electrical forces in the ancient universe, soon after the big bang, was slightly weaker than it is today. If correct, this would be the first example ever discovered in which the laws of physics changed over time. The difference was very small,  $5.7 \pm 1.0$  parts per million, but still highly statistically significant. Dividing, we get  $(5.7 - 0)/1.0 = 5.7$  for the number of standard deviations by which their measurement was different from the expected result of zero. Looking at the table above, we see that if the true value really was zero, the chances of this happening would be less than one in a million. In general, five standard deviations (“five sigma”) is considered the gold standard for statistical significance.

This is an example of how we test a hypothesis statistically, find a probability, and interpret the probability. The probability we find is the probability that our results

would differ this much from the hypothesis, if the hypothesis was true. It's not the probability that the hypothesis is true or false, nor is it the probability that our experiment is right or wrong.

However, there is a twist to this story that shows how statistics always have to be taken with a grain of salt. In 2004, Chand et al. redid the measurement by a more precise technique, and found that the change was  $0.6 \pm 0.6$  parts per million. This is only one standard deviation away from the expected value of 0, which should be interpreted as being statistically consistent with zero. If you measure something, and you think you know what the result is supposed to be theoretically, then one standard deviation is the amount you typically *expect* to be off by — that's why it's called the "standard" deviation. Moreover, the Chand result is wildly statistically inconsistent with the Webb result (see the example on page 89), which means that one experiment or the other is a mistake. Most likely Webb et al. underestimated their random errors, or perhaps there were systematic errors in their experiment that they didn't realize were there.

## Precision of an average

We decided that the standard deviation of our measurements of the length of the couch was 0.2 cm, i.e., the precision of each individual measurement was about 0.2 cm. But I told you that the average, 203.1 cm, was more precise than any individual measurement. How precise is the average? The answer is that the standard deviation of the average equals

$$\frac{\text{standard deviation of one measurement}}{\sqrt{\text{number of measurements}}}.$$

(An example on page 88 gives the reasoning that leads to the square root.) That means that you can theoretically measure anything to any desired precision, simply by averaging together enough measurements. In reality, no matter how small you make your random error, you can't get rid of systematic errors by averaging, so after a while it becomes pointless to take any more measurements.



## Appendix 3: Propagation of Errors

### Propagation of the error from a single variable

In the previous appendix we looked at techniques for estimating the random errors of raw data, but now we need to know how to evaluate the effects of those random errors on a final result calculated from the raw data. For instance, suppose you are given a cube made of some unknown material, and you are asked to determine its density. Density is defined as  $\rho = m/v$  ( $\rho$  is the Greek letter “rho”), and the volume of a cube with edges of length  $b$  is  $v = b^3$ , so the formula

$$\rho = m/b^3$$

will give you the density if you measure the cube’s mass and the length of its sides. Suppose you measure the mass very accurately as  $m = 1.658 \pm 0.003$  g, but you know  $b = 0.85 \pm 0.06$  cm with only two digits of precision. Your best value for  $\rho$  is  $1.658 \text{ g}/(0.85 \text{ cm})^3 = 2.7 \text{ g/cm}^3$ .

How can you figure out how precise this value for  $\rho$  is? We’ve already made sure not to keep more than two significant figures for  $\rho$ , since the less accurate piece of raw data had only two significant figures. We expect the last significant figure to be somewhat uncertain, but we don’t yet know how uncertain. A simple method for this type of situation is simply to change the raw data by one sigma, recalculate the result, and see how much of a change occurred. In this example, we add 0.06 cm to  $b$  for comparison.

$$\begin{aligned} b = 0.85 \text{ cm} & \quad \text{gave} \quad \rho = 2.7 \text{ g/cm}^3 \\ b = 0.91 \text{ cm} & \quad \text{gives} \quad \rho = 2.2 \text{ g/cm}^3 \end{aligned}$$

The resulting change in the density was  $0.5 \text{ g/cm}^3$ , so that is our estimate for how much it could have been off by:

$$\rho = 2.7 \pm 0.5 \text{ g/cm}^3$$

### Propagation of the error from several variables

What about the more general case in which no one piece of raw data is clearly the main source of error? For instance, suppose we get a more accurate measurement of the edge of the cube,  $b = 0.851 \pm 0.001$  cm. In percentage terms, the accuracies of  $m$  and

$b$  are roughly comparable, so both can cause significant errors in the density. The following more general method can be applied in such cases:

(1) Change one of the raw measurements, say  $m$ , by one standard deviation, and see by how much the final result,  $\rho$ , changes. Use the symbol  $Q_m$  for the absolute value of that change.

$$\begin{aligned} m = 1.658 \text{ g} & \quad \text{gave} \quad \rho = 2.690 \text{ g/cm}^3 \\ m = 1.661 \text{ g} & \quad \text{gives} \quad \rho = 2.695 \text{ g/cm}^3 \end{aligned}$$

$$Q_m = \text{change in } \rho = 0.005 \text{ g/cm}^3$$

(2) Repeat step (1) for the other raw measurements.

$$\begin{aligned} b = 0.851 \text{ cm} & \quad \text{gave} \quad \rho = 2.690 \text{ g/cm}^3 \\ b = 0.852 \text{ cm} & \quad \text{gives} \quad \rho = 2.681 \text{ g/cm}^3 \end{aligned}$$

$$Q_b = \text{change in } \rho = 0.009 \text{ g/cm}^3$$

(3) The error bars on  $\rho$  are given by the formula

$$\sigma_\rho = \sqrt{Q_m^2 + Q_b^2},$$

yielding  $\sigma_\rho = 0.01 \text{ g/cm}^3$ . Intuitively, the idea here is that if our result could be off by an amount  $Q_m$  because of an error in  $m$ , and by  $Q_b$  because of  $b$ , then if the two errors were in the same direction, we might be off by roughly  $|Q_m| + |Q_b|$ . However, it’s equally likely that the two errors would be in opposite directions, and at least partially cancel. The expression  $\sqrt{Q_m^2 + Q_b^2}$  gives an answer that’s smaller than  $Q_m + Q_b$ , representing the fact that the cancellation might happen.

The final result is  $\rho = 2.69 \pm 0.01 \text{ g/cm}^3$ .

#### Example: An average

On page 86 I claimed that averaging a bunch of measurements reduces the error bars by the square root of the number of measurements. We can now see that this is a special case of propagation of errors.

For example, suppose Alice measures the circumference  $c$  of a guinea pig’s waist to be 10 cm. Using the guess method, she estimates that her error bars are about  $\pm 1$  cm (worse than the normal normal  $\sim 1$  mm error bars for a tape measure, because the guinea pig was squirming). Bob then measures the same thing, and gets 12 cm. The average is computed as


$$c = \frac{A + B}{2},$$

where  $A$  is Alice’s measurement, and  $B$  is Bob’s, giving 11 cm. If Alice had been off by one standard deviation (1 cm), it would have changed the average by 0.5



cm, so we have  $Q_A = 0.5$  cm, and likewise  $Q_B = 0.5$  cm. Combining these, we find  $\sigma_c = \sqrt{Q_A^2 + Q_B^2} = 0.7$  cm, which is simply  $(1.0 \text{ cm})/\sqrt{2}$ . The final result is  $c = (11.0 \pm 0.7)$  cm. (This violates the usual rule for significant figures, which is that the final result should have no more sig figs than the least precise piece of data that went into the calculation. That's okay, because the sig fig rules are just a quick and dirty way of doing propagation of errors. We've done real propagation of errors in this example, and it turns out that the error is in the first decimal place, so the 0 in that place is entitled to hold its head high as a real sig fig, albeit a relatively imprecise one with an uncertainty of  $\pm 7$ .)

*Example: The difference between two measurements*  
In the example on page 85, we saw that two groups of scientists measured the same thing, and the results were  $W = 5.7 \pm 1.0$  for Webb et al. and  $C = 0.6 \pm 0.6$  for Chand et al. It's of interest to know whether the difference between their two results is small enough to be explained by random errors, or so big that it couldn't possibly have happened by chance, indicating that someone messed up. The figure shows each group's results, with error bars, on the number line. We see that the two sets of error bars don't overlap with one another, but error bars are not absolute limits, so it's perfectly possible to have non-overlapping error bars by chance, but the gap between the error bars is very large compared to the error bars themselves, so it looks implausible that the results could be statistically consistent with one another. I've tried to suggest this visually with the shading underneath the data-points.

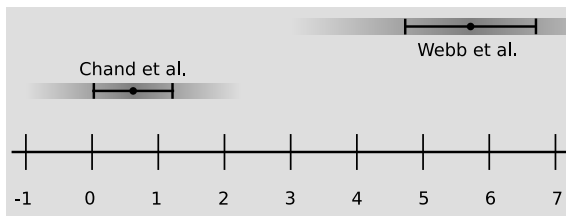

To get a sharper statistical test, we can calculate the difference  $d$  between the two results,


$$d = W - C$$

which is 5.1. Since the operation is simply the subtraction of the two numbers, an error in either input just causes an error in the output that is of the same size. Therefore we have  $Q_W = 1.0$  and  $Q_C = 0.6$ , resulting in  $\sigma_d = \sqrt{Q_W^2 + Q_C^2} = 1.2$ . We find that the difference between the two results is  $d = 5.1 \pm 1.2$ , which differs from zero by  $5.1/1.2 \approx 4$  standard deviations. Looking at the table on page 85, we see that the chances that  $d$  would be this big by chance are extremely small, less than about one in ten thousand. We can conclude to a high level of statistical confidence that the two groups' measurements are inconsistent with one another, and that one group is simply wrong.

# Appendix 4: Graphing

## Review of Graphing

Many of your analyses will involve making graphs. A graph can be an efficient way of presenting data visually, assuming you include all the information needed by the reader to interpret it. That means labeling the axes and indicating the units in parentheses, as in the example. A title is also helpful. Make sure that distances along the axes correctly represent the differences in the quantity being plotted. In the example, it would not have been correct to space the points evenly in the horizontal direction, because they were not actually measured at equally spaced points in time.

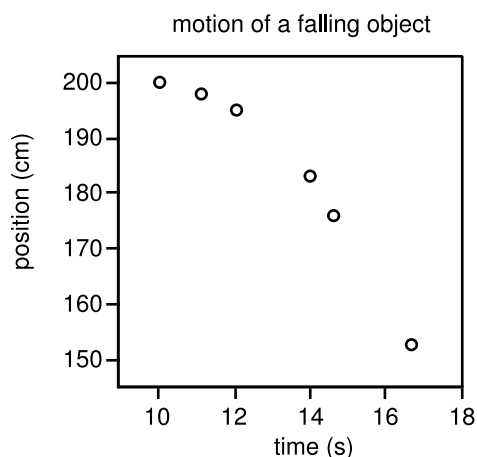

## Graphing on a Computer

Making graphs by hand in your lab notebook is fine, but in some cases you may find it saves you time to do graphs on a computer. For computer graphing, I recommend LibreOffice, which is free, open-source software. It's installed on the computers in rooms 416 and 418. Because LibreOffice is free, you can download it and put it on your own computer at home without paying money. If you already know Excel, it's very similar — you almost can't tell it's a different program.

Here's a brief rundown on using LibreOffice:

On Windows, go to the Start menu and choose All Programs, LibreOffice, and LibreOffice Calc. On Linux, do Applications, Office, OpenOffice, Spreadsheet.

Type in your x values in the first column, and your y values in the second column. For scientific notation, do, e.g., 5.2e-7 to represent  $5.2 \times 10^{-7}$ .

Select those two columns using the mouse.

From the Insert menu, do Object:Chart.


When it offers you various styles of graphs to choose from, choose the icon that shows a scatter plot, with dots on it (XY Chart).

Adjust the scales so the actual data on the plot is as big as possible, eliminating wasted space. To do this, double-click on the graph so that it's surrounded by a gray border. Then do Format, Axis, X Axis or Y Axis, Scale.

If you want error bars on your graph you can either draw them in by hand or put them in a separate column of your spreadsheet and doing Insert, Y Error Bars, Cell Range. Under Parameters, check "Same value for both." Click on the icon, and then use the mouse in the spreadsheet to select the cells containing the error bars.

## Fitting a Straight Line to a Graph by Hand

Often in this course you will end up graphing some data points, fitting a straight line through them with a ruler, and extracting the slope.

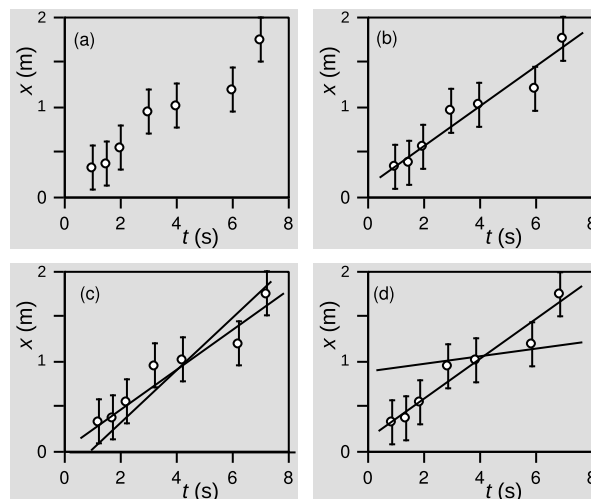

In this example, panel (a) shows the data, with error bars on each data point. Panel (b) shows a best fit, drawn by eye with a ruler. The slope of this best fit line is 100 cm/s. Note that the slope should be extracted from the line itself, not from two data points. The line is more reliable than any pair of individual data points.

In panel (c), a “worst believable fit” line has been drawn, which is as different in slope as possible from the best fit, while still pretty much staying consistent the data (going through or close to most of the error bars). Its slope is 60 cm/s. We can therefore estimate that the precision of our slope is +40 cm/s.

There is a tendency when drawing a “worst believable fit” line to draw instead an “unbelievably crazy fit” line, as in panel (d). The line in panel (d), with a very small slope, is just not believable compared to the data — it is several standard deviations away from most of the data points.

## Fitting a Straight Line to a Graph on a Computer

It’s also possible to fit a straight line to a graph using computer software such as LibreOffice.

To do this, first double-click on the graph so that a gray border shows up around it. Then right-click on a data-point, and a menu pops up. Choose Insert Trend Line.<sup>1</sup> choose Linear, and check the box for Show equation.

How accurate is your slope? A method for getting error bars on the slope is to artificially change one of your data points to reflect your estimate of how much it could have been off, and then redo the fit and find the new slope. The change in the slope tells you the error in the slope that results from the error in this data-point. You can then repeat this for the other points and proceed as in appendix 3.


An alternative method is to use the LINEST function that is available in many spreadsheet programs. For a description, see [tinyurl.com/ya7wmdft](http://tinyurl.com/ya7wmdft). Create the following formula in one cell of your spreadsheet: =Linest(y-values,x-value, True,True). Then, in excel, you need to press alt+ctrl+enter. In google sheets, press enter. A table with two columns and five rows will appear. The first number in the first column is the slope of the graph, and the second

number in the first column is the error in the slope.

In some cases, such as the absolute zero lab and the photoelectric effect lab, it’s very hard to tell how accurate your raw data are *a priori*; in these labs, you can use the typical amount of deviation of the points from the line as an estimate of their accuracy.

## Comparing Theory and Experiment

Figures (e) through (h) are examples of how we would compare theory and experiment on a graph. The convention is that theory is a line and experiment is points; this is because the theory is usually a prediction in the form of an equation, which can in principle be evaluated at infinitely many points, filling in all the gaps. One way to accomplish this with computer software is to graph both theory and experiment as points, but then print out the graph and draw a smooth curve through the theoretical points by hand.

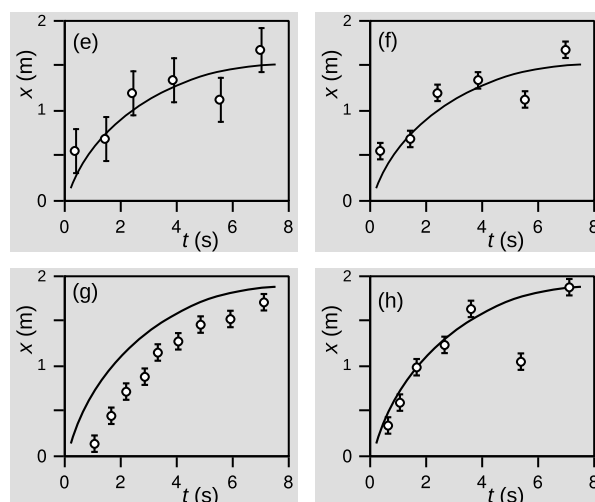

The point here is usually to compare theory and experiment, and arrive at a yes/no answer as to whether they agree. In (e), the theoretical curve goes through the error bars on four out of six of the data points. This is about what we expect statistically, since the probability of being within one standard deviation of the truth is about 2/3 for a standard bell curve. Given these data, we would conclude that theory and experiment agreed.

In graph (f), the points are exactly the same as in (e), but the conclusion is the opposite. The error bars are smaller, too small to explain the observed discrepancies between theory and experiment. The theoretical curve only goes through the error bars on two of the six points, and this is quite a bit less than

<sup>1</sup>“Trend line” is scientifically illiterate terminology that originates from Microsoft Office, which LibreOffice slavishly copies. If you don’t want to come off as an ignoramus, call it a “fit” or “line of best fit.”

we would expect statistically.

Graph (g) also shows disagreement between theory and experiment, but now we have a clear systematic error. In (h), the fifth data point looks like a mistake. Ideally you would notice during lab that something had gone wrong, and go back and check whether you could reproduce the result.



# Appendix 5: Finding Power Laws from Data

For many people, it is hard to imagine how scientists originally came up with all the equations that can now be found in textbooks. This appendix explains one method for finding equations to describe data from an experiment.

## Linear and nonlinear relationships

When two variables  $x$  and  $y$  are related by an equation of the form

$$y = cx \quad ,$$

where  $c$  is a constant (does not depend on  $x$  or  $y$ ), we say that a linear relationship exists between  $x$  and  $y$ . As an example, a harp has many strings of different lengths which are all of the same thickness and made of the same material. If the mass of a string is  $m$  and its length is  $L$ , then the equation

$$m = cL$$

will hold, where  $c$  is the mass per unit length, with units of kg/m. Many quantities in the physical world are instead related in a nonlinear fashion, i.e., the relationship does not fit the above definition of linearity. For instance, the mass of a steel ball bearing is related to its diameter by an equation of the form

$$m = cd^3 \quad ,$$

where  $c$  is the mass per unit volume, or density, of steel. Doubling the diameter does not double the mass, it increases it by a factor of eight.

## Power laws

Both examples above are of the general mathematical form

$$y = cx^p \quad ,$$

which is known as a power law. In the case of a linear relationship,  $p = 1$ . Consider the (made-up) experimental data shown in the table.

	$h$ =height of rodent at the shoulder (cm)	$f$ =food eaten per day (g)
shrew	1	3
rat	10	300
capybara	100	30,000

It's fairly easy to figure out what's going on just by staring at the numbers a little. Every time you increase the height of the animal by a factor of 10, its food consumption goes up by a factor of 100. This implies that  $f$  must be proportional to the square of  $h$ , or, displaying the proportionality constant  $k = 3$  explicitly,

$$f = 3h^2 \quad .$$

## Use of logarithms


Now we have found  $c = 3$  and  $p = 2$  by inspection, but that would be much more difficult to do if these weren't all round numbers. A more generally applicable method to use when you suspect a power-law relationship is to take logarithms of both variables. It doesn't matter at all what base you use, as long as you use the same base for both variables. Since the data above were increasing by powers of 10, we'll use logarithms to the base 10, but personally I usually just use natural logs for this kind of thing.

	$\log_{10} h$	$\log_{10} f$
shrew	0	0.48
rat	1	2.48
capybara	2	4.48

This is a big improvement, because differences are so much simpler to work mentally with than ratios. The difference between each successive value of  $h$  is 1, while  $f$  increases by 2 units each time. The fact that the logs of the  $f$ 's increase twice as quickly is the same as saying that  $f$  is proportional to the square of  $h$ .

## Log-log plots

Even better, the logarithms can be interpreted visually using a graph, as shown on the next page. The slope of this type of log-log graph gives the power  $p$ . Although it is also possible to extract the proportionality constant,  $c$ , from such a graph, the proportionality constant is usually much less interesting than  $p$ . For instance, we would suspect that if  $p = 2$  for rodents, then it might also equal 2 for frogs or ants. Also,  $p$  would be the same regardless of what units we used to measure the variables. The constant  $c$ , however, would be different if we used different units, and would also probably be different for other types of animals.

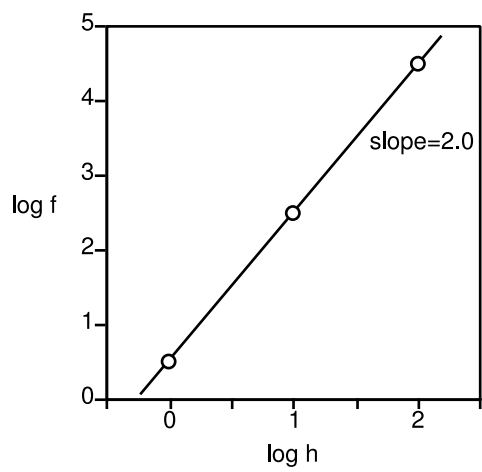

# Appendix 6: High Voltage Safety Checklist

Name: \_\_\_\_\_

----- Never work with high voltages by yourself.

----- Do not leave HV wires exposed - make sure there is insulation.

----- Turn the high-voltage supply off while working on the circuit.

----- When the voltage is on, avoid using both hands at once to touch the apparatus. Keep one hand in your pocket while using the other to touch the apparatus. That way, it is unlikely that you will get a shock across your chest.

----- It is possible for an electric current to cause your hand to clench involuntarily. If you observe this happening to your partner, do not try to pry their hand away, because you could become incapacitated as well — simply turn off the switch or pull the plug out of the wall.





# Appendix 7: Laser Safety Checklist

Name: \_\_\_\_\_

Before beginning a lab using lasers, make sure you understand these points, initial them, and show your safety checklist to your instructor. If you don't understand something, ask your instructor for clarification.

\_\_\_\_\_ The laser can damage your eyesight permanently if the beam goes in your eye.

\_\_\_\_\_ When you're not using the laser, turn it off or put something in front of it.

\_\_\_\_\_ Keep it below eye level and keep the beam horizontal. Don't bend or squat so that your eye is near the level of the beam.

\_\_\_\_\_ Keep the beam confined to your own lab bench. Whenever possible, orient your setup so that the beam is going toward the wall. If the beam is going to go off of your lab bench, use a backpack or a box to block the beam.

\_\_\_\_\_ Don't let the beam hit shiny surfaces such as faucets, because unpredictable reflections can result.

